

JENKEM TECHNOLOGY USA

INNOVATIVE MANUFACTURER
DEDICATED TO THE DEVELOPMENT
AND MANUFACTURE OF
HIGH-QUALITY PEG PRODUCTS

FROM R&D, CATALOG SCALE THROUGH
COMMERCIAL GMP PRODUCTION WITH
WORLDWIDE DISTRIBUTION

MOST ITEMS IN THIS CATALOG ARE OFFERED IN BULK
AT A SUBSTANTIAL DISCOUNT

CUSTOM SYNTHESIS CAPABILITIES ARE AVAILABLE
FOR PEG PRODUCTS NOT LISTED IN THIS CATALOG

PLEASE CONTACT US FOR PEGYLATION SERVICES

www.jenkemusa.com

Address:
4105 W. Spring Creek Pkwy, Suite 606B
Plano, TX 75024, U. S. A.

Phone:
866.953.6536 (866.9JENKEM)
972.673.0603

Fax:
469.327.2615

Email:
sales@jenkemusa.com

Prices are subject to change without notice.
Shipping and handling fees are extra.

TABLE OF CONTENTS

INTRODUCTION	5
NOVEL CELL PENETRATING PEPTIDE-PEG DERIVATIVES FOR GENE DELIVERY	6
NOVEL TECHNOLOGY FOR 3D-PRINTING OF PEG-BASED HYDROGELS	6
PRODUCTS FOR AMINE PEGYLATION	6
Y-shape PEG NHS Ester	7
Y-shape PEG Carboxyl	7
MPEG2 NHS Ester.....	7
Glucose PEG Succinimidyl Carboxymethyl Ester	8
Galactose PEG Succinimidyl Carboxymethyl Ester.....	8
Methoxy PEG Succinimidyl Carboxymethyl Ester	8
Methoxy PEG Acetic Acid	9
Methoxy PEG Succinimidyl Propionate.....	9
Methoxy PEG Propionic acid.....	10
Methoxy PEG Succinimidyl Butanoate	10
Methoxy PEG Butanoic acid	10
Methoxy PEG Succinimidyl Hexanoate.....	11
Methoxy PEG Hexanoic Acid.....	11
Methoxy PEG Succinimidyl Succinamide	11
Methoxy PEG Succinimidyl Glutaramide	12
Methoxy PEG Succinimidyl Carbonate	12
Methoxy PEG Nitrophenyl Carbonate.....	12
Methoxy PEG Succinimidyl Succinate	13
Methoxy PEG Succinimidyl Glutarate.....	13
PRODUCTS FOR THIOL PEGYLATION	14
Y-shape PEG Maleimide.....	15
MPEG2 Maleimide.....	15
Methoxy PEG Maleimide	15
Methoxy PEG Vinylsulfone	16
Methoxy PEG Thiol	16
PRODUCTS FOR N-TERMINAL PEGYLATION	18
Y-shape PEG Acetaldehyde	18
Y-shape PEG Propionaldehyde.....	18
Methoxy PEG Propionaldehyde	18
OTHER Y-SHAPE AND LINEAR METHOXY PEG DERIVATIVES	20
Y-shape PEG Alkyne	20
Y-shape PEG Amine	20
Y-shape PEG Fluorescein	20
Methoxy PEG Amine	21
Methoxy PEG Acrylate	21
Methoxy PEG Azide.....	21
Methoxy PEG Alkyne	22
Methoxy PEG Biotin	22
Methoxy PEG Hydrazide	22
Methoxy PEG DSPE.....	23
Methoxy PEG di-Glutamic Acid	23
Methoxy PEG tri-Glutamic Acid	23
Methoxy PEG Phosphate.....	24
Methoxy PEG Silane.....	24
Methoxy PEG Thiocctic acid	24
HOMOBIFUNCTIONAL PEG DERIVATIVES	26
Amine PEG Amine	26
Acrylate PEG Acrylate	26
Alkyne PEG Alkyne	27
Acetic Acid PEG Acetic Acid	27
Maleimide PEG Maleimide	27
PEG (Succinimidyl Carboxymethyl Ester)2	28
Thiol PEG Thiol.....	28
Vinylsulfone PEG Vinylsulfone	28
HETEROBIFUNCTIONAL PEG DERIVATIVES.....	30

Acrylate PEG Succinimidyl Carboxymethyl Ester	30
Acrylate PEG Succinimidyl Propionate	31
Amine PEG Acetic Acid	31
Alkyne PEG Maleimide	32
Azide PEG Succinimidyl Carboxymethyl Ester	32
Azide PEG Amine.....	32
Biotin PEG Maleimide.....	33
Biotin PEG Succinimidyl Carboxymethyl Ester	33
Biotin PEG Succinimidyl Glutaramide	34
Fmoc Amine PEG Amine	34
Fmoc Amine PEG Acetic Acid.....	34
Fmoc Amine PEG Hydroxyl	35
Fmoc Amine PEG Succinimidyl Carboxymethyl Ester	35
Hydroxyl PEG Amine.....	36
Hydroxyl PEG Acetic Acid	36
Hydroxyl PEG Succinimidyl Carboxymethyl Ester	37
Hydroxyl PEG Propionic Acid	37
Hydroxyl PEG Succinimidyl Propionate	38
Hydroxyl PEG Hexanoic Acid	38
Thiol PEG Amine	38
Thiol PEG Acetic Acid	39
Thiol PEG Succinimidyl Glutaramide.....	39
Thiol PEG Succinimidyl Propionate	39
t-Boc Amine PEG Hydroxyl	40
t-Boc Amine PEG Amine	40
t-Boc Amine PEG Acetic Acid.....	41
t-Boc Amine PEG Succinimidyl Carboxymethyl Ester	41
Maleimide PEG Amine	42
Maleimide PEG Acetic Acid.....	42
Maleimide PEG Hydroxyl	43
Maleimide PEG Succinimidyl Carboxymethyl Ester	43
OPSS PEG Amine	43
OPSS PEG Succinimidyl Carboxymethyl Ester.....	44
MONODISPERSE PEG DERIVATIVES.....	46
METHOXY MONODISPERSE PEG DERIVATIVES	46
Methoxy-PEG(5 or 6)-Amine	46
Methoxy-PEG(5 or 6)-Acrylate	47
Methoxy-PEG(5 or 6)-Azide	47
Methoxy-PEG(5 or 6)-Hydroxyl	47
Methoxy-PEG(5 or 6)-Propionic acid	47
Methoxy-PEG(5 or 6)-Thiol	48
HOMOBIFUNCTIONAL MONODISPERSE PEG MONODISPERSE DERIVATIVES	48
Amine PEG(8 or 12) Amine.....	48
Acetic Acid PEG(8 or 12) Acetic Acid	48
Hydroxyl-PEG(8 or 9 or 12)-Hydroxyl	49
Propionic Acid PEG(8 or 12) Propionic Acid	49
HETEROBIFUNCTIONAL MONODISPERSE PEG MONODISPERSE DERIVATIVES	49
Hydroxyl PEG(8 or 12) Amine	49
Hydroxyl PEG(8 or 12) Acetic Acid	50
Hydroxyl PEG(8 or 12) Propionic Acid	50
MULTI-ARM PEG DERIVATIVES	50
3ARM-PEG DERIVATIVES (GLYCEROL CORE)	51
3arm-PEG-Amine	51
4ARM-PEG DERIVATIVES (PENTAERYTHRITOL CORE)	51
4arm PEG Amine (pentaerythritol).....	51
4arm PEG Acrylate (pentaerythritol).....	51
4arm PEG Azide (pentaerythritol).....	52
4arm PEG Acetic Acid (pentaerythritol)	52
4arm PEG Isocyanate (pentaerythritol)	52
4arm PEG Maleimide (pentaerythritol)	53
4arm PEG Nitrophenyl Carbonate (pentaerythritol)	53

4arm PEG Succinimidyl Carboxymethyl Ester (pentaerythritol).....	53
4arm PEG Succinimidyl Glutamide (pentaerythritol).....	54
4arm PEG Succinimidyl Carbonate (pentaerythritol).....	54
4arm PEG Succinimidyl Succinate (pentaerythritol)	55
4arm PEG Succinimidyl Glutarate (pentaerythritol)	55
4arm PEG Thiol (pentaerythritol).....	55
4arm PEG Vinylsulfone (pentaerythritol)	56
6ARM-PEG DERIVATIVES (DIPENTAERYTHRITOL CORE)	56
6arm PEG Amine (dipentaerythritol)	56
8ARM-PEG DERIVATIVES (HEXAGLYCEROL CORE).....	57
8arm PEG Amine (hexaglycerol).....	57
8arm PEG Acrylate (hexaglycerol)	57
8arm PEG Acetic Acid (hexaglycerol)	57
8arm PEG Maleimide (hexaglycerol)	58
8arm PEG Succinimidyl Succinate (hexaglycerol).....	58
8arm PEG Succinimidyl Glutarate (hexaglycerol)	58
8arm PEG Thiol (hexaglycerol).....	59
8ARM-PEG DERIVATIVES (TRIPENTAERYTHRITOL CORE)	59
8arm PEG Acrylate (tripentaerythritol)	59
8arm PEG Amine (tripentaerythritol).....	59
8arm PEG Acetic Acid (tripentaerythritol)	60
8arm PEG Maleimide (tripentaerythritol)	60
8arm PEG Norbornene (tripentaerythritol)	60
8arm PEG Succinimidyl Succinate (tripentaerythritol)	61
8arm PEG Succinimidyl Glutarate (tripentaerythritol).....	61
8arm PEG Thiol (tripentaerythritol).....	62
8arm PEG Vinylsulfone (tripentaerythritol)	62
MULTI-ARM HETEROBIFUNCTIONAL PEG DERIVATIVES	64
(HO)3-4ARMPEG-CM.....	64
(HO)3-4ARMPEG-NH2HCl	64
(HO)2-4ARMPEG-(CM)2	65
(HO)2-4ARMPEG-(NH2HCl)2	65
HO-4ARMPEG-(CM)3.....	66
(CH ₃ O)3-4ARMPEG-SCM.....	66
(HO)7-8ARMPEG-CM.....	66
(HO)7-8ARMPEG(TP)-CM.....	67
BIODEGRADABLE PEG CO-POLYMERS	68
PEG-POLYLACTIDE BLOCK COPOLYMERS	68
4arm PEG Polylactide Block Co-Polymers.....	68
Methoxy PEG Polylactide Block Co-Polymers	68
Maleimide PEG Polylactide Block Co-Polymers	69
PEG-POLYCARPROLACTONE BLOCK COPOLYMERS	69
4arm PEG Polycaprolactone Block Co-Polymers	69
Methoxy PEG Polycaprolactone Block Co-Polymers	69
PEG RAW MATERIALS	70
LINEAR METHOXY PEG RAW MATERIALS	71
LINEAR BENZYL PEG RAW MATERIALS.....	71
MULTI-ARM PEG RAW MATERIALS	71
4ARM-PEG Raw Materials (pentaerythritol core)	71
6ARM-PEG Raw Materials (dipentaerythritol core).....	72
8ARM-PEG Raw Materials (hexaglycerol core)	72
8ARM-PEG Raw Materials (tripentaerythritol core)	72
PEG GPC STANDARDS	74
ORDERING INFORMATION.....	75
STORAGE AND HANDLING OF PEG DERIVATIVES.....	76

INTRODUCTION

Founded in 2001, JenKem Technology specializes in high-purity, low polydispersity polyethylene glycol (PEG) derivatives, PEG co-polymers, monodisperse PEGs, custom PEG derivative synthesis and PEGylation services. Specific strengths are supported by access to bulk ethylene oxide and the performance of all polymerizations and derivatizations in-house.

JenKem Technology supplies PEG derivatives worldwide, from R&D through GMP commercial quantities, for preclinical, clinical trials, and commercial products in the pharmaceutical, biotechnology, medical device, and diagnostics marketplaces.

JenKem Technology manufactures high quality reactive polyethylene glycol (PEG) derivatives with high purity and narrow polydispersity, from laboratory scale to commercial scale in both GMP and non-GMP grade. PEG raw materials are provided in bulk. JenKem Technology's PEG raw materials are back integrated to production directly from ethylene oxide. GPC calibration standards, necessary for accurate GPC analysis of PEG products, are provided in GPC calibration kits for linear PEGs and 4arm PEGs.

PEGs have been shown to provide improved water solubility and biocompatibility in drugs and other applications, and the ability to attach a variety of reactive functional groups to the terminal positions of these polymers has greatly increased their applications. Homo- or hetero-bifunctional PEGs are especially useful as cross-linking agents or spacers between two chemical entities, whereas mono-functional PEGs prevent bridging reactions which can affect the PEGylation of certain compounds. Branched PEG derivatives consist of linear methoxy PEG chains attached to a central core. JenKem® proprietary Y-shape PEGs are more selective, due to their sterically bulky structure. Multi-arm PEG derivatives, prepared by ethoxylation of different cores such as pentaerythritol, hexaglycerol, or tripentaerythritol, can be cross-linked in a variety of ways into hydrogels.

Bi-functional PEG derivatives with different terminal functional groups are widely used in the PEGylation of peptides, proteins, small molecules, like folate, mannose, pro-drugs or oligonucleotides, cells, nanoparticles and virus particles, or surfaces, and require extensive characterization with complex analytical techniques. Multi-arm PEG derivatives are largely employed in the formation of hydrogels. PEG hydrogels are used extensively in the controlled release of therapeutics, in medical devices, and regenerative medicine, and in various other applications, including wound sealing and healing, 3D printing of biosimilar materials, or as cartilage replacements.


Functionally, JenKem Technology's PEG products for PEGylation include activated PEG derivatives for amine PEGylation, activated PEG derivatives for N-terminal PEGylation, and activated peg derivatives for thiol PEGylation. Structurally, JenKem Technology's PEG products include JenKem® proprietary Y-shape branched PEG derivatives; linear M-PEG derivatives with low diol content; homobifunctional PEG derivatives; heterobifunctional PEG derivatives; multi-arm homofunctional reactive PEGs for hydrogels, multi-arm heterofunctional PEG derivatives for ADC development, and biodegradable PEG co-polymers.

The following innovative PEG product lines launched in 2018 are currently offered as made-to-order items. Please contact us for details on these novel products:

NOVEL CELL PENETRATING PEPTIDE-PEG DERIVATIVES FOR GENE DELIVERY

- Novel activated Cell Penetrating Peptide-PEG derivatives for covalent coupling
- Application for intracellular delivery of molecules that cannot penetrate the cell membrane on their own
- Presence of PEG eliminates aggregation and sedimentation of RNA or DNA
- Cleavable linkage for cytosolic release of delivered molecule
- In vitro demonstrated cellular uptake of siRNA and high gene-silencing efficacy of delivered siRNA

Schematic of novel CPP-PEG drug delivery system for RNA delivery:

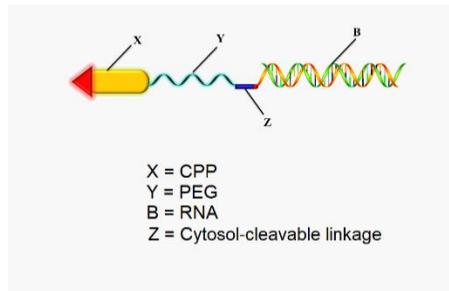

NOVEL TECHNOLOGY FOR 3D-PRINTING OF PEG-BASED HYDROGELS

- M-QTX novel photo-initiator with good water-solubility and biocompatibility, and low cytotoxicity
- Multiarm PEG and PEG-PCL acrylate and methacrylate derivatives for hydrogels
- Hydrogels 3D printed through SLA (stereolithography) at wavelength range of 380-420nm are transparent with well-defined shape
- Low toxicity, biocompatible hydrogels are suitable for tissue engineering and 3D cell culture

Pattern 3D printed from aqueous solution of 8ARM(TP)-ACLT-10K and M-QTX photoinitiator:




PRODUCTS FOR AMINE PEGYLATION

JenKem Technology provides high quality polyethylene glycol (PEG) derivatives for Amine PEGylation with high purity and low polydispersity.

JenKem Technology offers amine reactive PEGs with linear structure, Y-shape PEGs, U-shape PEGs, multiarm, and monodisperse PEGs. JenKem® proprietary Y-shape PEGs are more selective, due to their sterically bulky structure. JenKem Technology manufactures numerous NHS and Carboxyl PEGs for amine PEGylation. JenKem Technology's NHS PEGs are categorized as PEG NHS Esters with stable linkers (such as M-SCM, M-SPA, M-SBA, M-SHA, M-SGA, M-SSA); PEG NHS Esters with carbonate linker (such as M-SC and M-NPC); and degradable PEGs with cleavable NHS ester linker (such as M-SS and M-SG).


Amine PEGylation products with molecular weights, branching, and functional groups not listed in our catalog may be available by custom synthesis. GMP grade bulk products may also be available by custom synthesis. Please inquire at tech@jenkemusa.com about pricing and availability of custom activated PEGs.

Y-SHAPE PEG NHS ESTER

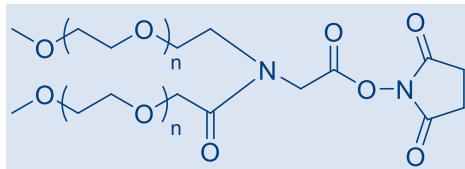


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A0001-1	Y-NHS-40K,1g	Y-shape PEG NHS Ester, MW 40000	\$130
A0001-10	Y-NHS-40K,10g	Y-shape PEG NHS Ester, MW 40000	\$1,040

Y-SHAPE PEG CARBOXYL

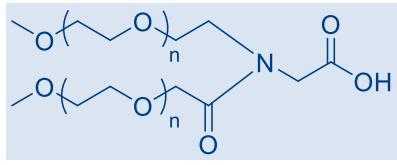

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A0019-1	Y-COOH-40K,1g	Y-shape PEG Carboxyl, MW 40000	\$120
A0019-10	Y-COOH-40K,10g	Y-shape PEG Carboxyl, MW 40000	\$960

MPEG2 NHS ESTER

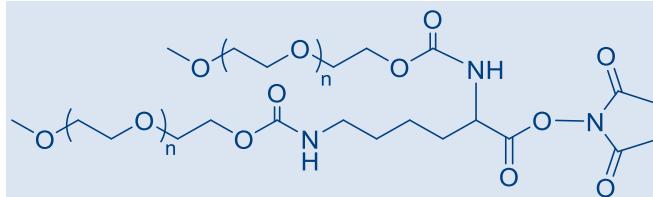


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A1003-1	MPEG2-NHS-40K,1g	MPEG2 NHS Ester, MW 40000	\$160
A1003-10	MPEG2-NHS-40K,10g	MPEG2 NHS Ester, MW 40000	\$1,280

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A1001-1	MPEG2-NHS-20K,1g	MPEG2 NHS Ester, MW 20000	\$170
*A1001-10	MPEG2-NHS-20K,10g	MPEG2 NHS, MW 20000	\$1,360


Please inquire if the PEG derivative that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.

GLUCOSE PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

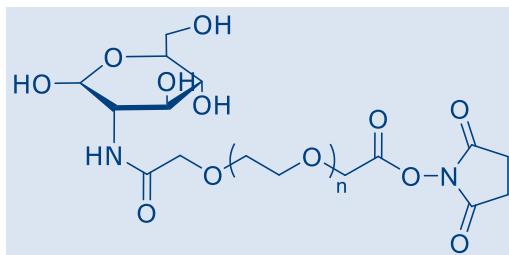

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A2001-1	GLUC-SCM-35K,1g	Glucose PEG Succinimidyl Carboxymethyl Ester, MW 35000	\$150
A2001-10	GLUC-SCM-35K,10g	Glucose PEG Succinimidyl Carboxymethyl Ester, MW 35000	\$1,200


GALACTOSE PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

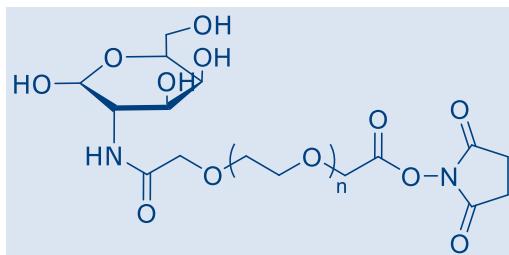

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A2002-1	GALA-SCM-35K,1g	Galactose PEG Succinimidyl Carboxymethyl Ester, MW 35000	\$200
*A2002-10	GALA-SCM-35K,10g	Galactose PEG Succinimidyl Carboxymethyl Ester, MW 35000	\$1,600

METHOXY PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

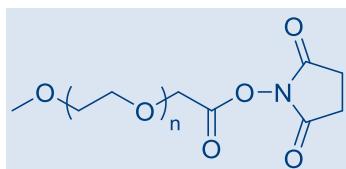


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3072-1	M-SCM-2000,1g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$110
A3072-10	M-SCM-2000,10g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$880
A3105-1	M-SCM-3000,1g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 3000	\$110
A3105-10	M-SCM-3000,10g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 3000	\$880
A3012-1	M-SCM-5000,1g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$110
A3012-10	M-SCM-5000,10g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$880
A3016-1	M-SCM-10K,1g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 10000	\$110
A3016-10	M-SCM-10K,10g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 10000	\$880
A3003-1	M-SCM-20K,1g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 20000	\$110
A3003-10	M-SCM-20K,10g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 20000	\$880

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3041-1	M-SCM-40K,1g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 40000	\$130
A3041-10	M-SCM-40K,10g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 40000	\$1,040

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3004-1	M-SCM-30K,1g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 30000	\$110
*A3004-10	M-SCM-30K,10g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 30000	\$880
*A3005-1	M-SCM-35K,1g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 35000	\$110
*A3005-10	M-SCM-35K,10g	Methoxy PEG Succinimidyl Carboxymethyl Ester, MW 35000	\$880

METHOXY PEG ACETIC ACID

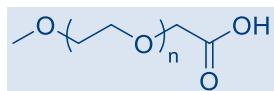

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3070-1	M-CM-5000,1g	Methoxy PEG Acetic Acid, MW 5000	\$100
A3070-10	M-CM-5000,10g	Methoxy PEG Acetic Acid, MW 5000	\$800
A3081-1	M-CM-10K,1g	Methoxy PEG Acetic Acid, MW 10000	\$100
A3081-10	M-CM-10K,10g	Methoxy PEG Acetic Acid, MW 10000	\$800
A3089-1	M-CM-20K,1g	Methoxy PEG Acetic Acid, MW 20000	\$100
A3089-10	M-CM-20K,10g	Methoxy PEG Acetic Acid, MW 20000	\$800

METHOXY PEG SUCCINIMIDYL PROPIONATE

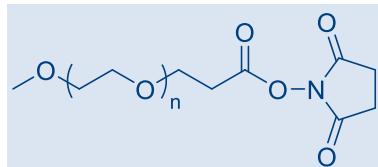


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3013-1	M-SPA-5000,1g	Methoxy PEG Succinimidyl Propionate, MW 5000	\$160
A3013-10	M-SPA-5000,10g	Methoxy PEG Succinimidyl Propionate, MW 5000	\$1,280

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3019-1	M-SPA-20K,1g	Methoxy PEG Succinimidyl Propionate, MW 20000	\$160
*A3019-10	M-SPA-20K,10g	Methoxy PEG Succinimidyl Propionate, MW 20000	\$1,280
*A3107-1	M-SPA-30K,1g	Methoxy PEG Succinimidyl Propionate, MW 30000	\$170
*A3107-10	M-SPA-30K,10g	Methoxy PEG Succinimidyl Propionate, MW 30000	\$1,360

METHOXY PEG PROPIONIC ACID

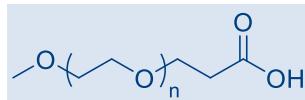


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3225-1	M-PA-5000,1g	Methoxy PEG Propionic Acid, MW 5000	\$150
A3225-10	M-PA-5000,10g	Methoxy PEG Propionic Acid, MW 5000	\$1,200

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A6023-1	M-PA-20K,1g	Methoxy PEG Propionic Acid, MW 20000	\$150
*A6023-10	M-PA-20K,10g	Methoxy PEG Propionic Acid, MW 20000	\$1,200
*A3226-1	M-PA-30K,1g	Methoxy PEG Propionic Acid, MW 30000	\$160
*A3226-10	M-PA-30K,10g	Methoxy PEG Propionic Acid, MW 30000	\$1,280


METHOXY PEG SUCCINIMIDYL BUTANOATE

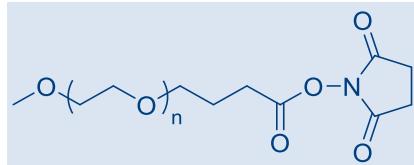

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3144-1	M-SBA-20K,1g	Methoxy PEG Succinimidyl Butanoate, MW 20000	\$140
*A3144-10	M-SBA-20K,10g	Methoxy PEG Succinimidyl Butanoate, MW 20000	\$1,120


METHOXY PEG BUTANOIC ACID

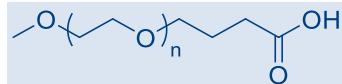

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3227-1	M-BA-20K,1g	Methoxy PEG Butanoic Acid, MW 20000	\$130
*A3227-10	M-BA-20K,10g	Methoxy PEG Butanoic Acid, MW 20000	\$1,040

METHOXY PEG SUCCINIMIDYL HEXANOATE

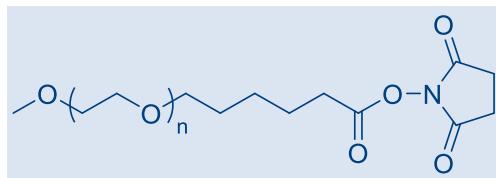


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3148-1	M-SHA-20K,1g	Methoxy PEG Succinimidyl Hexanoate, MW 20000	\$120
A3148-10	M-SHA-20K,10g	Methoxy PEG Succinimidyl Hexanoate, MW 20000	\$960
A3103-1	M-SHA-30K,1g	Methoxy PEG Succinimidyl Hexanoate, MW 30000	\$120
A3103-10	M-SHA-30K,10g	Methoxy PEG Succinimidyl Hexanoate, MW 30000	\$960

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A3146-1	M-SHA-5000,1g	Methoxy PEG Succinimidyl Hexanoate, MW 5000	\$120
* A3146-10	M-SHA-5000,10g	Methoxy PEG Succinimidyl Hexanoate, MW 5000	\$960

METHOXY PEG HEXANOIC ACID

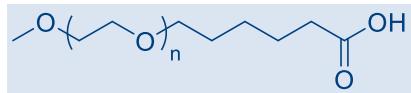


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3147-1	M-HA-20K,1g	Methoxy PEG Hexanoic Acid, MW 20000	\$100
A3147-10	M-HA-20K,10g	Methoxy PEG Hexanoic Acid, MW 20000	\$800
A3259-1	M-HA-30K,1g	Methoxy PEG Hexanoic Acid, MW 30000	\$100
A3259-10	M-HA-30K,10g	Methoxy PEG Hexanoic Acid, MW 30000	\$800

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A3145-1	M-HA-5000,1g	Methoxy PEG Hexanoic Acid, MW 5000	\$100
* A3145-10	M-HA-5000,10g	Methoxy PEG Hexanoic Acid, MW 5000	\$800


METHOXY PEG SUCCINIMIDYL SUCCINAMIDE

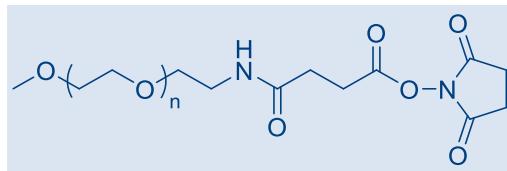

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3143-1	M-SSA-5000,1g	Methoxy PEG Succinimidyl Succinamide, MW 5000	\$120
*A3143-10	M-SSA-5000,10g	Methoxy PEG Succinimidyl Succinamide, MW 5000	\$960
*A3141-1	M-SSA-20K,1g	Methoxy PEG Succinimidyl Succinamide, MW 20000	\$120
*A3141-10	M-SSA-20K,10g	Methoxy PEG Succinimidyl Succinamide, MW 20000	\$960


METHOXY PEG SUCCINIMIDYL GLUTARAMIDE

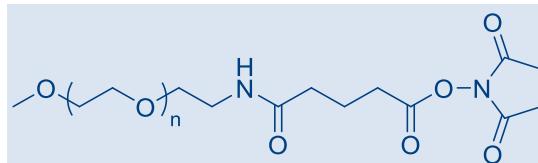

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3050-1	M-SGA-5000,1g	Methoxy PEG Succinimidyl Glutaramide, MW 5000	\$120
*A3050-10	M-SGA-5000,10g	Methoxy PEG Succinimidyl Glutaramide, MW 5000	\$960
*A3040-1	M-SGA-20K,1g	Methoxy PEG Succinimidyl Glutaramide, MW 20000	\$120
*A3040-10	M-SGA-20K,10g	Methoxy PEG Succinimidyl Glutaramide, MW 20000	\$960

METHOXY PEG SUCCINIMIDYL CARBONATE

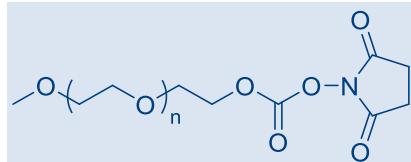

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3083-1	M-SC-10K,1g	Methoxy PEG Succinimidyl Carbonate, MW 10000	\$80
A3083-10	M-SC-10K,10g	Methoxy PEG Succinimidyl Carbonate, MW 10000	\$640
A3061-1	M-SC-12K,1g	Methoxy PEG Succinimidyl Carbonate, MW 12000	\$80
A3061-10	M-SC-12K,10g	Methoxy PEG Succinimidyl Carbonate, MW 12000	\$640
A3021-1	M-SC-20K,1g	Methoxy PEG Succinimidyl Carbonate, MW 20000	\$80
A3021-10	M-SC-20K,10g	Methoxy PEG Succinimidyl Carbonate, MW 20000	\$640


METHOXY PEG NITROPHENYL CARBONATE

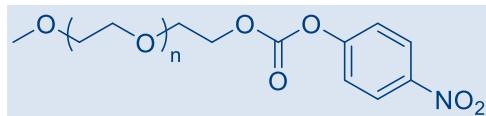

Please inquire if the PEG derivative that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3094-1	M-NPC-10K,1g	Methoxy PEG Nitrophenyl Carbonate, MW 10000	\$90
*A3094-10	M-NPC-10K,10g	Methoxy PEG Nitrophenyl Carbonate, MW 10000	\$720
*A3038-1	M-NPC-20K,1g	Methoxy PEG Nitrophenyl Carbonate, MW 20000	\$90
*A3038-10	M-NPC-20K,10g	Methoxy PEG Nitrophenyl Carbonate, MW 20000	\$720

METHOXY PEG SUCCINIMIDYL SUCCINATE

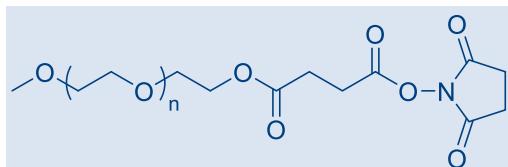

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3011-1	M-SS-5000,1g	Methoxy PEG Succinimidyl Succinate, MW 5000	\$80
A3011-10	M-SS-5000,10g	Methoxy PEG Succinimidyl Succinate, MW 5000	\$640
A3022-1	M-SS-20K,1g	Methoxy PEG Succinimidyl Succinate, MW 20000	\$80
A3022-10	M-SS-20K,10g	Methoxy PEG Succinimidyl Succinate, MW 20000	\$640

METHOXY PEG SUCCINIMIDYL GLUTARATE

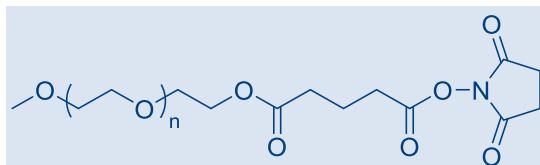

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3101-1	M-SG-5000,1g	Methoxy PEG Succinimidyl Glutarate, MW 5000	\$80
A3101-10	M-SG-5000,10g	Methoxy PEG Succinimidyl Glutarate, MW 5000	\$640
A3031-1	M-SG-20K,1g	Methoxy PEG Succinimidyl Glutarate, MW 20000	\$80
A3031-10	M-SG-20K,10g	Methoxy PEG Succinimidyl Glutarate, MW 20000	\$640

REFERENCES:

1. Xiong, Q., et al., Facile fabrication of reduction-responsive supramolecular nanoassemblies for co-delivery of doxorubicin and sorafenib towards hepatoma cells, *Frontiers in pharmacology*, 2018, 9, p.61.
2. Mastorakos, P., et al., Biodegradable brain-penetrating DNA nanocomplexes and their use to treat malignant brain tumors, *Journal of Controlled Release*, 2017, 262, P. 37-46.
3. Vaillard, V.A., et al., MPEG-NHS carbonates: Effect of alkyl spacers on the reactivity: Kinetic and mechanistic insights, *Journal of Applied Polymer Science*, 2019, 136(5):47028.
4. Luan, J., et al., GSDMD membrane pore is critical for IL-1 β release and antagonizing IL-1 β by hepatocyte-specific nanobiologics is a promising therapeutics for murine alcoholic steatohepatitis, *Biomaterials*, 2020, V. 227.
5. Sousa, A., et al., Design of experiments to select triphenylphosphonium-polyplexes with suitable physicochemical properties for mitochondrial gene therapy, *Journal of Molecular Liquids*, 2020, V. 302.
6. Chen, L., Ultra-small MoS₂ nanodots-incorporated mesoporous silica nanospheres for pH-sensitive drug delivery and CT imaging, *Journal of Materials Science & Technology*, 2021, V. 63, P. 91-96.

7. Yan, J., et al., Tumor Contrast Imaging with Gas Vesicles by Circumventing the Reticuloendothelial System, Ultrasound in Medicine & Biology, 2020, V. 46(2); P. 359-368.
8. Mesken, J., et al., Modifying plasmid-loaded HSA-nanoparticles with cell penetrating peptides–Cellular uptake and enhanced gene delivery, International journal of pharmaceutics, 2017, 522(1):198-209.
9. Li, C., et al., Real-Time Monitoring Surface Chemistry-Dependent In Vivo Behaviors of Protein Nanocages via Encapsulating an NIR-II Ag2S Quantum Dot, ACS Nano, 2015, 9 (12), 12255-12263.
10. Suarez, et I., Intramyocardial injection of hydrogel with high interstitial spread does not impact action potential propagation, Acta Biomaterialia, 2015, V. 26, P. 13-22.
11. Pfister, D., et al., Integrated process for high conversion and high yield protein PEGylation, Biotechnol. Bioeng., 2016.
12. Goffin V, et al., Pegvisomant Pfizer/Sensus. Curr Opin Investig Drugs, 2004, 5:463-468.
13. Liu, Y., et al., Cathodic protected Mn2+ by NaxWO3 nanorods for stable magnetic resonance imaging-guided tumor photothermal therapy, Biomaterials, 2020, V. 234.
14. Javanmardi, S., et al., Redox-Sensitive, PEG-Shielded Carboxymethyl PEI Nanogels Silencing MicroRNA-21, Sensitizes Resistant Ovarian Cancer Cells to Cisplatin, Asian Journal of Pharmaceutical Sciences, 2018.
15. Jones, S.K., et al., Revisiting the value of competition assays in folate receptor-mediated drug delivery, Biomaterials, 2017.
16. Zhang, C., Impact of Large Aggregated Uricases and PEG Diol on Accelerated Blood Clearance of PEGylated Canine Uricase, PLoS One, 2012, 7(6): e39659.
17. Wang, C., et al., Synthesis and formation mechanism of bone mineral, whitlockite nanocrystals in tri-solvent system, Journal of Colloid and Interface Science, 2020, V. 569, P. 1-11.
18. Ji, F., et al., A Dual pH/Magnetic Responsive Nanocarrier Based on PEGylated Fe3O4 Nanoparticles for Doxorubicin Delivery, Journal of Nanoscience and Nanotechnology, 2018, 18.7: 4464-4470.
19. Xiao, S., et al., Aptamer-mediated gene therapy enhanced antitumor activity against human hepatocellular carcinoma in vitro and in vivo, Journal of Controlled Release, 2017.
20. Jing, P., et al, Enhanced growth inhibition of prostate cancer in vitro and in vivo by a recombinant adenovirus-mediated dual-aptamer modified drug delivery system, Cancer Letters, 2016, V. 383(2), P. 230-242.
21. Lu, C., et al., Poly(ethylene glycol) crosslinked multi-armed poly(ϵ -benzyloxycarbonyl-L-lysine)s as super-amphiphiles: Synthesis, self-assembly, and evaluation as efficient delivery systems for poorly water-soluble drugs, Colloids and Surfaces B: Biointerfaces, 2019, 182.
22. AlQahtani, A.D., et al., Production of “biobetter” glucarpidase variants to improve drug detoxification and antibody directed enzyme prodrug therapy for cancer treatment, European Journal of Pharmaceutical Sciences, 2019, V. 127, P. 79-91.
23. Haruta, K., et al., A Novel PEGylation Method for Improving the Pharmacokinetic Properties of Anti-Interleukin-17A RNA Aptamers, Nucleic acid therapeutics, 2017, 27(1):36-44.
24. Ma, G., PEGylated drugs: Concept, Design and Application, Science Press, 2016, 314.
25. Ma, S.S., et al., The pharmacokinetic and pharmacodynamic properties of site-specific pegylated genetically modified recombinant human interleukin-11 in normal and thrombocytopenic monkeys, European Journal of Pharmaceutics and Biopharmaceutics, 2017.
26. Oddone, N., et al., ROS-responsive “smart” polymeric conjugate: Synthesis, characterization and proof-of-concept study, International Journal of Pharmaceutics, 2019, 570.


PRODUCTS FOR THIOL PEGYLATION

JenKem Technology provides high quality activated polyethylene glycol derivatives (PEGs) for Thiol PEGylation with high purity and low polydispersity.

Methoxy PEG Maleimide is selective for thiol groups on cysteine side chains. Y-shape PEG Maleimide is more selective, due to its sterically bulky structure. Methoxy PEG Vinylsulfone is a sulphydryl reactive PEG that reacts at high pH. Methoxy PEG Thiol polydisperse or monodisperse PEGs are reactive towards Thiol groups on cysteine side chains under mild reaction conditions. PolyLactide PEG Maleimide (PLA PEG Maleimide) is a biodegradable thiol reactive PEG copolymer, usually employed for controlled drug release.


Thiol PEGylation products with molecular weights, branching, and functional groups not listed in our catalog may be available by custom synthesis. GMP grade bulk products may also be available by custom synthesis. Please inquire at tech@jenkemusa.com about pricing and availability of custom activated PEGs.

Y-SHAPE PEG MALEIMIDE

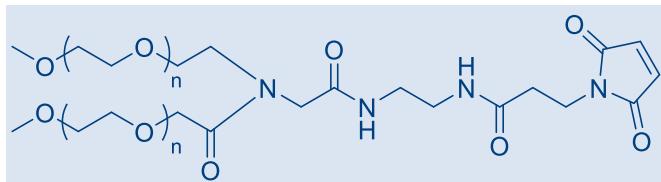


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A0002-1	Y-MAL-40K,1g	Y-shape PEG Maleimide, MW 40000	\$170
A0002-10	Y-MAL-40K,10g	Y-shape PEG Maleimide, MW 40000	\$1,360

MPEG2 MALEIMIDE

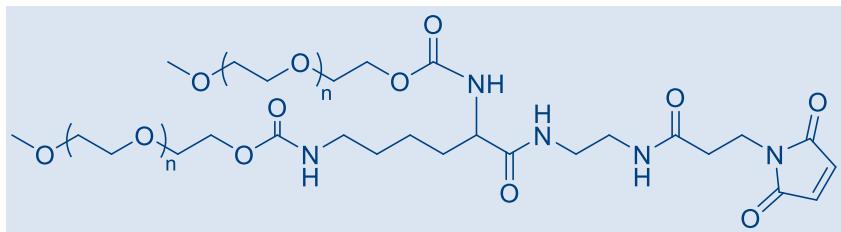


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A1009-1	MPEG2-MAL-40K,1g	MPEG2 Maleimide, MW 40000	\$180
A1009-10	MPEG2-MAL-40K,10g	MPEG2 Maleimide, MW 40000	\$1,440

METHOXY PEG MALEIMIDE

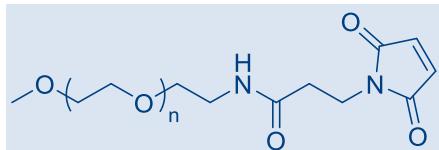

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3124-1	M-MAL-2000,1g	Methoxy PEG Maleimide, MW 2000	\$130
A3124-10	M-MAL-2000,10g	Methoxy PEG Maleimide, MW 2000	\$1,040
A3125-1	M-MAL-5000,1g	Methoxy PEG Maleimide, MW 5000	\$130
A3125-10	M-MAL-5000,10g	Methoxy PEG Maleimide, MW 5000	\$1,040
A3119-1	M-MAL-10K,1g	Methoxy PEG Maleimide, MW 10000	\$130
A3119-10	M-MAL-10K,10g	Methoxy PEG Maleimide, MW 10000	\$1,040
A3115-1	M-MAL-20K,1g	Methoxy PEG Maleimide, MW 20000	\$130
A3115-10	M-MAL-20K,10g	Methoxy PEG Maleimide, MW 20000	\$1,040
A3116-1	M-MAL-30K,1g	Methoxy PEG Maleimide, MW 30000	\$130
A3116-10	M-MAL-30K,10g	Methoxy PEG Maleimide, MW 30000	\$1,040
A3123-1	M-MAL-40K,1g	Methoxy PEG Maleimide, MW 40000	\$160
A3123-10	M-MAL-40K,10g	Methoxy PEG Maleimide, MW 40000	\$1,280


METHOXY PEG VINYLSULFONE

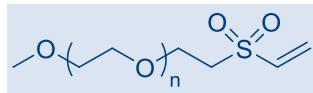

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3034-1	M-VS-5000,1g	Methoxy PEG Vinylsulfone, MW 5000	\$130
*A3034-10	M-VS-5000,10g	Methoxy PEG Vinylsulfone, MW 5000	\$1,040
*A3006-1	M-VS-20K,1g	Methoxy PEG Vinylsulfone, MW 20000	\$130
*A3006-10	M-VS-20K,10g	Methoxy PEG Vinylsulfone, MW 20000	\$1,040

METHOXY PEG THIOL

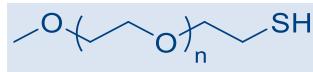

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3029-1	M-SH-5000,1g	Methoxy PEG Thiol, MW 5000	\$120
A3029-10	M-SH-5000,10g	Methoxy PEG Thiol, MW 5000	\$960

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3007-1	M-SH-20K,1g	Methoxy PEG Thiol, MW 20000	\$120
*A3007-10	M-SH-20K,10g	Methoxy PEG Thiol, MW 20000	\$960
*A3095-1	M-SH-30K,1g	Methoxy PEG Thiol, MW 30000	\$150
*A3095-10	M-SH-30K,10g	Methoxy PEG Thiol, MW 30000	\$1,200
*A3051-1	M-SH-40K,1g	Methoxy PEG Thiol, MW 40000	\$170
*A3051-10	M-SH-40K,10g	Methoxy PEG Thiol, MW 40000	\$1,360

REFERENCES:

- Li, H., et al., Dual MMP7-Proximity-Activated and Folate Receptor-Targeted Nanoparticles for siRNA Delivery, *Biomacromolecules*, 2015, 16 (1), p: 192–201.
- Umeshappa, C. S., et al., Liver-specific T regulatory type-1 cells program local neutrophils to suppress hepatic autoimmunity via CRAMP, *Cell Reports*, 2021, V. 34(13).
- Shen, L., et al., Facile synthesis of organosilica-capped mesoporous silica nanocarriers with selective redox-triggered drug release properties for safe tumor chemotherapy, *Biomaterials science*, 2019.
- Kaminskas, L.M., et al., A 30 kDa polyethylene glycol-enfuvirtide complex enhances the exposure of enfuvirtide in lymphatic viral reservoirs in rats, *European Journal of Pharmaceutics and Biopharmaceutics*, 2019.
- Mahou, R., et al., Injectable and inherently vascularizing semi-interpenetrating polymer network for delivering cells to the subcutaneous space, *Biomaterials*, 2017.
- Wang, X.-M., et al., Exposure-time-dependent subcellular staging of gold nanoparticles deposition and vesicle destruction in mice livers, *Nanomedicine: Nanotechnology, Biology and Medicine*, 2021, V. 34.

7. Han, L., et al., Surface-enhanced resonance Raman scattering guided brain tumor surgery showing prognostic benefit in rat models, *ACS applied materials & interfaces*, 2019.
8. Liu, Y., et al., Cathodic protected Mn²⁺ by NaxWO₃ nanorods for stable magnetic resonance imaging-guided tumor photothermal therapy, *Biomaterials*, 2020, V. 234.
9. Guo, T., et al., Highly-selective detection of EGFR mutation gene in lung cancer based on surface enhanced Raman spectroscopy and asymmetric PCR, *Journal of Pharmaceutical and Biomedical Analysis*, 2020, 190,113522.
10. Mei, L., et al., Increased tumor targeted delivery using a multistage liposome system functionalized with RGD, TAT and cleavable PEG. *International Journal of Pharmaceutics*, 2014, 468(1–2): p. 26–38.
11. Day, J.R., et al., The impact of functional groups of poly(ethylene glycol) macromers on the physical properties of photo-polymerized hydrogels and the local inflammatory response in the host, *Acta Biomaterialia*, 2018, V. 67, P. 42-52.
12. Bae, S.W., et al., Dynamic, Bioresponsive Hydrogels via Changes in DNA Aptamer Conformation, *Macromolecular bioscience*, 2019, 19(2):1800353.
13. Li, H., et al., Bioengineered three-dimensional scaffolds to elucidate the effects of material biodegradability on cell behavior using POSS-PEG hybrid hydrogels, *Polymer Degradation and Stability*, 2019, 164:118-26.
14. Day, J.R., et al., The impact of functional groups of poly(ethylene glycol) macromers on the physical properties of photo-polymerized hydrogels and the local inflammatory response in the host, *Acta Biomaterialia*, 2018, 67, P. 42-52.
15. DeGroot, A.C., et al., Entropic control of receptor recycling using engineered ligands, *Biophysical journal*, 2018, 114(6), p.1377-1388.
16. Pereira Gomes, C., et al., Fine tuning neuronal targeting of nanoparticles by adjusting the ligand grafting density and combining PEG spacers of different length, *Acta Biomaterialia*, 2018, V. 78, P. 247-259.
17. Wan, X., et al., Effect of protein immunogenicity and PEG size and branching on the anti-PEG immune response to PEGylated proteins, *Process Biochemistry*, 2017, 52:183-91.
18. Zhao, Y.Z., et al., PEGylation with the thiosuccinimido butylamine linker significantly increases the stability of haloalkane dehalogenase DhaA, *Journal of Biotechnology*, 2017.
19. Xu, L., et al., Oxidized catechol-derived poly (ethylene glycol) for thiol-specific conjugation. *Reactive and Functional Polymers*, 2017.
20. Zhang, C., et al., Development of long-acting ciliary neurotrophic factor by site-specific conjugation with different-sized polyethylene glycols and transferrin, *International Journal of Pharmaceutics*, 2017.
21. Rong, L., et al., Iron chelated melanin-like nanoparticles for tumor-associated macrophage repolarization and cancer therapy, *Biomaterials*, 2019.
22. Adibnia, V., et al., Chitosan hydrogel micro-bio-devices with complex capillary patterns via reactive-diffusive self-assembly, *Acta Biomaterialia*, 2019.
23. Yue, Q., et al., An EGFRvIII targeted dual-modal gold nanoprobe for imaging-guided brain tumor surgery, *Nanoscale*, 2017.
24. Kampert, T., et al., Phenotypically Screened Carbon Nanoparticles for Enhanced Combinatorial Therapy in Triple Negative Breast Cancer, *Cellular and Molecular Bioengineering*, 2017, 1-6.
25. Billingsley, M.M., et al., Antibody-nanoparticle conjugates to enhance the sensitivity of ELISA-based detection methods, *PloS one*, 2017, 12(5).
26. Liu, S., et al., Effects of the proportion of two different cross-linkers on the material and biological properties of enzymatically degradable PEG hydrogels, *Polymer Degradation and Stability*, 2020, V. 172.
27. Rutz, A. L., et al., Employing PEG crosslinkers to optimize cell viability in gel phase bioinks and tailor post printing mechanical properties, *Acta Biomaterialia*, 2019.
28. DiVito, KA, et al., Hydrodynamic Focusing-Enabled Blood Vessel Fabrication for in Vitro Modeling of Neural Surrogates. *Journal of Medical and Biological Engineering*. 2021, 1-4.
29. Nam, K, et al., Anisotropically Functionalized Aptamer-DNA Nanostructures for Enhanced Cell Proliferation and Target-Specific Adhesion in 3D Cell Cultures. *Biomacromolecules*. 2021.
30. Peng, H, et al., Polymer Cell Surface Coating Enhances Mesenchymal Stem Cell Retention and Cardiac Protection. *ACS Applied Bio Materials*. 2021, 4(2):1655-67.
31. Tan, JJ, et al., Investigating PEGDA and GelMA Microgel Models for Sustained 3D Heterotypic Dermal Papilla and Keratinocyte Co-Cultures. *International Journal of Molecular Sciences*. 2021, 22(4):2143.
32. Chae, JJ, et al., Drug-Free, Nonsurgical Reduction of Intraocular Pressure for Four Months after Suprachoroidal Injection of Hyaluronic Acid Hydrogel. *Advanced Science*. 2021; 8(2):2001908.


PRODUCTS FOR N-TERMINAL PEGYLATION

JenKem Technology provides high quality polyethylene glycol (PEG) derivatives for N-Terminal PEGylation with high purity and low polydispersity.

JenKem Technology's activated PEG aldehydes for N-Terminal PEGylation, are reactive towards the N-Terminal amines influenced by the lower pKa of the N-Terminal amine in proteins. JenKem Technology's proprietary Y-shape PEGs are more selective, due to their sterically bulky structure. Y-shape Propionaldehyde PEG is reactive towards N-Terminal amines in the presence of a reducing reagent at pH 5-8; less selective but more reactive compared with Y-AALD-40K. Y-shape Acetaldehyde PEG is also an N-Terminal amine reactive PEG in the presence of a reducing reagent at pH 5-8; less reactive but more selective compared to linear PEG aldehydes. Methoxy PEG Propionaldehyde reacts with the N-Terminal at pH 5-8 in the presence of a reducing reagent.

N-Terminal PEGylation products with molecular weights, branching, and functional groups not listed in our catalog may be available by custom synthesis. GMP grade bulk products may also be available by custom synthesis. Please inquire at tech@jenkemusa.com about pricing and availability of custom activated PEGs.


Y-SHAPE PEG ACETALDEHYDE

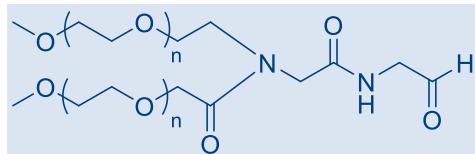

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A0003-1	Y-AALD-40K,1g	Y-shape PEG Acetaldehyde, MW 40000	\$170
*A0003-10	Y-AALD-40K,10g	Y-shape PEG Acetaldehyde, MW 40000	\$1,360

Y-SHAPE PEG PROPIONALDEHYDE

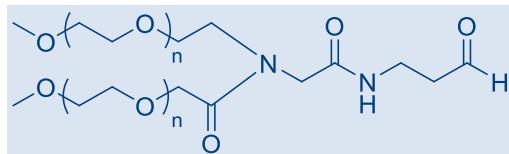

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A0017-1	Y-PALD-40K,1g	Y-shape PEG Propionaldehyde, MW 40000	\$170
A0017-10	Y-PALD-40K,10g	Y-shape PEG Propionaldehyde, MW 40000	\$1,360

METHOXY PEG PROPIONALDEHYDE

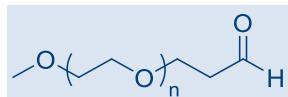

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3039-1	M-ALD-5000,1g	Methoxy PEG Propionaldehyde, MW 5000	\$120


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3039-10	M-ALD-5000,10g	Methoxy PEG Propionaldehyde, MW 5000	\$960
A3033-1	M-ALD-10K,1g	Methoxy PEG Propionaldehyde, MW 10000	\$120
A3033-10	M-ALD-10K,10g	Methoxy PEG Propionaldehyde, MW 10000	\$960
A3001-1	M-ALD-20K,1g	Methoxy PEG Propionaldehyde, MW 20000	\$120
A3001-10	M-ALD-20K,10g	Methoxy PEG Propionaldehyde, MW 20000	\$960
A3026-1	M-ALD-30K,1g	Methoxy PEG Propionaldehyde, MW 30000	\$130
A3026-10	M-ALD-30K,10g	Methoxy PEG Propionaldehyde, MW 30000	\$1,040
A3053-1	M-ALD-40K,1g	Methoxy PEG Propionaldehyde, MW 40000	\$150
A3053-10	M-ALD-40K,10g	Methoxy PEG Propionaldehyde, MW 40000	\$1,200

REFERENCES:

1. Yu, K.-M., et al., Preclinical evaluation of the mono-PEGylated recombinant human interleukin-11 in cynomolgus monkeys, *Toxicology and Applied Pharmacology*, 2018, V. 342, P. 39-49.
2. Lorey, S., et al., Novel Ubiquitin-derived High Affinity Binding Proteins with Tumor Targeting Properties, *J. Biol. Chem.*, 2014, 289, 8493-8507.
3. Liebner, R., et al., Head to Head Comparison of the Formulation and Stability of Concentrated Solutions of HESylated versus PEGylated Anakinra, *Journal of Pharmaceutical Sciences*, 2015, 104(2): 515-526.
4. Mayolo-Deloisa, K., et al., Aldehyde PEGylation of laccase from *Trametes versicolor* in route to increase its stability: effect on enzymatic activity, *Journal of Molecular Recognition*, 2015, 28(3): 173-179.
5. Zaykov, A. N., et al., Insulin-like peptide 5 fails to improve metabolism or body weight in obese mice, *Peptides*, 2019, 120.
6. Zaghami, A., et al., Mechanisms of activity loss for a multi-PEGylated protein by experiment and simulation, *Materials Today Chemistry*, 2019, 12:121-31.
7. Hernandez-Vargas, G., et al., Thermo-separating polymer-based aqueous two-phase systems for the recovery of PEGylated lysozyme species, *Journal of Chromatography B.*, 2019, 1105:120-8.
8. Behi, J., et al., Optimization of PEGylation reaction time and molar ratio of rhG-CSF toward increasing bioactive potency of monoPEGylated protein, *International Journal of Biological Macromolecules*, 2018, V. 109, P. 888-895.
9. Kateja, N., et al., Development of an integrated continuous PEGylation and purification Process for granulocyte colony stimulating factor, *Journal of Biotechnology*, 2020, 322, p. 79-89.
10. Zhang, Y., et al., Co-delivery of doxorubicin and curcumin by pH-sensitive prodrug nanoparticle for combination therapy of cancer. *Scientific Reports*, 2016, 6:21225.
11. Mata-Gomez, M.A., et al., Modelling of electrokinetic phenomena for capture of PEGylated ribonuclease A in a microdevice with insulating structures, *Biomicrofluidics*, 2016, 10(3): 033106.
12. Abbasi, S., et al., Design and cell cytotoxicity assessment of palmitoylated polyethylene glycol-grafted chitosan as nanomicelle carrier for paclitaxel. *J. Appl. Polym. Sci.*, 2015, 133, 43233.
13. Mayolo-Deloisa, K., et al., Aldehyde PEGylation of laccase from *Trametes versicolor* in route to increase its stability: effect on enzymatic activity, *Journal of Molecular Recognition*, 2015, 28(3): 173-179.
14. Wu, L., et al., N-Terminal Modification with Pseudo-Bifunctional PEG-Hexadecane Markedly Improves the Pharmacological Profile of Human Growth Hormone, *Molecular Pharmaceutics*, 2015.
15. Zaghami, A., et al., Determination of the degree of PEGylation of protein bioconjugates using data from proton nuclear magnetic resonance spectroscopy, *Data in Brief*, V. 25, 2019, 104037.
16. Zhao, Y.Z., et al., PEGylation with the thiosuccinimido butylamine linker significantly increases the stability of haloalkane dehalogenase DhaA, *Journal of Biotechnology*, 2017.
17. Mejia-Manzano, L.A., et al., Recovery of PEGylated and native lysozyme using an in situ aqueous two-phase system directly from the PEGylation reaction, *Journal of Chemical Technology and Biotechnology*, 2017.
18. Zhang, L., et al., Suppression for lung metastasis by depletion of collagen I and lysyl oxidase via losartan assisted with paclitaxel-loaded pH-sensitive liposomes in breast cancer, *Drug Delivery*, 2016.
19. Mayolo-Deloisa, K., et al., PEGylated protein separation using different hydrophobic interaction supports: Conventional and monolithic supports. *Biotechnology progress*, 2016.

OTHER Y-SHAPE AND LINEAR METHOXY PEG DERIVATIVES


Y-SHAPE PEG ALKYNE

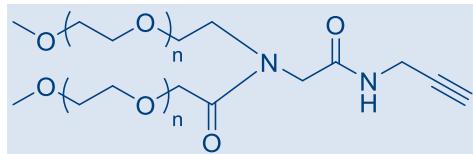

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A0026-1	Y-ALKYNE-40K,1g	Y-shape PEG Alkyne, MW 40000	\$160
*A0026-10	Y-ALKYNE-40K,10g	Y-shape PEG Alkyne, MW 40000	\$1,280


Y-SHAPE PEG AMINE

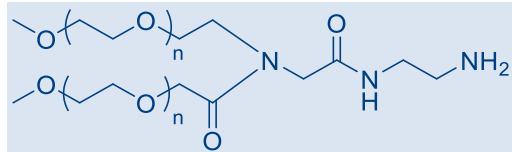

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A0010F-1	Y-NH2-40K,1g	Y-shape PEG Amine, MW 40000	\$160
*A0010F-10	Y-NH2-40K,10g	Y-shape PEG Amine, MW 40000	\$1,280


Y-SHAPE PEG FLUORESCIN

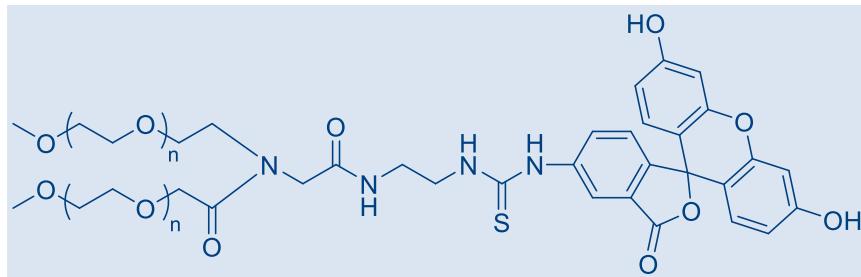

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A0014-1	Y-FITC-40K,1g	Y-shape PEG Fluorescein Isothiocyanate, MW 40000	\$250
*A0014-10	Y-FITC-40K,10g	Y-shape PEG Fluorescein Isothiocyanate, MW 40000	\$2,000

METHOXY PEG AMINE

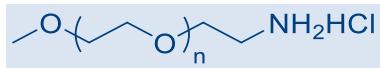

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3071-1	M-NH2HCl-2000,1g	Methoxy PEG Amine, HCl Salt, MW 2000	\$50
A3071-10	M-NH2HCl-2000,10g	Methoxy PEG Amine, HCl Salt, MW 2000	\$400
A3085-1	M-NH2HCl-3000,1g	Methoxy PEG Amine, HCl Salt, MW 3000	\$50
A3085-10	M-NH2HCl-3000,10g	Methoxy PEG Amine, HCl Salt, MW 3000	\$400
A3035-1	M-NH2HCl-5000,1g	Methoxy PEG Amine, HCl Salt, MW 5000	\$50
A3035-10	M-NH2HCl-5000,10g	Methoxy PEG Amine, HCl Salt, MW 5000	\$400
A3052-1	M-NH2HCl-10K,1g	Methoxy PEG Amine, HCl Salt, MW 10000	\$50
A3052-10	M-NH2HCl-10K,10g	Methoxy PEG Amine, HCl Salt, MW 10000	\$400
A3008-1	M-NH2HCl-20K,1g	Methoxy PEG Amine, HCl Salt, MW 20000	\$50
A3008-10	M-NH2HCl-20K,10g	Methoxy PEG Amine, HCl Salt, MW 20000	\$400
A3043-1	M-NH2HCl-30K,1g	Methoxy PEG Amine, HCl Salt, MW 30000	\$60
A3043-10	M-NH2HCl-30K,10g	Methoxy PEG Amine, HCl Salt, MW 30000	\$480
A3049-1	M-NH2HCl-40K,1g	Methoxy PEG Amine, HCl Salt, MW 40000	\$80
A3049-10	M-NH2HCl-40K,10g	Methoxy PEG Amine, HCl Salt, MW 40000	\$640


METHOXY PEG ACRYLATE

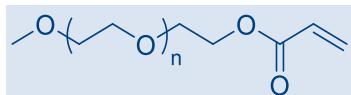

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A3098-1	M-ACLT-5000,1g	Methoxy PEG Acrylate, MW 5000	\$80
* A3098-10	M-ACLT-5000,10g	Methoxy PEG Acrylate, MW 5000	\$640
* A3099-1	M-ACLT-20K,1g	Methoxy PEG Acrylate, MW 20000	\$80
* A3099-10	M-ACLT-20K,10g	Methoxy PEG Acrylate, MW 20000	\$640

METHOXY PEG AZIDE

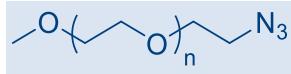


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3082-1	M-AZIDE-5000,1g	Methoxy PEG Azide, MW 5000	\$100
A3082-10	M-AZIDE-5000,10g	Methoxy PEG Azide, MW 5000	\$800
A3093-1	M-AZIDE-20K,1g	Methoxy PEG Azide, MW 20000	\$100
A3093-10	M-AZIDE-20K,10g	Methoxy PEG Azide, MW 20000	\$800

METHOXY PEG ALKyne

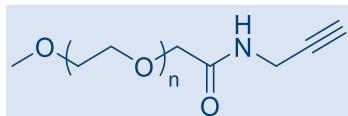

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3177-1	M-ALKYNE-5000,1g	Methoxy PEG Alkyne, MW 5000	\$140
A3177-10	M-ALKYNE-5000,10g	Methoxy PEG Alkyne, MW 5000	\$1,120


METHOXY PEG BIOTIN

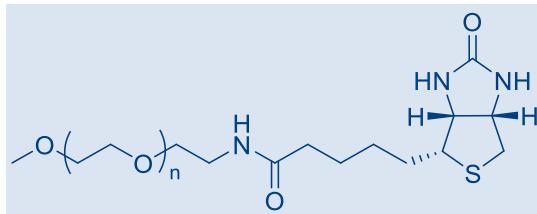

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3097-1	M-BIOTIN-5000,1g	Methoxy PEG Biotin, MW 5000	\$180
*A3097-10	M-BIOTIN-5000,10g	Methoxy PEG Biotin, MW 5000	\$1,440
*A3048-1	M-BIOTIN-20K,1g	Methoxy PEG Biotin, MW 20000	\$180
*A3048-10	M-BIOTIN-20K,10g	Methoxy PEG Biotin, MW 20000	\$1,440


METHOXY PEG HYDRAZIDE

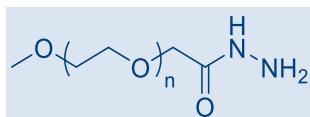

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3060-1	M-HZ-5000,1g	Methoxy PEG Hydrazide, MW 5000	\$120
*A3060-10	M-HZ-5000,10g	Methoxy PEG Hydrazide, MW 5000	\$960
*A3096-1	M-HZ-10K,1g	Methoxy PEG Hydrazide, MW 10000	\$120
*A3096-10	M-HZ-10K,10g	Methoxy PEG Hydrazide, MW 10000	\$960
*A3059-1	M-HZ-20K,1g	Methoxy PEG Hydrazide, MW 20000	\$120
*A3059-10	M-HZ-20K,10g	Methoxy PEG Hydrazide, MW 20000	\$960
*A3065-1	M-HZ-30K,1g	Methoxy PEG Hydrazide, MW 30000	\$140
*A3065-10	M-HZ-30K,10g	Methoxy PEG Hydrazide, MW 30000	\$1,120
*A3165F-1	M-HZ-40K,1g	Methoxy PEG Hydrazide, MW 40000	\$140
*A3165F-10	M-HZ-40K,10g	Methoxy PEG Hydrazide, MW 40000	\$1,120


METHOXY PEG DSPE

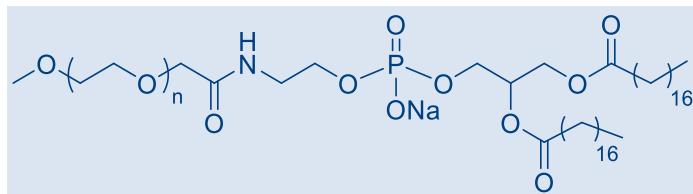

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3140-1	M-DSPE-5000,1g	Methoxy PEG DSPE, MW 5000	\$130
*A3140-10	M-DSPE-5000,10g	Methoxy PEG DSPE, MW 5000	\$1,040


METHOXY PEG DI-GLUTAMIC ACID

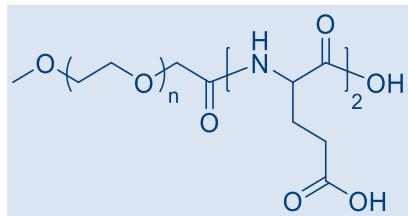

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3156-1	M-GLU2-20K,1g	Methoxy PEG di-Glutamic Acid, MW 20000	\$120
*A3156-10	M-GLU2-20K,10g	Methoxy PEG di-Glutamic Acid, MW 20000	\$960


METHOXY PEG TRI-GLUTAMIC ACID

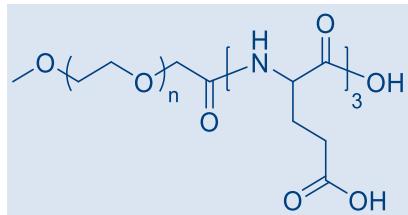

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3157-1	M-GLU3-20K,1g	Methoxy PEG tri-Glutamic Acid, MW 20000	\$150
*A3157-10	M-GLU3-20K,10g	Methoxy PEG tri-Glutamic Acid, MW 20000	\$1,200

METHOXY PEG PHOSPHATE

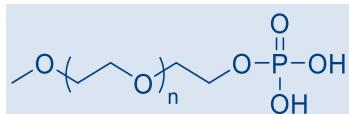

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3106-1	M-PHOS-5000,1g	Methoxy PEG Phosphate, MW 5000	\$120
A3106-10	M-PHOS-5000,10g	Methoxy PEG Phosphate, MW 5000	\$960

METHOXY PEG SILANE

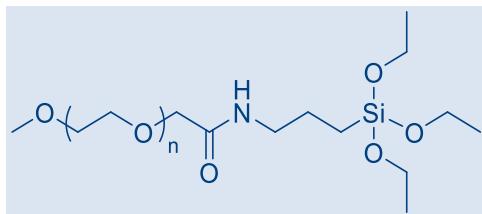


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A3037-1	M-SLN-5000,1g	Methoxy PEG Silane, MW 5000	\$130
A3037-10	M-SLN-5000,10g	Methoxy PEG Silane, MW 5000	\$1,040

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3100-1	M-SLN-20K,1g	Methoxy PEG Silane, MW 20000	\$130
*A3100-10	M-SLN-20K,10g	Methoxy PEG Silane, MW 20000	\$1,040

METHOXY PEG THIOCTIC ACID

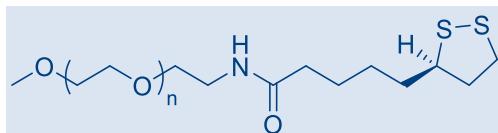

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A3076-1	M-TA-2000,1g	Methoxy PEG Thioctic Acid, MW 2000	\$200
*A3076-10	M-TA-2000,10g	Methoxy PEG Thioctic Acid, MW 2000	\$1,600

Please inquire if the PEG derivative that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.

REFERENCES:

1. Yu, K.-M., et al., Preclinical evaluation of the mono-PEGylated recombinant human interleukin-11 in cynomolgus monkeys, *Toxicology and Applied Pharmacology*, 2018, V. 342, P. 39-49.
2. Muddineti, O.S., et al., Vitamin-E/lipid based PEGylated polymeric micellar doxorubicin to sensitize doxorubicin-resistant cells towards treatment, *Reactive and Functional Polymers*, 2019, 134:49-57.
3. Yu, S., et al., A density gradient of VAPG peptides on a cell-resisting surface achieves selective adhesion and directional migration of smooth muscle cells over fibroblasts. *Acta biomaterialia*, 2018, 72, pp.70-81.
4. Muddineti, O.S., et al., Transferrin-Modified Vitamin-E/Lipid Based Polymeric Micelles for Improved Tumor Targeting and Anticancer Effect of Curcumin, *Pharmaceutical research*, 2018; 35(5):97.
5. Zhang, X., et al., Pharmaceutical micelles featured with singlet oxygen-responsive cargo release and mitochondrial targeting for enhanced photodynamic therapy, *Nanotechnology*, 2018, 29(25).
6. Gao, D., et al., Polymeric micelles encapsulating pH-responsive doxorubicin prodrug and glutathione-activated zinc (II) phthalocyanine for combined chemotherapy and photodynamic therapy, *Journal of Controlled Release*, 2018.
7. Shu, F., et al., Fabrication of a hyaluronic acid conjugated metal organic framework for targeted drug delivery and magnetic resonance imaging, *RSC Advances*, 2018, 8(12), pp.6581-6589.
8. Wang, L., et al., Delivery of radix ophiopogonis polysaccharide via sucrose acetateisobutyrate-based in situ forming systems alone or combined with itsmono-PEGylation, *Drug delivery*, 2018, 25(1), pp.267-277.
9. Wojnilowicz, M., et al., Glycogen-nucleic acid constructs for gene silencing in multicellular tumor spheroids, *Biomaterials*, 2018.
10. Khondee, S., et al., Doxorubicin-Loaded Micelle Targeting MUC1: A Potential Therapeutic for MUC1 Triple Negative Breast Cancer Treatment, *Current drug delivery*, 2018, 15(3), pp.406-416.
11. Sun, T., et al., A targeting theranostics nanomedicine as an alternative approach for hyperthermia perfusion, *Biomaterials*, 2018.
12. Pan, C., et al., Theranostic pH-sensitive nanoparticles for highly efficient targeted delivery of doxorubicin for breast tumor treatment. *International Journal of Nanomedicine*, 2018, 13.
13. Wang, X., et al., Polycation-telodendrimer nanocomplexes for intracellular protein delivery. *Colloids and Surfaces B: Biointerfaces*, 2018, 162, pp.405-414.
14. Wu, T., et al., Synthesis and characterization of poly (ethylene oxide)/polylactide/polylysine tri-arm star copolymers for gene delivery, *Journal of Polymer Science Part A: Polymer Chemistry*, 2018, 56(6).
15. Li, H., et al., Lactoferrin functionalized PEG-PLGA nanoparticles of shikonin for brain targeting therapy of glioma, *International journal of biological macromolecules* 2018, 107, pp.204-211.
16. Wang, X., et al., Polycation-telodendrimer nanocomplexes for intracellular protein delivery, *Colloids and Surfaces B: Biointerfaces* 2018, V. 162, p. 405-414.
17. DiVito, K.A., et al., Microfabricated blood vessels undergo neoangiogenesis, *Biomaterials*, 2017.
18. Hu, Y., et al., Facile Construction of Mitochondria-Targeting Nanoparticles for Enhanced Phototherapeutic Effects, *Biomaterials Science*, 2017.
19. Li, J., et al., Smart Asymmetric Vesicles with Triggered Availability of Inner Cell-Penetrating Shells for Specific Intracellular Drug Delivery, *ACS Applied Materials & Interfaces*, 2017.
20. Wu, Y., et al., Fabrication of thermo-sensitive complex micelles for reversible cell targeting, *Journal of Materials Science: Materials in Medicine*, 2015, 26:255.
21. Chen, Y., et al., An upconversion nanoparticle/Ru (ii) polypyridyl complex assembly for NIR-activated release of a DNA covalent-binding agent. *RSC Advances*. 2016, 6(28):23804-8.
22. Vila-Caballer, M., et al., A pH-sensitive stearoyl-PEG-poly(methacryloyl sulfadimethoxine)-decorated liposome system for protein delivery: An application for bladder cancer treatment, *Journal of Controlled Release* 2016.
23. Gonzalez-Villegas, J., et al., Poly (ethylene glycol)-modified zirconium phosphate nanoplatelets for improved doxorubicin delivery, *Inorganica Chimica Acta*, 2017.
24. Wang, L., et al., Iron-engineered mesoporous silica nanocatalyst with biodegradable and catalytic framework for tumor-specific therapy, *Biomaterials* 2018, V. 163, P. 1-13.
25. Yu, L., et al., Ultrasmall mesoporous organosilica nanoparticles: Morphology modulations and redox-responsive biodegradability for tumor-specific drug delivery, *Biomaterials* 2018, V. 161.


HOMOBIFUNCTIONAL PEG DERIVATIVES

JenKem Technology provides high quality activated homobifunctional polyethylene glycol derivatives (PEGs) with high purity and low polydispersity.

JenKem Technology homobifunctional PEG derivatives have numerous applications as cross-linkers, including for PEGylation of proteins and peptides, or nanoparticle and surface modifications. Homobifunctional PEG derivatives have the general structure: X—PEG—X, where X is a functional or reactive group.

Homobifunctional PEG products with molecular weights, branching, and functional groups not listed in our catalog may be available by custom synthesis. GMP grade bulk products may also be available by custom synthesis. Please inquire at tech@jenkemusa.com about pricing and availability of custom activated PEGs.

AMINE PEG AMINE

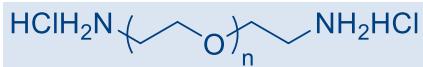


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A4015-1	NH2HCl-PEG2000-NH2HCl,1g	PEG (Amine)2, HCl Salt, MW 2000	\$60
A4015-5	NH2HCl-PEG2000-NH2HCl,5g	PEG (Amine)2, HCl Salt, MW 2000	\$240
A4016-1	NH2HCl-PEG3500-NH2HCl,1g	PEG (Amine)2, HCl Salt, MW 3500	\$60
A4016-5	NH2HCl-PEG3500-NH2HCl,5g	PEG (Amine)2, HCl Salt, MW 3500	\$240
A4017-1	NH2HCl-PEG5000-NH2HCl,1g	PEG (Amine)2, HCl Salt, MW 5000	\$60
A4017-5	NH2HCl-PEG5000-NH2HCl,5g	PEG (Amine)2, HCl Salt, MW 5000	\$240

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A4018-1	NH2HCl-PEG7500-NH2HCl,1g	PEG (Amine)2, HCl Salt, MW 7500	\$60
*A4018-5	NH2HCl-PEG7500-NH2HCl,5g	PEG (Amine)2, HCl Salt, MW 7500	\$240

ACRYLATE PEG ACRYLATE

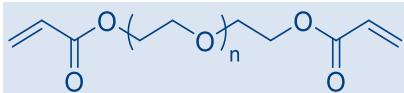


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A4047-1	ACLT-PEG2000-ACLT,1g	PEG (Acrylate)2, MW 2000	\$100
A4047-5	ACLT-PEG2000-ACLT,5g	PEG (Acrylate)2, MW 2000	\$400
A4048-1	ACLT-PEG3500-ACLT,1g	PEG (Acrylate)2, MW 3500	\$100
A4048-5	ACLT-PEG3500-ACLT,5g	PEG (Acrylate)2, MW 3500	\$400

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A4049-1	ACLT-PEG5000-ACLT,1g	PEG (Acrylate)2, MW 5000	\$100
*A4049-5	ACLT-PEG5000-ACLT,5g	PEG (Acrylate)2, MW 5000	\$400


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A4050-1	ACLT-PEG7500-ACLT,1g	PEG (Acrylate)2, MW 7500	\$100
*A4050-5	ACLT-PEG7500-ACLT,5g	PEG (Acrylate)2, MW 7500	\$400

ALKYNE PEG ALKYNE

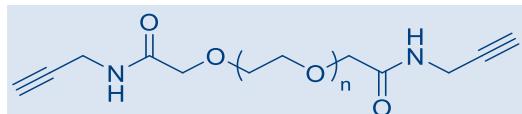

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A4095-1	ALKYNE-PEG5000-ALKYNE,1g	PEG (Alkyne)2, MW 5000	\$140
A4095-5	ALKYNE-PEG5000-ALKYNE,5g	PEG (Alkyne)2, MW 5000	\$560

ACETIC ACID PEG ACETIC ACID

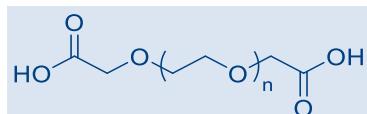


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A4003-1	CM-PEG2000-CM,1g	PEG (Acetic Acid)2, MW 2000	\$80
A4003-5	CM-PEG2000-CM,5g	PEG (Acetic Acid)2, MW 2000	\$320
A4004-1	CM-PEG3500-CM,1g	PEG (Acetic Acid)2, MW 3500	\$80
A4004-5	CM-PEG3500-CM,5g	PEG (Acetic Acid)2, MW 3500	\$320
A4005-1	CM-PEG5000-CM,1g	PEG (Acetic Acid)2, MW 5000	\$80
A4005-5	CM-PEG5000-CM,5g	PEG (Acetic Acid)2, MW 5000	\$320

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A4006-1	CM-PEG7500-CM,1g	PEG (Acetic Acid)2, MW 7500	\$80
*A4006-5	CM-PEG7500-CM,5g	PEG (Acetic Acid)2, MW 7500	\$320

MALEIMIDE PEG MALEIMIDE

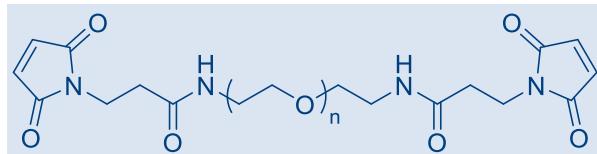


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A4011-1	MAL-PEG2000-MAL,1g	PEG (Maleimide)2, MW 2000	\$140
A4011-5	MAL-PEG2000-MAL,5g	PEG (Maleimide)2, MW 2000	\$560
A4010-1	MAL-PEG3500-MAL,1g	PEG (Maleimide)2, MW 3500	\$140
A4010-5	MAL-PEG3500-MAL,5g	PEG (Maleimide)2, MW 3500	\$560
A4009-1	MAL-PEG5000-MAL,1g	PEG (Maleimide)2, MW 5000	\$140

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A4009-5	MAL-PEG5000-MAL,5g	PEG (Maleimide)2, MW 5000	\$560

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A4014-1	MAL-PEG7500-MAL,1g	PEG (Maleimide)2, MW 7500	\$140
*A4014-5	MAL-PEG7500-MAL,5g	PEG (Maleimide)2, MW 7500	\$560

PEG (SUCCINIMIDYL CARBOXYMETHYL ESTER)2

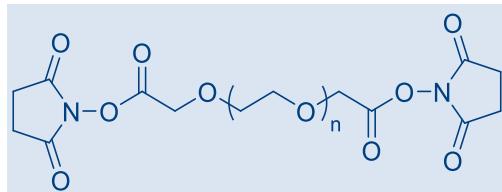


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A4051-1	SCM-PEG2000-SCM,1g	PEG (Succinimidyl Carboxymethyl Ester)2, MW 2000	\$100
A4051-5	SCM-PEG2000-SCM,5g	PEG (Succinimidyl Carboxymethyl Ester)2, MW 2000	\$400
A4052-1	SCM-PEG3500-SCM,1g	PEG (Succinimidyl Carboxymethyl Ester)2, MW 3500	\$100
A4052-5	SCM-PEG3500-SCM,5g	PEG (Succinimidyl Carboxymethyl Ester)2, MW 3500	\$400
A4013-1	SCM-PEG5000-SCM,1g	PEG (Succinimidyl Carboxymethyl Ester)2, MW 5000	\$100
A4013-5	SCM-PEG5000-SCM,5g	PEG (Succinimidyl Carboxymethyl Ester)2, MW 5000	\$400

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A4053-1	SCM-PEG7500-SCM,1g	PEG (Succinimidyl Carboxymethyl Ester)2, MW 7500	\$100
*A4053-5	SCM-PEG7500-SCM,5g	PEG (Succinimidyl Carboxymethyl Ester)2, MW 7500	\$400

THIOL PEG THIOL

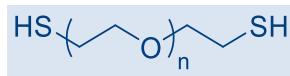

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A4001-1	HS-PEG3500-SH,1g	PEG (Thiol)2, MW 3500	\$120
A4001-5	HS-PEG3500-SH,5g	PEG (Thiol)2, MW 3500	\$480
A4075-1	HS-PEG5000-SH,1g	PEG (Thiol)2, MW 5000	\$120
A4075-5	HS-PEG5000-SH,5g	PEG (Thiol)2, MW 5000	\$480

VINYLSULFONE PEG VINYL SULFONE

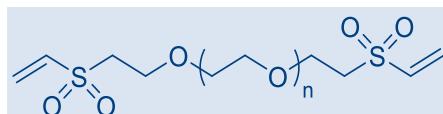

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A4076-1	VS-PEG3500-VS,1g	PEG (Vinylsulfone)2, MW 3500	\$150
*A4076-5	VS-PEG3500-VS,5g	PEG (Vinylsulfone)2, MW 3500	\$600
*A4077-1	VS-PEG5000-VS,1g	PEG (Vinylsulfone)2, MW 5000	\$150
*A4077-5	VS-PEG5000-VS,5g	PEG (Vinylsulfone)2, MW 5000	\$600

REFERENCES:

1. Wang, C., et al., Synthesis and formation mechanism of bone mineral, whitlockite nanocrystals in tri-solvent system, *Journal of Colloid and Interface Science*, 2020, V. 569, P. 1-11.
2. Ji, F., et al., A Dual pH/Magnetic Responsive Nanocarrier Based on PEGylated Fe3O4 Nanoparticles for Doxorubicin Delivery, *Journal of Nanoscience and Nanotechnology*, 2018, 18:7: 4464-4470.
3. Xiao, S., et al., Aptamer-mediated gene therapy enhanced antitumor activity against human hepatocellular carcinoma in vitro and in vivo, *Journal of Controlled Release*, 2017.
4. Jing, P., et al, Enhanced growth inhibition of prostate cancer in vitro and in vivo by a recombinant adenovirus-mediated dual-aptamer modified drug delivery system, *Cancer Letters*, 2016, V. 383(2).
5. Inostroza-Brito, K.E., et al., Cross-linking of a Biopolymer-Peptide Co-Assembling System, *Acta Biomaterialia*, 2017.
6. Lu, C., et al., Poly(ethylene glycol) crosslinked multi-armed poly(ϵ -benzyloxycarbonyl-L-lysine)s as super-amphiphiles: Synthesis, self-assembly, and evaluation as efficient delivery systems for poorly water-soluble drugs, *Colloids and Surfaces B: Biointerfaces*, 2019, 182.
7. Yang, S.-B., et al., The colossal role of H-MnO2-PEG in ischemic stroke, *Nanomedicine: Nanotechnology, Biology and Medicine*, 2021, 33, 102362.
8. Zhang, T., et al., Magnetic/pH dual-responsive nanocomposites loading doxorubicin hydrochloride for cancer therapy, *Micro & Nano Letters*, 2019.
9. Goh, S.C., et al., Polydopamine-polyethylene glycol-albumin antifouling coatings on multiple substrates: variations, *Journal of Materials Chemistry B*, 2018.
10. Kawanami, T., et al., A novel diclofenac-hydrogel conjugate system for intraarticular sustained release: Development of 2-pyridylamino-substituted 1-phenylethanol (PAPE) and its derivatives as tunable traceless linkers, *International Journal of Pharmaceutics*, 2020, 585, 119519.
11. Goor, O., et al., Introduction of anti-fouling coatings at the surface of supramolecular elastomeric materials via post-modification of reactive supramolecular additives, *Polymer Chemistry*, 2017.
12. Kim, H.C., et al., Highly stable and reduction responsive micelles from a novel polymeric surfactant with a repeating disulfide-based gemini structure for efficient drug delivery, *Polymer*, 2017.
13. Yu, H., et al., Polyvinylpyrrolidone functionalization induces deformable structure of graphene oxide nanosheets for lung-targeting delivery, *Nano Today*, 2021, V. 38.
14. Liao, H. T., et al, A bioactive multi-functional heparin-grafted aligned poly(lactide-co-glycolide)/curcumin nanofiber membrane to accelerate diabetic wound healing, *Materials Science and Engineering: C*, 2021.
15. Mou, J., et al., A New Green Titania with Enhanced NIR Absorption for Mitochondria-Targeted Cancer Therapy, *Theranostics*, 2017, 7(6):1531-1542.
16. Day, J.R., et al., The impact of functional groups of poly(ethylene glycol) macromers on the physical properties of photo-polymerized hydrogels and the local inflammatory response in the host, *Acta Biomaterialia*, 2018, V. 67, P. 42-52.
17. Ozdemir, T., et al., Bottom-up assembly of salivary gland microtissues for assessing myoepithelial cell function, *Biomaterials*, 2017, V. 142, P. 124-135.
18. Gottipati, A., et al., Gelatin Based Polymer Cell Coating Improves Bone Marrow-Derived Cell Retention in the Heart after Myocardial Infarction, *Stem Cell Reviews and Reports*, 2019.
19. Day, J.R., et al., The impact of functional groups of poly(ethylene glycol) macromers on the physical properties of photo-polymerized hydrogels and the local inflammatory response in the host, *Acta Biomaterialia*, 2018, V. 67, P. 42-52.
20. Tan, J.J., et al., Impact of substrate stiffness on dermal papilla aggregates in microgels, *Biomaterials science*, 2018.
21. Marozas, I. A., et al., Adaptable boronate ester hydrogels with tunable viscoelastic spectra to probe timescale dependent mechanotransduction, *Biomaterials*, 2019, 223.

22. Pedron, S., et al., Patterning Three-Dimensional Hydrogel Microenvironments Using Hyperbranched Polyglycerols for Independent Control of Mesh Size and Stiffness. *Biomacromolecules*, 2017, 18(4):1393-400.
23. DiVito, K.A., et al., Data characterizing microfabricated human blood vessels created via hydrodynamic focusing, *Data in Brief*, 2017, 14, P. 156-162.
24. Liu, S., et al., Effects of the proportion of two different cross-linkers on the material and biological properties of enzymatically degradable PEG hydrogels, *Polymer Degradation and Stability*, 2020, V. 172.
25. Truong, N.F., et al., Microporous annealed particle hydrogel stiffness, void space size, and adhesion properties impact cell proliferation, cell spreading, and gene transfer, *Acta Biomaterialia*, 2019.
26. Bae, S.W., et al., Dynamic, Bioresponsive Hydrogels via Changes in DNA Aptamer Conformation, *Macromolecular bioscience*, 2019, 19(2):1800353.
27. Stoichevska, V., et al., Engineering specific chemical modification sites into a collagen-like protein from *Streptococcus pyogenes*, *Journal of Biomedical Materials Research Part A*, 2017.
28. Li, H., et al., Bioengineered three-dimensional scaffolds to elucidate the effects of material biodegradability on cell behavior using POSS-PEG hybrid hydrogels, *Polymer Degradation and Stability*, 2019, 164:118-26.
29. Day, J.R., et al., The impact of functional groups of poly(ethylene glycol) macromers on the physical properties of photo-polymerized hydrogels and the local inflammatory response in the host, *Acta Biomaterialia*, 2018, 67, P. 42-52.
30. Rutz, A. L., et al., Employing PEG crosslinkers to optimize cell viability in gel phase bioinks and tailor post printing mechanical properties, *Acta Biomaterialia*, 2019.

HETEROBIFUNCTIONAL PEG DERIVATIVES


JenKem Technology provides high quality heterobifunctional glycol derivatives (PEGs) with high purity and low polydispersity. JenKem Technology's heterobifunctional PEG derivatives are generally employed as crosslinking agents or as spacers between two different chemical entities. The PEG moiety in the heterobifunctional PEG derivatives provides water solubility, biocompatibility, and flexibility to the linker. Heterobifunctional PEG derivatives have the following general structure:

X—PEG—Y, where X and Y are two different functional or reactive groups.

JenKem Technology's heterofunctional PEG Derivatives are useful as linkers for Antibody-Drug Conjugates or ADCs. ADCs conjugated via heterobifunctional PEGs exhibit improved water solubility and PK/PD profile.

Heterobifunctional PEG products with molecular weights, branching, and functional groups not listed in our catalog may be available by custom synthesis. GMP grade bulk products may also be available by custom synthesis. Please inquire at tech@jenkemusa.com about pricing and availability of custom activated PEGs.

ACRYLATE PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

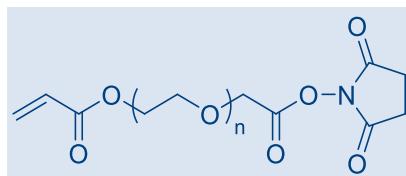


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5021-1	ACLT-PEG2000-SCM,1g	Acrylate PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$350
A5021-5	ACLT-PEG2000-SCM,5g	Acrylate PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$1,400
A5022-1	ACLT-PEG3500-SCM,1g	Acrylate PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$330
A5022-5	ACLT-PEG3500-SCM,5g	Acrylate PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$1,320
A5023-1	ACLT-PEG5000-SCM,1g	Acrylate PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$330

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5023-5	ACLT-PEG5000-SCM,5g	Acrylate PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$1,320

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A5024-1	ACLT-PEG7500-SCM,1g	Acrylate PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$330
* A5024-5	ACLT-PEG7500-SCM,5g	Acrylate PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$1,320


ACRYLATE PEG SUCCINIMIDYL PROPIONATE

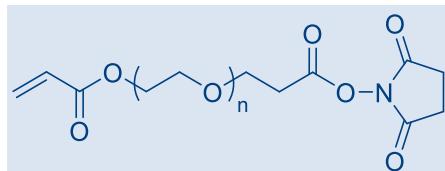

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A5320-1	ACLT-PEG3500-SPA,1g	Acrylate PEG Succinimidyl Propionate, MW 3500	\$396
* A5320-5	ACLT-PEG3500-SPA,5g	Acrylate PEG Succinimidyl Propionate, MW 3500	\$1,584
* A5321-1	ACLT-PEG5000-SPA,1g	Acrylate PEG Succinimidyl Propionate, MW 5000	\$396
* A5321-5	ACLT-PEG5000-SPA,5g	Acrylate PEG Succinimidyl Propionate, MW 5000	\$1,584

AMINE PEG ACETIC ACID

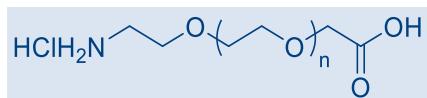


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5017-1	NH2HCl-PEG2000-CM,1g	Amine PEG Acetic Acid, HCl Salt, MW 2000	\$280
A5017-5	NH2HCl-PEG2000-CM,5g	Amine PEG Acetic Acid, HCl Salt, MW 2000	\$1,120
A5018-1	NH2HCl-PEG3500-CM,1g	Amine PEG Acetic Acid, HCl Salt, MW 3500	\$280
A5018-5	NH2HCl-PEG3500-CM,5g	Amine PEG Acetic Acid, HCl Salt, MW 3500	\$1,120
A5019-1	NH2HCl-PEG5000-CM,1g	Amine PEG Acetic Acid, HCl Salt, MW 5000	\$280
A5019-5	NH2HCl-PEG5000-CM,5g	Amine PEG Acetic Acid, HCl Salt, MW 5000	\$1,120

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A5236-1	NH2HCl-PEG1000-CM,1g	Amine PEG Acetic Acid, HCl Salt, MW 1000	\$350
* A5236-5	NH2HCl-PEG1000-CM,5g	Amine PEG Acetic Acid, HCl Salt, MW 1000	\$1,400
* A5020-1	NH2HCl-PEG7500-CM,1g	Amine PEG Acetic Acid, HCl Salt, MW 7500	\$280
* A5020-5	NH2HCl-PEG7500-CM,5g	Amine PEG Acetic Acid, HCl Salt, MW 7500	\$1,120
* A5159-1	NH2HCl-PEG10K-CM,1g	Amine PEG Acetic Acid, HCl Salt, MW 10000	\$350

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5159-5	NH2HCl-PEG10K-CM,5g	Amine PEG Acetic Acid, HCl Salt, MW 10000	\$1,400
*A5160-1	NH2HCl-PEG20K-CM,1g	Amine PEG Acetic Acid, HCl Salt, MW 20000	\$350
*A5160-5	NH2HCl-PEG20K-CM,5g	Amine PEG Acetic Acid, HCl Salt, MW 20000	\$1,400


ALKYNE PEG MALEIMIDE

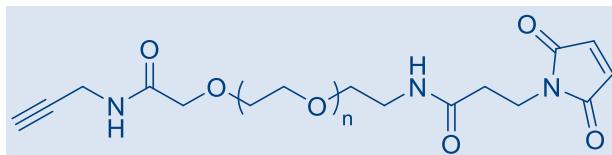

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5290-1	ALKYNE-PEG3500-MAL,1g	Alkyne PEG Maleimide, MW 3500	\$500
*A5290-5	ALKYNE-PEG3500-MAL,5g	Alkyne PEG Maleimide, MW 3500	\$2,000

AZIDE PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

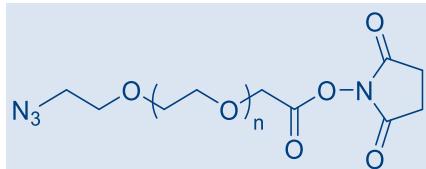

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5086-1	AZIDE-PEG3500-SCM,1g	Azide PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$500
A5086-5	AZIDE-PEG3500-SCM,5g	Azide PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$2,000
A5088-1	AZIDE-PEG5000-SCM,1g	Azide PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$500
A5088-5	AZIDE-PEG5000-SCM,5g	Azide PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$2,000


AZIDE PEG AMINE

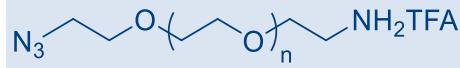

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5238-1	AZIDE-PEG3500-NH2TFA,1g	Azide PEG Amine, TFA Salt, MW 3500	\$400
A5238-5	AZIDE-PEG3500-NH2TFA,5g	Azide PEG Amine, TFA Salt, MW 3500	\$1,600
A5239-1	AZIDE-PEG5000-NH2TFA,1g	Azide PEG Amine, TFA Salt, MW 5000	\$400
A5239-5	AZIDE-PEG5000-NH2TFA,5g	Azide PEG Amine, TFA Salt, MW 5000	\$1,600

Please inquire if the PEG derivative that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.


BIOTIN PEG MALEIMIDE

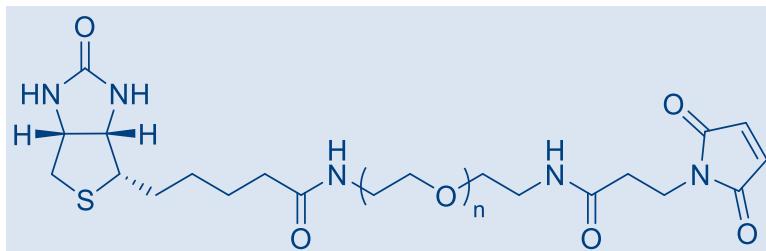

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5029-1	BIOTIN-PEG2000-MAL,1g	Biotin PEG Maleimide, MW 2000	\$550
*A5029-5	BIOTIN-PEG2000-MAL,5g	Biotin PEG Maleimide, MW 2000	\$2,200
*A5030-1	BIOTIN-PEG3500-MAL,1g	Biotin PEG Maleimide, MW 3500	\$550
*A5030-5	BIOTIN-PEG3500-MAL,5g	Biotin PEG Maleimide, MW 3500	\$2,200
*A5031-1	BIOTIN-PEG5000-MAL,1g	Biotin PEG Maleimide, MW 5000	\$550
*A5031-5	BIOTIN-PEG5000-MAL,5g	Biotin PEG Maleimide, MW 5000	\$2,200
*A5032-1	BIOTIN-PEG7500-MAL,1g	Biotin PEG Maleimide, MW 7500	\$550
*A5032-5	BIOTIN-PEG7500-MAL,5g	Biotin PEG Maleimide, MW 7500	\$2,200

BIOTIN PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

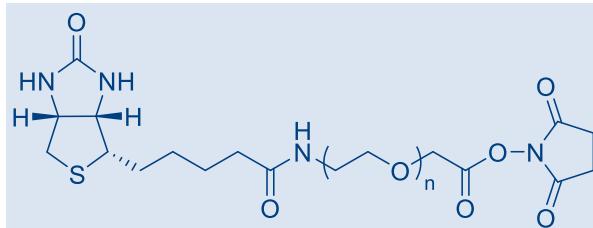


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5025-1	BIOTIN-PEG2000-SCM,1g	Biotin PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$500
A5025-5	BIOTIN-PEG2000-SCM,5g	Biotin PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$2,000

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5026-1	BIOTIN-PEG3500-SCM,1g	Biotin PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$500
*A5026-5	BIOTIN-PEG3500-SCM,5g	Biotin PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$2,000
*A5027-1	BIOTIN-PEG5000-SCM,1g	Biotin PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$500
*A5027-5	BIOTIN-PEG5000-SCM,5g	Biotin PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$2,000
*A5028-1	BIOTIN-PEG7500-SCM,1g	Biotin PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$500
*A5028-5	BIOTIN-PEG7500-SCM,5g	Biotin PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$2,000


BIOTIN PEG SUCCINIMIDYL GLUTARAMIDE

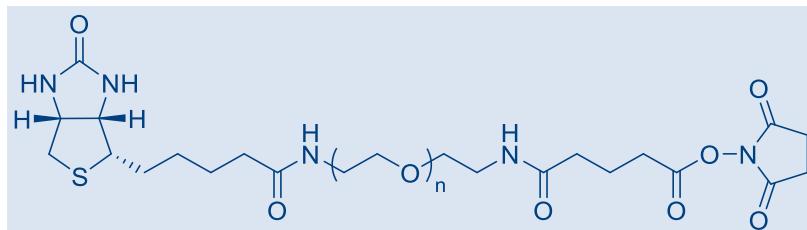

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5057-1	BIOTIN-PEG3500-SGA,1g	Biotin PEG Succinimidyl Glutaramide, MW 3500	\$550
*A5057-5	BIOTIN-PEG3500-SGA,5g	Biotin PEG Succinimidyl Glutaramide, MW 3500	\$2,200


FMOC AMINE PEG AMINE

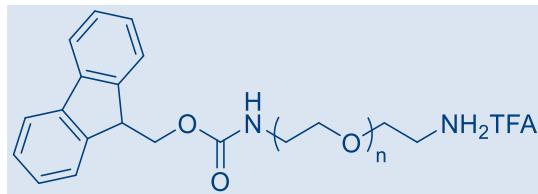

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5138-1	FMOC-PEG2000-NH2TFA,1g	Fmoc Amine PEG Amine, TFA Salt, MW 2000	\$300
*A5138-5	FMOC-PEG2000-NH2TFA,5g	Fmoc Amine PEG Amine, TFA Salt, MW 2000	\$1,200
*A5139-1	FMOC-PEG3500-NH2TFA,1g	Fmoc Amine PEG Amine, TFA Salt, MW 3500	\$300
*A5139-5	FMOC-PEG3500-NH2TFA,5g	Fmoc Amine PEG Amine, TFA Salt, MW 3500	\$1,200
*A5140-1	FMOC-PEG5000-NH2TFA,1g	Fmoc Amine PEG Amine, TFA Salt, MW 5000	\$300
*A5140-5	FMOC-PEG5000-NH2TFA,5g	Fmoc Amine PEG Amine, TFA Salt, MW 5000	\$1,200
*A5141-1	FMOC-PEG7500-NH2TFA,1g	Fmoc Amine PEG Amine, TFA Salt, MW 7500	\$300
*A5141-5	FMOC-PEG7500-NH2TFA,5g	Fmoc Amine PEG Amine, TFA Salt, MW 7500	\$1,200

FMOC AMINE PEG ACETIC ACID

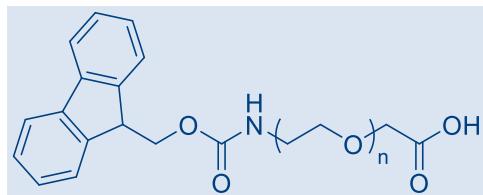


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5081-1	FMOC-PEG2000-CM,1g	Fmoc Amine PEG Acetic Acid, MW 2000	\$375
A5081-5	FMOC-PEG2000-CM,5g	Fmoc Amine PEG Acetic Acid, MW 2000	\$1,500

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5135-1	FMOC-PEG3500-CM,1g	Fmoc Amine PEG Acetic Acid, MW 3500	\$375
*A5135-5	FMOC-PEG3500-CM,5g	Fmoc Amine PEG Acetic Acid, MW 3500	\$1,500
*A5136-1	FMOC-PEG5000-CM,1g	Fmoc Amine PEG Acetic Acid, MW 5000	\$375
*A5136-5	FMOC-PEG5000-CM,5g	Fmoc Amine PEG Acetic Acid, MW 5000	\$1,500
*A5137-1	FMOC-PEG7500-CM,1g	Fmoc Amine PEG Acetic Acid, MW 7500	\$375
*A5137-5	FMOC-PEG7500-CM,5g	Fmoc Amine PEG Acetic Acid, MW 7500	\$1,500


FMOC AMINE PEG HYDROXYL

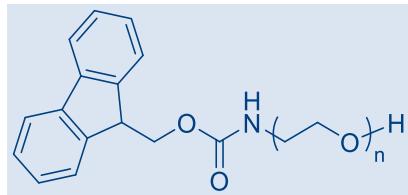

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5131-1	FMOC-PEG2000-OH,1g	Fmoc Amine PEG Hydroxyl, MW 2000	\$250
*A5131-5	FMOC-PEG2000-OH,5g	Fmoc Amine PEG Hydroxyl, MW 2000	\$1,000
*A5129-1	FMOC-PEG3500-OH,1g	Fmoc Amine PEG Hydroxyl, MW 3500	\$250
*A5129-5	FMOC-PEG3500-OH,5g	Fmoc Amine PEG Hydroxyl, MW 3500	\$1,000
*A5128-1	FMOC-PEG5000-OH,1g	Fmoc Amine PEG Hydroxyl, MW 5000	\$250
*A5128-5	FMOC-PEG5000-OH,5g	Fmoc Amine PEG Hydroxyl, MW 5000	\$1,000
*A5130-1	FMOC-PEG7500-OH,1g	Fmoc Amine PEG Hydroxyl, MW 7500	\$250
*A5130-5	FMOC-PEG7500-OH,5g	Fmoc Amine PEG Hydroxyl, MW 7500	\$1,000

FMOC AMINE PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

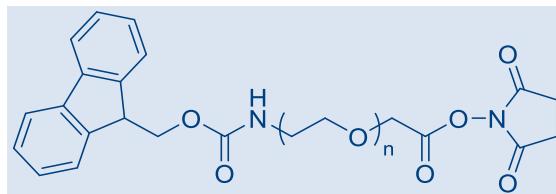


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5049-1	FMOC-PEG2000-SCM,1g	Fmoc Amine PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$400
A5049-5	FMOC-PEG2000-SCM,5g	Fmoc Amine PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$1,600

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5050-1	FMOC-PEG3500-SCM,1g	Fmoc Amine PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$400

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5050-5	FMOC-PEG3500-SCM,5g	Fmoc Amine PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$1,600
*A5051-1	FMOC-PEG5000-SCM,1g	Fmoc Amine PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$400
*A5051-5	FMOC-PEG5000-SCM,5g	Fmoc Amine PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$1,600
*A5052-1	FMOC-PEG7500-SCM,1g	Fmoc Amine PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$400
*A5052-5	FMOC-PEG7500-SCM,5g	Fmoc Amine PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$1,600

HYDROXYL PEG AMINE

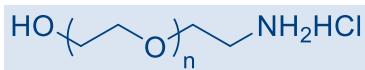


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5013-1	HO-PEG2000-NH2HCl,1g	Hydroxyl PEG Amine, HCl Salt, MW 2000	\$200
A5013-5	HO-PEG2000-NH2HCl,5g	Hydroxyl PEG Amine, HCl Salt, MW 2000	\$800
A5014-1	HO-PEG3500-NH2HCl,1g	Hydroxyl PEG Amine, HCl Salt, MW 3500	\$200
A5014-5	HO-PEG3500-NH2HCl,5g	Hydroxyl PEG Amine, HCl Salt, MW 3500	\$800
A5015-1	HO-PEG5000-NH2HCl,1g	Hydroxyl PEG Amine, HCl Salt, MW 5000	\$200
A5015-5	HO-PEG5000-NH2HCl,5g	Hydroxyl PEG Amine, HCl Salt, MW 5000	\$800

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5234-1	HO-PEG1000-NH2HCl,1g	Hydroxyl PEG Amine, HCl Salt, MW 1000	\$300
*A5234-5	HO-PEG1000-NH2HCl,5g	Hydroxyl PEG Amine, HCl Salt, MW 1000	\$1,200
*A5016-1	HO-PEG7500-NH2HCl,1g	Hydroxyl PEG Amine, HCl Salt, MW 7500	\$200
*A5016-5	HO-PEG7500-NH2HCl,5g	Hydroxyl PEG Amine, HCl Salt, MW 7500	\$800
*A5157-1	HO-PEG10K-NH2HCl,1g	Hydroxyl PEG Amine, HCl Salt, MW 10000	\$200
*A5157-5	HO-PEG10K-NH2HCl,5g	Hydroxyl PEG Amine, HCl Salt, MW 10000	\$800
*A5158-1	HO-PEG20K-NH2HCl,1g	Hydroxyl PEG Amine, HCl Salt, MW 20000	\$200
*A5158-5	HO-PEG20K-NH2HCl,5g	Hydroxyl PEG Amine, HCl Salt, MW 20000	\$800

HYDROXYL PEG ACETIC ACID

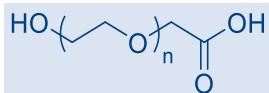


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5009-1	HO-PEG2000-CM,1g	Hydroxyl PEG Acetic Acid, MW 2000	\$200
A5009-5	HO-PEG2000-CM,5g	Hydroxyl PEG Acetic Acid, MW 2000	\$800
A5010-1	HO-PEG3500-CM,1g	Hydroxyl PEG Acetic Acid, MW 3500	\$200
A5010-5	HO-PEG3500-CM,5g	Hydroxyl PEG Acetic Acid, MW 3500	\$800
A5011-1	HO-PEG5000-CM,1g	Hydroxyl PEG Acetic Acid, MW 5000	\$200
A5011-5	HO-PEG5000-CM,5g	Hydroxyl PEG Acetic Acid, MW 5000	\$800

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5233-1	HO-PEG1000-CM,1g	Hydroxyl PEG Acetic Acid, MW 1000	\$300
*A5233-5	HO-PEG1000-CM,5g	Hydroxyl PEG Acetic Acid, MW 1000	\$1,200
*A5012-1	HO-PEG7500-CM,1g	Hydroxyl PEG Acetic Acid, MW 7500	\$200
*A5012-5	HO-PEG7500-CM,5g	Hydroxyl PEG Acetic Acid, MW 7500	\$800
*A5155-1	HO-PEG10K-CM,1g	Hydroxyl PEG Acetic Acid, MW 10000	\$250
*A5155-5	HO-PEG10K-CM,5g	Hydroxyl PEG Acetic Acid, MW 10000	\$1,000
*A5156-1	HO-PEG20K-CM,1g	Hydroxyl PEG Acetic Acid, MW 20000	\$300
*A5156-5	HO-PEG20K-CM,5g	Hydroxyl PEG Acetic Acid, MW 20000	\$1,200

HYDROXYL PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

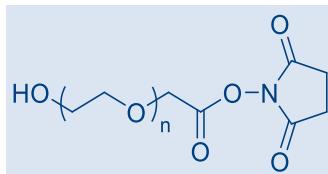


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5078-1	HO-PEG5000-SCM,1g	Hydroxyl PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$250
A5078-5	HO-PEG5000-SCM,5g	Hydroxyl PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$1,000

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5076-1	HO-PEG2000-SCM,1g	Hydroxyl PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$250
*A5076-5	HO-PEG2000-SCM,5g	Hydroxyl PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$1,000
*A5077-1	HO-PEG3500-SCM,1g	Hydroxyl PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$250
*A5077-5	HO-PEG3500-SCM,5g	Hydroxyl PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$1,000
*A5079-1	HO-PEG7500-SCM,1g	Hydroxyl PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$250
*A5079-5	HO-PEG7500-SCM,5g	Hydroxyl PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$1,000


HYDROXYL PEG PROPIONIC ACID

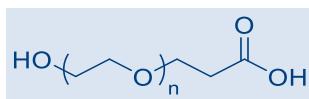

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5324-1	HO-PEG3500-PA,1g	Hydroxyl PEG Propionic Acid, MW 3500	\$300
*A5324-5	HO-PEG3500-PA,5g	Hydroxyl PEG Propionic Acid, MW 3500	\$1,200
*A5325-1	HO-PEG5000-PA,1g	Hydroxyl PEG Propionic Acid, MW 5000	\$300
*A5325-5	HO-PEG5000-PA,5g	Hydroxyl PEG Propionic Acid, MW 5000	\$1,200


HYDROXYL PEG SUCCINIMIDYL PROPIONATE

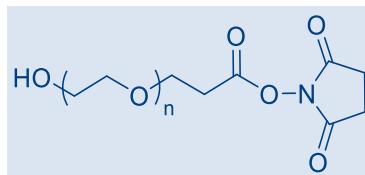

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5326-1	HO-PEG3500-SPA,1g	Hydroxyl PEG Succinimidyl Propionate, MW 3500	\$350
*A5326-5	HO-PEG3500-SPA,5g	Hydroxyl PEG Succinimidyl Propionate, MW 3500	\$1,400
*A5327-1	HO-PEG5000-SPA,1g	Hydroxyl PEG Succinimidyl Propionate, MW 5000	\$350
*A5327-5	HO-PEG5000-SPA,5g	Hydroxyl PEG Succinimidyl Propionate, MW 5000	\$1,400


HYDROXYL PEG HEXANOIC ACID

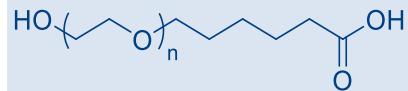

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5187-1	HO-PEG5000-HA,1g	Hydroxyl PEG Hexanoic Acid, MW 5000	\$300
*A5187-5	HO-PEG5000-HA,5g	Hydroxyl PEG Hexanoic Acid, MW 5000	\$1,200

THIOL PEG AMINE

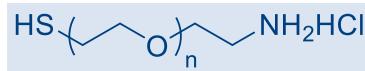


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5143-1	HS-PEG2000-NH2HCl,1g	Thiol PEG Amine, HCl Salt, MW 2000	\$400
A5143-5	HS-PEG2000-NH2HCl,5g	Thiol PEG Amine, HCl Salt, MW 2000	\$1,600
A5144-1	HS-PEG3500-NH2HCl,1g	Thiol PEG Amine, HCl Salt, MW 3500	\$350
A5144-5	HS-PEG3500-NH2HCl,5g	Thiol PEG Amine, HCl Salt, MW 3500	\$1,400
A5145-1	HS-PEG5000-NH2HCl,1g	Thiol PEG Amine, HCl Salt, MW 5000	\$350
A5145-5	HS-PEG5000-NH2HCl,5g	Thiol PEG Amine, HCl Salt, MW 5000	\$1,400

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5146-1	HS-PEG7500-NH2HCl,1g	Thiol PEG Amine, HCl Salt, MW 7500	\$350
*A5146-5	HS-PEG7500-NH2HCl,5g	Thiol PEG Amine, HCl Salt, MW 7500	\$1,400

THIOL PEG ACETIC ACID

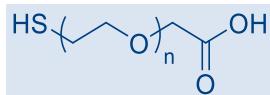


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5098-1	HS-PEG2000-CM,1g	Thiol PEG Acetic Acid, MW 2000	\$350
A5098-5	HS-PEG2000-CM,5g	Thiol PEG Acetic Acid, MW 2000	\$1,400
A5099-1	HS-PEG3500-CM,1g	Thiol PEG Acetic Acid, MW 3500	\$300
A5099-5	HS-PEG3500-CM,5g	Thiol PEG Acetic Acid, MW 3500	\$1,200
A5100-1	HS-PEG5000-CM,1g	Thiol PEG Acetic Acid, MW 5000	\$300
A5100-5	HS-PEG5000-CM,5g	Thiol PEG Acetic Acid, MW 5000	\$1,200

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5101-1	HS-PEG7500-CM,1g	Thiol PEG Acetic Acid, MW 7500	\$300
*A5101-5	HS-PEG7500-CM,5g	Thiol PEG Acetic Acid, MW 7500	\$1,200


THIOL PEG SUCCINIMIDYL GLUTARAMIDE

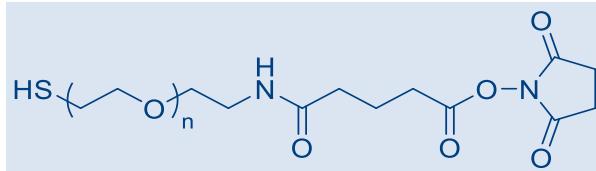

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5216-1	HS-PEG3500-SGA,1g	Thiol PEG Succinimidyl Glutaramide, MW 3500	\$600
*A5216-5	HS-PEG3500-SGA,5g	Thiol PEG Succinimidyl Glutaramide, MW 3500	\$2,400
*A5240-1	HS-PEG5000-SGA,1g	Thiol PEG Succinimidyl Glutaramide, MW 5000	\$600
*A5240-5	HS-PEG5000-SGA,5g	Thiol PEG Succinimidyl Glutaramide, MW 5000	\$2,400

THIOL PEG SUCCINIMIDYL PROPIONATE

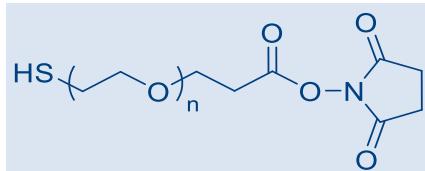


PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5330-1	HS-PEG3500-SPA,1g	Thiol PEG Succinimidyl Propionate, MW 3500	\$660

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5330-5	HS-PEG3500-SPA,5g	Thiol PEG Succinimidyl Propionate, MW 3500	\$2,640
*A5331-1	HS-PEG5000-SPA,1g	Thiol PEG Succinimidyl Propionate, MW 5000	\$720
*A5331-5	HS-PEG5000-SPA,5g	Thiol PEG Succinimidyl Propionate, MW 5000	\$2,880


T-BOC AMINE PEG HYDROXYL

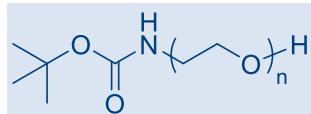

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5045-1	TBOC-PEG2000-OH,1g	t-Boc Amine PEG Hydroxyl, MW 2000	\$250
*A5045-5	TBOC-PEG2000-OH,5g	t-Boc Amine PEG Hydroxyl, MW 2000	\$1,000
*A5046-1	TBOC-PEG3500-OH,1g	t-Boc Amine PEG Hydroxyl, MW 3500	\$250
*A5046-5	TBOC-PEG3500-OH,5g	t-Boc Amine PEG Hydroxyl, MW 3500	\$1,000
*A5047-1	TBOC-PEG5000-OH,1g	t-Boc Amine PEG Hydroxyl, MW 5000	\$250
*A5047-5	TBOC-PEG5000-OH,5g	t-Boc Amine PEG Hydroxyl, MW 5000	\$1,000
*A5048-1	TBOC-PEG7500-OH,1g	t-Boc Amine PEG Hydroxyl, MW 7500	\$250
*A5048-5	TBOC-PEG7500-OH,5g	t-Boc Amine PEG Hydroxyl, MW 7500	\$1,000

T-BOC AMINE PEG AMINE

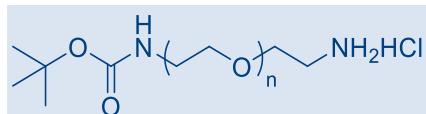


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5042-1	TBOC-PEG2000-NH2HCl,1g	t-Boc Amine PEG Amine, HCl Salt, MW 2000	\$280
A5042-5	TBOC-PEG2000-NH2HCl,5g	t-Boc Amine PEG Amine, HCl Salt, MW 2000	\$1,120

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5043-1	TBOC-PEG3500-NH2HCl,1g	t-Boc Amine PEG Amine, HCl Salt, MW 3500	\$280
*A5043-5	TBOC-PEG3500-NH2HCl,5g	t-Boc Amine PEG Amine, HCl Salt, MW 3500	\$1,120
*A5040-1	TBOC-PEG5000-NH2HCl,1g	t-Boc Amine PEG Amine, HCl Salt, MW 5000	\$280
*A5040-5	TBOC-PEG5000-NH2HCl,5g	t-Boc Amine PEG Amine, HCl Salt, MW 5000	\$1,120
*A5044-1	TBOC-PEG7500-NH2HCl,1g	t-Boc Amine PEG Amine, HCl Salt, MW 7500	\$280
*A5044-5	TBOC-PEG7500-NH2HCl,5g	t-Boc Amine PEG Amine, HCl Salt, MW 7500	\$1,120

T-BOC AMINE PEG ACETIC ACID

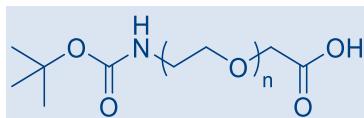


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5071-1	TBOC-PEG2000-CM,1g	t-Boc Amine PEG Acetic Acid, MW 2000	\$375
A5071-5	TBOC-PEG2000-CM,5g	t-Boc Amine PEG Acetic Acid, MW 2000	\$1,500

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5116-1	TBOC-PEG3500-CM,1g	t-Boc Amine PEG Acetic Acid, MW 3500	\$375
*A5116-5	TBOC-PEG3500-CM,5g	t-Boc Amine PEG Acetic Acid, MW 3500	\$1,500
*A5163-1	TBOC-PEG5000-CM,1g	t-Boc Amine PEG Acetic Acid, MW 5000	\$375
*A5163-5	TBOC-PEG5000-CM,5g	t-Boc Amine PEG Acetic Acid, MW 5000	\$1,500
*A5064-1	TBOC-PEG7500-CM,1g	t-Boc Amine PEG Acetic Acid, MW 7500	\$375
*A5064-5	TBOC-PEG7500-CM,5g	t-Boc Amine PEG Acetic Acid, MW 7500	\$1,500

T-BOC AMINE PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

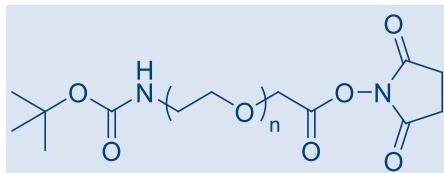


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5095-1	TBOC-PEG2000-SCM,1g	t-Boc Amine PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$400
A5095-5	TBOC-PEG2000-SCM,5g	t-Boc Amine PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$1,600

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5056-1	TBOC-PEG3500-SCM,1g	t-Boc Amine PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$400
*A5056-5	TBOC-PEG3500-SCM,5g	t-Boc Amine PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$1,600
*A5096-1	TBOC-PEG5000-SCM,1g	t-Boc Amine PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$400
*A5096-5	TBOC-PEG5000-SCM,5g	t-Boc Amine PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$1,600
*A5097-1	TBOC-PEG7500-SCM,1g	t-Boc Amine PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$400
*A5097-5	TBOC-PEG7500-SCM,5g	t-Boc Amine PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$1,600

Please inquire if the PEG derivative that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.

MALEIMIDE PEG AMINE

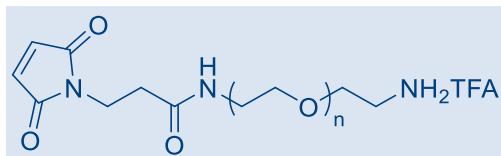


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5005-1	MAL-PEG2000-NH2TFA,1g	Maleimide PEG Amine, TFA Salt, MW 2000	\$350
A5005-5	MAL-PEG2000-NH2TFA,5g	Maleimide PEG Amine, TFA Salt, MW 2000	\$1,400
A5006-1	MAL-PEG3500-NH2TFA,1g	Maleimide PEG Amine, TFA Salt, MW 3500	\$350
A5006-5	MAL-PEG3500-NH2TFA,5g	Maleimide PEG Amine, TFA Salt, MW 3500	\$1,400
A5007-1	MAL-PEG5000-NH2TFA,1g	Maleimide PEG Amine, TFA Salt, MW 5000	\$350
A5007-5	MAL-PEG5000-NH2TFA,5g	Maleimide PEG Amine, TFA Salt, MW 5000	\$1,400

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5008-1	MAL-PEG7500-NH2TFA,1g	Maleimide PEG Amine, TFA Salt, MW 7500	\$350
*A5008-5	MAL-PEG7500-NH2TFA,5g	Maleimide PEG Amine, TFA Salt, MW 7500	\$1,400

MALEIMIDE PEG ACETIC ACID

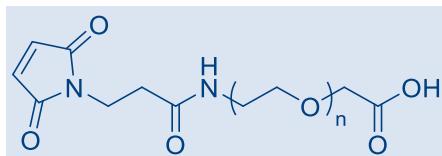


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5082-1	MAL-PEG2000-CM,1g	Maleimide PEG Acetic Acid, MW 2000	\$350
A5082-5	MAL-PEG2000-CM,5g	Maleimide PEG Acetic Acid, MW 2000	\$1,400
A5080-1	MAL-PEG3500-CM,1g	Maleimide PEG Acetic Acid, MW 3500	\$350
A5080-5	MAL-PEG3500-CM,5g	Maleimide PEG Acetic Acid, MW 3500	\$1,400
A5083-1	MAL-PEG5000-CM,1g	Maleimide PEG Acetic Acid, MW 5000	\$350
A5083-5	MAL-PEG5000-CM,5g	Maleimide PEG Acetic Acid, MW 5000	\$1,400

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5084-1	MAL-PEG7500-CM,1g	Maleimide PEG Acetic Acid, MW 7500	\$350
*A5084-5	MAL-PEG7500-CM,5g	Maleimide PEG Acetic Acid, MW 7500	\$1,400

Please inquire if the PEG derivative that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.

MALEIMIDE PEG HYDROXYL

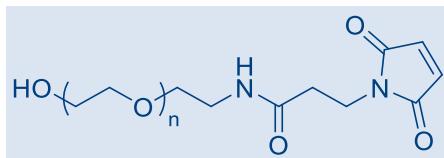


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5036-1	MAL-PEG5000-OH,1g	Maleimide PEG Hydroxyl, MW 5000	\$320
A5036-5	MAL-PEG5000-OH,5g	Maleimide PEG Hydroxyl, MW 5000	\$1,280

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5034-1	MAL-PEG2000-OH,1g	Maleimide PEG Hydroxyl, MW 2000	\$320
*A5034-5	MAL-PEG2000-OH,5g	Maleimide PEG Hydroxyl, MW 2000	\$1,280
*A5035-1	MAL-PEG3500-OH,1g	Maleimide PEG Hydroxyl, MW 3500	\$320
*A5035-5	MAL-PEG3500-OH,5g	Maleimide PEG Hydroxyl, MW 3500	\$1,280
*A5037-1	MAL-PEG7500-OH,1g	Maleimide PEG Hydroxyl, MW 7500	\$320
*A5037-5	MAL-PEG7500-OH,5g	Maleimide PEG Hydroxyl, MW 7500	\$1,280

MALEIMIDE PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

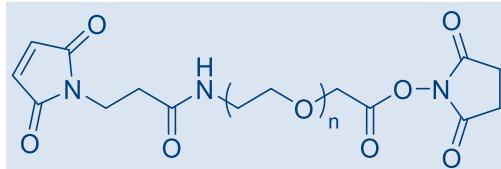

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5001-1	MAL-PEG2000-SCM,1g	Maleimide PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$450
A5001-5	MAL-PEG2000-SCM,5g	Maleimide PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$1,800
A5002-1	MAL-PEG3500-SCM,1g	Maleimide PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$425
A5002-5	MAL-PEG3500-SCM,5g	Maleimide PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$1,700
A5003-1	MAL-PEG5000-SCM,1g	Maleimide PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$425
A5003-5	MAL-PEG5000-SCM,5g	Maleimide PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$1,700

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5004-1	MAL-PEG7500-SCM,1g	Maleimide PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$425
*A5004-5	MAL-PEG7500-SCM,5g	Maleimide PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$1,700


OPSS PEG AMINE



PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5245-1	OPSS-PEG3500-NH2HCl,1g	OPSS PEG Amine, HCl Salt, MW 3500	\$600
*A5245-5	OPSS-PEG3500-NH2HCl,5g	OPSS PEG Amine, HCl Salt, MW 3500	\$2,400

OPSS PEG SUCCINIMIDYL CARBOXYMETHYL ESTER

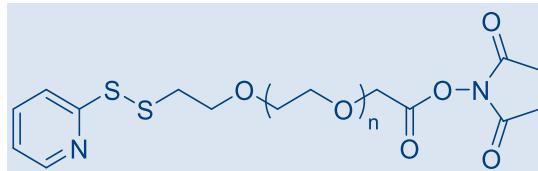

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A5107-1	OPSS-PEG2000-SCM,1g	OPSS PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$500
A5107-5	OPSS-PEG2000-SCM,5g	OPSS PEG Succinimidyl Carboxymethyl Ester, MW 2000	\$2,000

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A5108-1	OPSS-PEG3500-SCM,1g	OPSS PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$500
*A5108-5	OPSS-PEG3500-SCM,5g	OPSS PEG Succinimidyl Carboxymethyl Ester, MW 3500	\$2,000
*A5109-1	OPSS-PEG5000-SCM,1g	OPSS PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$500
*A5109-5	OPSS-PEG5000-SCM,5g	OPSS PEG Succinimidyl Carboxymethyl Ester, MW 5000	\$2,000
*A5110-1	OPSS-PEG7500-SCM,1g	OPSS PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$500
*A5110-5	OPSS-PEG7500-SCM,5g	OPSS PEG Succinimidyl Carboxymethyl Ester, MW 7500	\$2,000

REFERENCES:

- Wang, C., et al., Synthesis and formation mechanism of bone mineral, whitlockite nanocrystals in tri-solvent system, Journal of Colloid and Interface Science, 2020, V. 569, P. 1-11.
- Abstiens, K., et al., Gold-tagged Polymeric Nanoparticles with Spatially Controlled Composition for Enhanced Detectability in Biological Environments, ACS Applied Nano Materials, 2019.
- Abstiens, K., et al., Interaction of functionalized nanoparticles with serum proteins and its impact on colloidal stability and cargo leaching, Soft matter., 2019.
- Dutta, R., et al., Pharmacokinetics and Biodistribution of GDC-0449 Loaded Micelles in Normal and Liver Fibrotic Mice, Pharmaceutical research, 2017, 34(3):564-78.
- Wang, Y., et al., A pH-responsive silica–metal–organic framework hybrid nanoparticle for the delivery of hydrophilic drugs, nucleic acids, and CRISPR-Cas9 genome-editing machineries, Journal of Controlled Release, 2020.
- Tapia-Arellano, A., et al., Functionalization with PEG/Angiopep-2 peptide to improve the delivery of gold nanoparticles to central nervous system: in vitro and in vivo studies, Materials Science and Engineering: C, 2021.
- Zhao, B., et al., Single-step, wash-free digital immunoassay for rapid quantitative analysis of serological antibody against SARS-CoV-2 by photonic resonator absorption microscopy, Talanta, 2021, 225, 122004.
- Barros, D., et al., An affinity-based approach to engineer laminin-presenting cell instructive microenvironments, Biomaterials, 2019, 192:601-11.

9. Peters, M., et al., PEGylating poly(p-phenylene vinylene)-based bioimaging nanoprobes, *Journal of Colloid and Interface Science*, 2021, 581B, P. 566-575.
10. Feng, S., et al., Sorafenib encapsulated in nanocarrier functionalized with glycan-3 specific peptide for targeted therapy of hepatocellular carcinoma, *Colloids and Surfaces B: Biointerfaces*, 2019, 184.
11. Lumen, D., et al., Investigation of silicon nanoparticles produced by centrifuge chemical vapor deposition for applications in therapy and diagnostics, *European Journal of Pharmaceutics and Biopharmaceutics*, 2021, V. 158.
12. Poojari, R. et al., Distinct stratification of normal liver, hepatocellular carcinoma (HCC), and anticancer nanomedicine-treated tumor tissues by Raman fingerprinting for HCC therapeutic monitoring, *Nanomedicine: Nanotechnology, Biology and Medicine*, 2021, V. 33.
13. Feldmann, D.P., et al., The impact of microfluidic mixing of triblock micelleplexes on in vitro in vivo gene silencing and intracellular trafficking, *Nanotechnology*, 2017, 28(22):224001.
14. Han, Y., et al., Effective oral delivery of Exenatide-Zn²⁺ complex through distal ileum-targeted double layers nanocarriers modified with deoxycholic acid and glycocholic acid in diabetes therapy, *Biomaterials*, 2021, V. 275.
15. Soni, K.S., et al., Tuning polypeptide-based micellar carrier for efficient combination therapy of ErbB2-positive breast cancer, *Journal of Controlled Release*, 2017, V. 264, P. 276-287.
16. Chen, F., et al., Glycyrhetic acid-decorated and reduction-sensitive micelles to enhance the bioavailability and anti-hepatocellular carcinoma efficacy of tanshinone IIA, *Biomater. Sci.*, 2016, 4, 167-182.
17. Guarnieri, D., et al., Tumor-activated prodrug (TAP)-conjugated nanoparticles with cleavable domains for safe doxorubicin delivery, *Biotechnology and Bioengineering*, 2015, 112(3): 601-611.
18. Chen, Q., et al., Biodegradable nanoparticles decorated with different carbohydrates for efficient macrophage-targeted gene therapy, *Journal of Controlled Release*, 2020; 323:179-90.
19. Post, A., et al., Elucidation of Endothelial Cell Hemostatic Regulation with Integrin-Targeting Hydrogels, *Annals of biomedical engineering*, 2019.
20. Ghuman, H., et al., ECM hydrogel improves the delivery of PEG microsphere-encapsulated neural stem cells and endothelial cells into tissue cavities caused by stroke, *Brain Research Bulletin*, 2021, V. 168, P. 120-137.
21. Post, A., et al., Introduction of sacrificial bonds to hydrogels to increase defect tolerance during suturing of multilayer vascular grafts, *Acta Biomaterialia*, 2018, V. 69, p. 313-322.
22. Stukel, J.M., et al., The interplay of peptide affinity and scaffold stiffness on neuronal differentiation of neural stem cells. *Biomedical Materials*, 2018, 13(2), p.024102.
23. Xu, X., et al., Efficient and targeted drug/siRNA co-delivery mediated by reversibly crosslinked polymersomes toward anti-inflammatory treatment of ulcerative colitis (UC), *Nano Research*, 2019, 1-9.
24. Poli, H., et al., Optimisation of alendronate conjugation to polyethylene glycol for functionalisation of biopolymers and nanoparticles, *European Polymer Journal*, 2021, V. 155.
25. Babaei, M., et al., Targeted rod-shaped mesoporous silica nanoparticles for the co-delivery of camptothecin and survivin shRNA in to colon adenocarcinoma in vitro and in vivo, *European Journal of Pharmaceutics and Biopharmaceutics*, 2020.
26. Steinbauer, P., et al., Biomimetic adhesion motifs based on RAFT polymers with phosphonate groups, *European Polymer Journal*, 2021, V. 143.
27. JNi, J., et al., PSMA-targeted nanoparticles for specific penetration of blood-brain tumor barrier and combined therapy of brain metastases, *Journal of Controlled Release*, 2021, V. 329, P. 934-947.
28. Keuth, J., et al., Reversion of arterial calcification by elastin-targeted DTPA-HSA nanoparticles, *European Journal of Pharmaceutics and Biopharmaceutics*, 2020, V. 150, P. 108-119.
29. Xu, Y., et al., Cell-Penetrating Peptide Enhanced Insulin Buccal Absorption, *International Journal of Pharmaceutics*, 2020, 119469.
30. Khan, N.U., et al., Escape from abluminal LRP1-mediated clearance for boosted nanoparticle brain delivery and brain metastasis treatment, *Acta Pharmaceutica Sinica B*, 2020.
31. Jiang, H., et al., Multimodal theranostics augmented by transmembrane polymer-sealed nano-enzymatic porous MoS₂ nanoflowers, *International Journal of Pharmaceutics*, 2020, 586, 119606.
32. Balzer, C. J., et al., Single-Turnover Activation of Arp2/3 Complex by Dip1 May Balance Nucleation of Linear versus Branched Actin Filaments, *Current Biology*, 2019.
33. Wen, M., et al., Performance of TMC-g-PEG-VAPG/miRNA-145 complexes in electrospun membranes for target-regulating vascular SMCs, *Colloids and Surfaces B: Biointerfaces*, 2019, 182.
34. Kim, K.L., et al., Systematic detection of m6A-modified transcripts at single-molecule and single-cell resolution, *Cell Reports Methods*, 2021, V.1 (5).
35. Sun, T., et al., Co-delivery of Cu(I) chelator and chemotherapeutics as a new strategy for tumor theranostic, *Journal of Controlled Release*, 2020, V. 321, P. 483-496.

36. Badkas, A., et al., Modulation of in vitro phagocytic uptake and immunogenicity potential of modified Herceptin®-conjugated PLGA-PEG nanoparticles for drug delivery, *Colloids and Surfaces B: Biointerfaces*, 2018, V. 162.
37. Ju, L., et al., Von Willebrand factor-A1 domain binds platelet glycoprotein Iba in multiple states with distinctive force-dependent dissociation kinetics, *Thrombosis Research*, 2015, V. 136:3, P. 606-612.
38. Alpsoy, L., et al., Synthesis and Characterization of Carboxylated Luteolin (CL)-Functionalized SPION, *Journal of Superconductivity and Novel Magnetism*, 2017.
39. Li, Z., et al., Nanoparticle depots for controlled and sustained gene delivery, *Journal of Controlled Release*, 2020.
40. Lu, J., et al., Fabrication of thermo-and pH-sensitive cellulose nanofibrils-reinforced hydrogel with biomass nanoparticles, *Carbohydrate Polymers*, 2019.
41. Li, H., et al., Combination of active targeting, enzyme-triggered release and fluorescent dye into gold nanoclusters for endomicroscopy-guided photothermal/photodynamic therapy to pancreatic ductal adenocarcinoma, *Biomaterials*, 2017.
42. Stefanick, J.F., et al., Dual-receptor targeted strategy in nanoparticle design achieves tumor cell selectivity through cooperativity, *Nanoscale*, 2019, 11(10):4414-27.
43. Lim, YGet al., Calcium-binding near-infrared fluorescent nanoprobe for bone tissue imaging. *Journal of Industrial and Engineering Chemistry*, 2020, 89:442-7.
44. Lv, F., et al., Enhanced mucosal penetration and efficient inhibition efficacy against cervical cancer of PEGylated docetaxel nanocrystals by TAT modification, *Journal of Controlled Release*, 2021, V. 336, P. 572-582.
45. Chen, Q., et al., Deep tumor-penetrated nanosystem eliminates cancer stem cell for highly efficient liver cancer therapy, *Chemical Engineering Journal*, 2021, V. 421 (2).
46. Xia, C., et al., Redox-responsive nanoassembly restrained myeloid-derived suppressor cells recruitment through autophagy-involved lactate dehydrogenase A silencing for enhanced cancer immunochemotherapy, *Journal of Controlled Release*, 2021, V. 335, P. 557-574.


MONODISPERSE PEG DERIVATIVES

JenKem Technology's monodisperse PEG products are synthesized via very reproducible chemical reactions. JenKem Technology's monodisperse PEGs retain the polyethylene glycol's good solubility in water, offering to PEGylated drugs same PEG benefits of reduced aggregation and immunogenicity, and increased hydrodynamic volume.

JenKem Technology catalog monodisperse PEG products consist of linear chains of five to twelve –CH₂-CH₂-O- units, and various terminal functional and reactive groups. Monodisperse PEG products with molecular weights, functional groups and branching not listed in our catalog may be available by custom synthesis. Please inquire at tech@jenkemusa.com about pricing and availability.


METHOXY MONODISPERSE PEG DERIVATIVES

METHOXY-PEG(5 OR 6)-AMINE

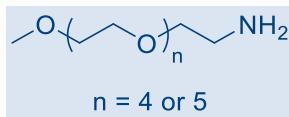


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A13020-1	M-PEG5-NH2,1g	Methoxy PEG5 Amine, MW 251.4	\$235
A13020-5	M-PEG5-NH2,5g	Methoxy PEG5 Amine, MW 251.4	\$500
A13001-1	M-PEG6-NH2,1g	Methoxy PEG6 Amine, MW 295.3	\$240
A13001-5	M-PEG6-NH2,5g	Methoxy PEG6 Amine, MW 295.3	\$525

METHOXY-PEG(5 OR 6)-ACRYLATE

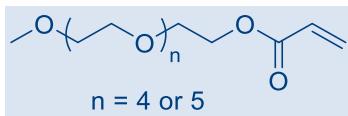


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A13005-1	M-PEG5-ACLT,1g	Methoxy PEG5 Acrylate, MW 306.3	\$400
A13005-5	M-PEG5-ACLT,5g	Methoxy PEG5 Acrylate, MW 306.3	\$1,600
A13007-1	M-PEG6-ACLT,1g	Methoxy PEG6 Acrylate, MW 350.4	\$400
A13007-5	M-PEG6-ACLT,5g	Methoxy PEG6 Acrylate, MW 350.4	\$1,600

METHOXY-PEG(5 OR 6)-AZIDE

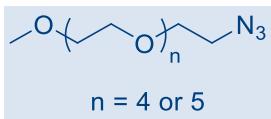


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A13016-1	M-PEG5-AZIDE,1g	Methoxy PEG5 Azide, MW 277.3	\$400
A13016-5	M-PEG5-AZIDE,5g	Methoxy PEG5 Azide, MW 277.3	\$1,600
A13017-1	M-PEG6-AZIDE,1g	Methoxy PEG6 Azide, MW 321.4	\$400
A13017-5	M-PEG6-AZIDE,5g	Methoxy PEG6 Azide, MW 321.4	\$1,600

METHOXY-PEG(5 OR 6)-HYDROXYL

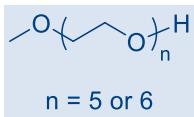

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
P12002-5	M-PEG5-OH,5g	Methoxy PEG5 Hydroxyl, MW 252.2	\$184
P12002-25	M-PEG5-OH,25g	Methoxy PEG5 Hydroxyl, MW 252.2	\$460
P12001-5	M-PEG6-OH,5g	Methoxy PEG6 Hydroxyl, MW 296.2	\$196
P12001-25	M-PEG6-OH,25g	Methoxy PEG6 Hydroxyl, MW 296.2	\$490


METHOXY-PEG(5 OR 6)-PROPIONIC ACID

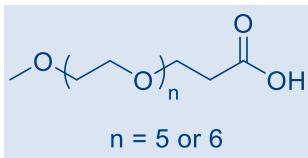

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A13010-1	M-PEG5-PA,1g	Methoxy PEG5 Propionic Acid, MW 324.4	\$220
A13010-5	M-PEG5-PA,5g	Methoxy PEG5 Propionic Acid, MW 324.4	\$460
A13011-1	M-PEG6-PA,1g	Methoxy PEG6 Propionic Acid, MW 368.2	\$235

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A13011-5	M-PEG6-PA,5g	Methoxy PEG6 Propionic Acid, MW 368.2	\$500


METHOXY-PEG(5 OR 6)-THIOL

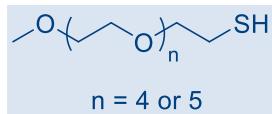

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A13012-1	M-PEG5-SH,1g	Methoxy PEG5 Thiol, MW 268.4	\$300
A13012-5	M-PEG5-SH,5g	Methoxy PEG5 Thiol, MW 268.4	\$1,000
A13013-1	M-PEG6-SH,1g	Methoxy PEG6 Thiol, MW 312.4	\$339
A13013-5	M-PEG6-SH,5g	Methoxy PEG6 Thiol, MW 312.4	\$1,130


HOMOBIFUNCTIONAL MONODISPERSE PEG MONODISPERSE DERIVATIVES

AMINE PEG(8 OR 12) AMINE

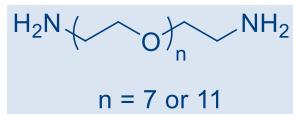


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A12011-1	NH2-PEG8-NH2,1g	Amine PEG8 Amine, MW 368.5	\$384
A12011-5	NH2-PEG8-NH2,5g	Amine PEG8 Amine, MW 368.5	\$1,280
A12012-1	NH2-PEG12-NH2,1g	Amine PEG12 Amine, MW 544.8	\$519
A12012-5	NH2-PEG12-NH2,5g	Amine PEG12 Amine, MW 544.8	\$1,730

ACETIC ACID PEG(8 OR 12) ACETIC ACID

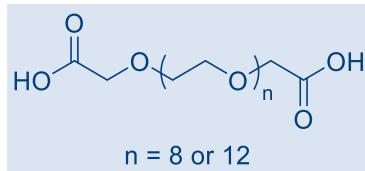


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A12019-1	CM-PEG8-CM,1g	Acetic Acid PEG8 Acetic Acid, MW 486.5	\$444
A12019-5	CM-PEG8-CM,5g	Acetic Acid PEG8 Acetic Acid, MW 486.5	\$1,480
A12020-1	CM-PEG12-CM,1g	Acetic Acid PEG12 Acetic Acid, MW 662.7	\$468
A12020-5	CM-PEG12-CM,5g	Acetic Acid PEG12 Acetic Acid, MW 662.7	\$1,560

HYDROXYL-PEG(8 OR 9 OR 12)-HYDROXYL

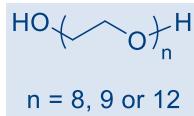

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
P11002-5	HO-PEG8-OH,5g	Octaethylene Glycol, MW 370.5	\$225
P11002-25	HO-PEG8-OH,25g	Octaethylene Glycol, MW 370.5	\$476
P11003-5	HO-PEG9-OH,5g	Enneethylene glycol, MW 414.5	\$235
P11003-25	HO-PEG9-OH,25g	Enneethylene glycol, MW 414.5	\$490
P11004-5	HO-PEG12-OH,5g	Dodecaethylene glycol, MW 546.7	\$250
P11004-25	HO-PEG12-OH,25g	Dodecaethylene glycol, MW 546.7	\$625


PROPIONIC ACID PEG(8 OR 12) PROPIONIC ACID

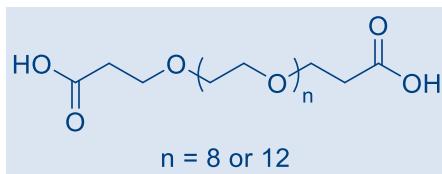

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A12026-1	PA-PEG8-PA,1g	Propionic Acid PEG8 Propionic Acid, MW 514.6	\$372
A12026-5	PA-PEG8-PA,5g	Propionic Acid PEG8 Propionic Acid, MW 514.6	\$1,240
A12018-1	PA-PEG12-PA,1g	Propionic Acid PEG12 Propionic Acid, MW 690.8	\$390
A12018-5	PA-PEG12-PA,5g	Propionic Acid PEG12 Propionic Acid, MW 690.8	\$1,300


HETEROBIFUNCTIONAL MONODISPERSE PEG MONODISPERSE DERIVATIVES

HYDROXYL PEG(8 OR 12) AMINE

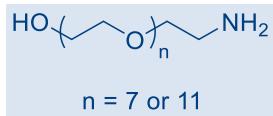


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A12008-1	HO-PEG8-NH2,1g	Hydroxyl PEG8 Amine, MW 369.5	\$324
A12008-5	HO-PEG8-NH2,5g	Hydroxyl PEG8 Amine, MW 369.5	\$1,080
A12009-1	HO-PEG12-NH2,1g	Hydroxyl PEG12 Amine, MW 545.8	\$519
A12009-5	HO-PEG12-NH2,5g	Hydroxyl PEG12 Amine, MW 545.8	\$1,730

Please inquire if the PEG derivative that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.

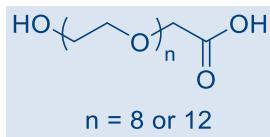HYDROXYL PEG(8 OR 12) ACETIC ACID

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A12014-1	HO-PEG8-CM,1g	Hydroxyl PEG8 Acetic Acid, MW 428.5	\$550
A12014-5	HO-PEG8-CM,5g	Hydroxyl PEG8 Acetic Acid, MW 428.5	\$2,200
A12015-1	HO-PEG12-CM,1g	Hydroxyl PEG12 Acetic Acid, MW 604.8	\$675
A12015-5	HO-PEG12-CM,5g	Hydroxyl PEG12 Acetic Acid, MW 604.8	\$2,700

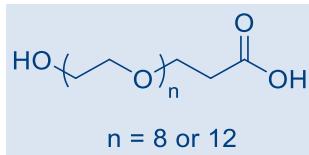HYDROXYL PEG(8 OR 12) PROPIONIC ACID

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A12002-1	HO-PEG8-PA,1g	Hydroxyl PEG8 Propionic Acid, MW 442.5	\$279
A12002-5	HO-PEG8-PA,5g	Hydroxyl PEG8 Propionic Acid, MW 442.5	\$930
A12003-1	HO-PEG12-PA,1g	Hydroxyl PEG12 Propionic Acid, MW 618.8	\$438
A12003-5	HO-PEG12-PA,5g	Hydroxyl PEG12 Propionic Acid, MW 618.8	\$1,460

MULTI-ARM PEG DERIVATIVES

JenKem Technology provides high quality activated multi-arm polyethylene glycol derivatives (PEGs) with high purity and low polydispersity.


JenKem Technology's multi-arm homo-and hetero-functional PEG derivatives can be cross-linked into hydrogels. PEG hydrogels have a variety of applications in medical devices and regenerative medicine, and are especially of interest for controlled release of drugs, for 3D cell culture, and for wound sealing and healing. JenKem Technology's multi-arm heterofunctional PEG Derivatives are designed as linkers for Antibody-Drug Conjugates or ADCs, as they facilitate the coupling of high ratios of drug molecules per antibody. ADCs conjugated via multiarm PEGs deliver a higher concentration of drugs to targets for each antibody-antigen interaction, and exhibit improved solubility and PK/PD profile.

JenKem Technology's multi-arm star PEGs are synthesized by ethoxylation of glycerol (3ARM PEG), pentaerythritol (4ARM PEG), tripentaerythritol (8ARM(TP) PEG), hexaglycerol (8ARM PEG), or dipentaerythritol (6ARM PEG). The number of ethylene oxide units in the PEG chain may not be equal for all arms. The total molecular weight reported for the JenKem multi-arm PEGs is the sum of the molecular weights of the PEG chains for each arm.

Multi-arm star PEG products with molecular weights and functional groups not listed in our catalog may be available by custom synthesis. Please inquire at tech@jenkemusa.com about pricing and availability.

3ARM-PEG DERIVATIVES (GLYCEROL CORE)


3ARM-PEG-AMINE

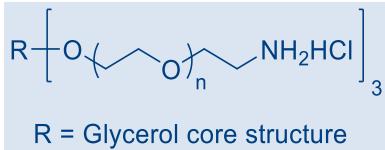

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A9027-1	3ARM-NH2HCl-15K,1g	3arm PEG Amine (glycerol), HCl salt, MW 15000	\$60
A9027-10	3ARM-NH2HCl-15K,10g	3arm PEG Amine (glycerol), HCl salt, MW 15000	\$480
A9028-1	3ARM-NH2HCl-20K,1g	3arm PEG Amine (glycerol), HCl salt, MW 20000	\$60
A9028-10	3ARM-NH2HCl-20K,10g	3arm PEG Amine (glycerol), HCl salt, MW 20000	\$480


4ARM-PEG DERIVATIVES (PENTAERYTHRITOL CORE)

4ARM PEG AMINE (PENTAERYTHRITOL)

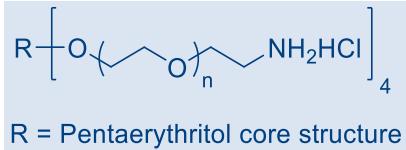

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7032F-1	4ARM-NH2-2000,1g	4arm PEG Amine (pentaerythritol), Free Amine, MW 2000	\$70
A7032F-10	4ARM-NH2-2000,10g	4arm PEG Amine (pentaerythritol), Free Amine, MW 2000	\$560
A7020-1	4ARM-NH2HCl-5000,1g	4arm PEG Amine (pentaerythritol), HCl Salt, MW 5000	\$50
A7020-10	4ARM-NH2HCl-5000,10g	4arm PEG Amine (pentaerythritol), HCl Salt, MW 5000	\$400
A7011-1	4ARM-NH2HCl-10K,1g	4arm PEG Amine (pentaerythritol), HCl Salt, MW 10000	\$50
A7011-10	4ARM-NH2HCl-10K,10g	4arm PEG Amine (pentaerythritol), HCl Salt, MW 10000	\$400
A7026-1	4ARM-NH2HCl-20K,1g	4arm PEG Amine (pentaerythritol), HCl Salt, MW 20000	\$50
A7026-10	4ARM-NH2HCl-20K,10g	4arm PEG Amine (pentaerythritol), HCl Salt, MW 20000	\$400
A7024-1	4ARM-NH2HCl-40K,1g	4arm PEG Amine (pentaerythritol), HCl Salt, MW 40000	\$60
A7024-10	4ARM-NH2HCl-40K,10g	4arm PEG Amine (pentaerythritol), HCl Salt, MW 40000	\$480


4ARM PEG ACRYLATE (PENTAERYTHRITOL)

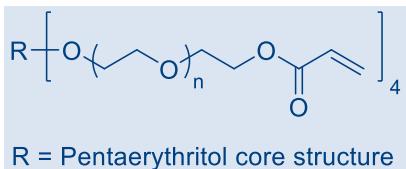

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7068-1	4ARM-ACLT-10K,1g	4arm PEG Acrylate (pentaerythritol), MW 10000	\$80
A7068-10	4ARM-ACLT-10K,10g	4arm PEG Acrylate (pentaerythritol), MW 10000	\$640

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7034-1	4ARM-ACLT-20K,1g	4arm PEG Acrylate (pentaerythritol), MW 20000	\$80
A7034-10	4ARM-ACLT-20K,10g	4arm PEG Acrylate (pentaerythritol), MW 20000	\$640


4ARM PEG AZIDE (PENTAERYTHRITOL)

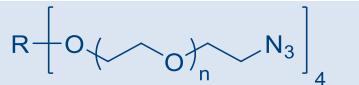

R = Pentaerythritol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7163-1	4ARM-AZIDE-10K,1g	4arm PEG Azide (pentaerythritol), MW 10000	\$130
A7163-10	4ARM-AZIDE-10K,10g	4arm PEG Azide (pentaerythritol), MW 10000	\$1,040
A7185-1	4ARM-AZIDE-20K,1g	4arm PEG Azide (pentaerythritol), MW 20000	\$130
A7185-10	4ARM-AZIDE-20K,10g	4arm PEG Azide (pentaerythritol), MW 20000	\$1,040


4ARM PEG ACETIC ACID (PENTAERYTHRITOL)

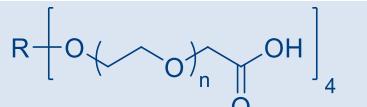

R = Pentaerythritol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7109-1	4ARM-CM-5000,1g	4arm PEG Acetic Acid (pentaerythritol), MW 5000	\$100
A7109-10	4ARM-CM-5000,10g	4arm PEG Acetic Acid (pentaerythritol), MW 5000	\$800
A7027-1	4ARM-CM-10K,1g	4arm PEG Acetic Acid (pentaerythritol), MW 10000	\$100
A7027-10	4ARM-CM-10K,10g	4arm PEG Acetic Acid (pentaerythritol), MW 10000	\$800
A7028-1	4ARM-CM-20K,1g	4arm PEG Acetic Acid (pentaerythritol), MW 20000	\$100
A7028-10	4ARM-CM-20K,10g	4arm PEG Acetic Acid (pentaerythritol), MW 20000	\$800
A7066-1	4ARM-CM-40K,1g	4arm PEG Acetic Acid (pentaerythritol), MW 40000	\$110
A7066-10	4ARM-CM-40K,10g	4arm PEG Acetic Acid (pentaerythritol), MW 40000	\$880

4ARM PEG ISOCYANATE (PENTAERYTHRITOL)

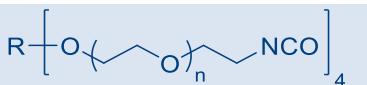


R = Pentaerythritol core structure

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A7106-1	4ARM-NCO-10K,1g	4arm PEG Isocyanate (pentaerythritol), MW 10000	\$150
* A7106-10	4ARM-NCO-10K,10g	4arm PEG Isocyanate (pentaerythritol), MW 10000	\$1,200
* A7111-1	4ARM-NCO-20K,1g	4arm PEG Isocyanate (pentaerythritol), MW 20000	\$150
* A7111-10	4ARM-NCO-20K,10g	4arm PEG Isocyanate (pentaerythritol), MW 20000	\$1,200

4ARM PEG MALEIMIDE (PENTAERYTHRITOL)

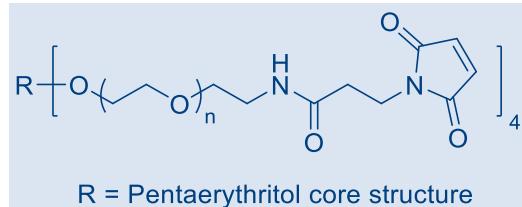

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7018-1	4ARM-MAL-10K,1g	4arm PEG Maleimide (pentaerythritol), MW 10000	\$130
A7018-10	4ARM-MAL-10K,10g	4arm PEG Maleimide (pentaerythritol), MW 10000	\$1040
A7029-1	4ARM-MAL-20K,1g	4arm PEG Maleimide (pentaerythritol), MW 20000	\$130
A7029-10	4ARM-MAL-20K,10g	4arm PEG Maleimide (pentaerythritol), MW 20000	\$1040
A7067-1	4ARM-MAL-40K,1g	4arm PEG Maleimide (pentaerythritol), MW 40000	\$140
A7067-10	4ARM-MAL-40K,10g	4arm PEG Maleimide (pentaerythritol), MW 40000	\$1,120


4ARM PEG NITROPHENYL CARBONATE (PENTAERYTHRITOL)

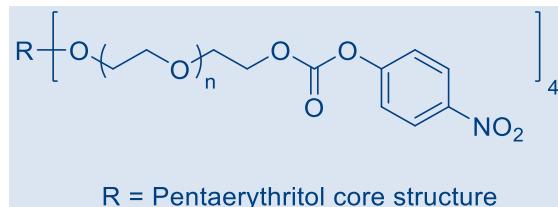

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A7080-1	4ARM-NPC-10K,1g	4arm PEG Nitrophenyl Carbonate (pentaerythritol), MW 10000	\$90
* A7080-10	4ARM-NPC-10K,10g	4arm PEG Nitrophenyl Carbonate (pentaerythritol), MW 10000	\$720


4ARM PEG SUCCINIMIDYL CARBOXYMETHYL ESTER (PENTAERYTHRITOL)

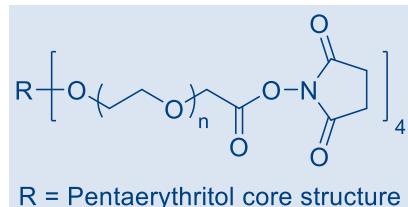

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7014-1	4ARM-SCM-5000,1g	4arm PEG Succinimidyl Carboxymethyl Ester (pentaerythritol), MW 5000	\$110
A7014-10	4ARM-SCM-5000,10g	4arm PEG Succinimidyl Carboxymethyl Ester (pentaerythritol), MW 5000	\$880
A7015-1	4ARM-SCM-10K,1g	4arm PEG Succinimidyl Carboxymethyl Ester (pentaerythritol), MW 10000	\$110
A7015-10	4ARM-SCM-10K,10g	4arm PEG Succinimidyl Carboxymethyl Ester (pentaerythritol), MW 10000	\$880

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7038-1	4ARM-SCM-20K,1g	4arm PEG Succinimidyl Carboxymethyl Ester (pentaerythritol), MW 20000	\$110
A7038-10	4ARM-SCM-20K,10g	4arm PEG Succinimidyl Carboxymethyl Ester (pentaerythritol), MW 20000	\$880
A7070-1	4ARM-SCM-40K,1g	4arm PEG Succinimidyl Carboxymethyl Ester (pentaerythritol), MW 40000	\$120
A7070-10	4ARM-SCM-40K,10g	4arm PEG Succinimidyl Carboxymethyl Ester (pentaerythritol), MW 40000	\$960

4ARM PEG SUCCINIMIDYL GLUTARAMIDE (PENTAERYTHRITOL)

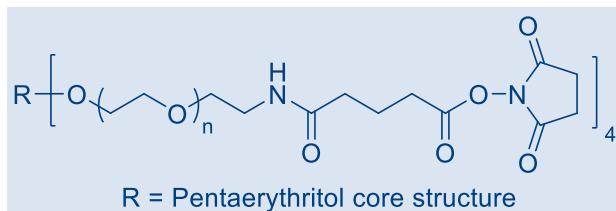


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7013-1	4ARM-SGA-20K,1g	4arm PEG Succinimidyl Glutaramide (pentaerythritol), MW 20000	\$120
A7013-10	4ARM-SGA-20K,10g	4arm PEG Succinimidyl Glutaramide (pentaerythritol), MW 20000	\$960

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A7110-1	4ARM-SGA-10K,1g	4arm PEG Succinimidyl Glutaramide (pentaerythritol), MW 10000	\$120
* A7110-10	4ARM-SGA-10K,10g	4arm PEG Succinimidyl Glutaramide (pentaerythritol), MW 10000	\$960

4ARM PEG SUCCINIMIDYL CARBONATE (PENTAERYTHRITOL)

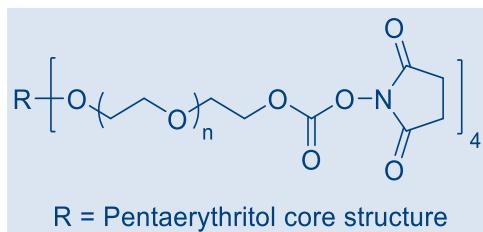


PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A7184-1	4ARM-SC-10K,1g	4ARM PEG Succinimidyl Carbonate (pentaerythritol), MW 10000	\$120
* A7184-10	4ARM-SC-10K,10g	4ARM PEG Succinimidyl Carbonate (pentaerythritol), MW 10000	\$960

Please inquire if the PEG derivative that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.

4ARM PEG SUCCINIMIDYL SUCCINATE (PENTAERYTHRITOL)

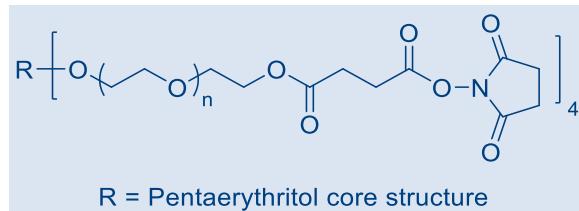


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7006-1	4ARM-SS-10K,1g	4arm PEG Succinimidyl Succinate (pentaerythritol), MW 10000	\$80
A7006-10	4ARM-SS-10K,10g	4arm PEG Succinimidyl Succinate (pentaerythritol), MW 10000	\$640
A7030-1	4ARM-SS-20K,1g	4arm PEG Succinimidyl Succinate (pentaerythritol), MW 20000	\$80
A7030-10	4ARM-SS-20K,10g	4arm PEG Succinimidyl Succinate (pentaerythritol), MW 20000	\$640

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A7069-1	4ARM-SS-40K,1g	4arm PEG Succinimidyl Succinate (pentaerythritol), MW 40000	\$90
* A7069-10	4ARM-SS-40K,10g	4arm PEG Succinimidyl Succinate (pentaerythritol), MW 40000	\$720

4ARM PEG SUCCINIMIDYL GLUTARATE (PENTAERYTHRITOL)

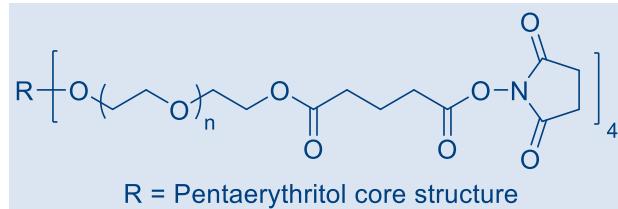


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7031-1	4ARM-SG-10K,1g	4arm PEG Succinimidyl Glutarate (pentaerythritol), MW 10000	\$80
A7031-10	4ARM-SG-10K,10g	4arm PEG Succinimidyl Glutarate (pentaerythritol), MW 10000	\$640
A7010-1	4ARM-SG-20K,1g	4arm PEG Succinimidyl Glutarate (pentaerythritol), MW 20000	\$80
A7010-10	4ARM-SG-20K,10g	4arm PEG Succinimidyl Glutarate (pentaerythritol), MW 20000	\$640

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
* A7017-1	4ARM-SG-40K,1g	4arm PEG Succinimidyl Glutarate (pentaerythritol), MW 40000	\$90
* A7017-10	4ARM-SG-40K,10g	4arm PEG Succinimidyl Glutarate (pentaerythritol), MW 40000	\$720

4ARM PEG THIOL (PENTAERYTHRITOL)

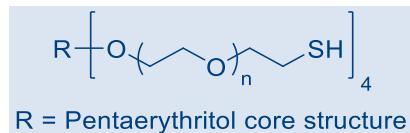

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7002-1	4ARM-SH-5000,1g	4arm PEG Thiol (pentaerythritol), MW 5000	\$80
A7002-10	4ARM-SH-5000,10g	4arm PEG Thiol (pentaerythritol), MW 5000	\$640
A7008-1	4ARM-SH-10K,1g	4arm PEG Thiol (pentaerythritol), MW 10000	\$80
A7008-10	4ARM-SH-10K,10g	4arm PEG Thiol (pentaerythritol), MW 10000	\$640
A7039-1	4ARM-SH-20K,1g	4arm PEG Thiol (pentaerythritol), MW 20000	\$80
A7039-10	4ARM-SH-20K,10g	4arm PEG Thiol (pentaerythritol), MW 20000	\$640

4ARM PEG VINYLSULFONE (PENTAERYTHRITOL)

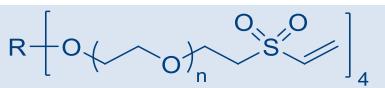


R = Pentaerythritol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A7005-1	4ARM-VS-10K,1g	4arm PEG Vinylsulfone (pentaerythritol), MW 10000	\$150
A7005-10	4ARM-VS-10K,10g	4arm PEG Vinylsulfone (pentaerythritol), MW 10000	\$1,200
A7025-1	4ARM-VS-20K,1g	4arm PEG Vinylsulfone (pentaerythritol), MW 20000	\$150
A7025-10	4ARM-VS-20K,10g	4arm PEG Vinylsulfone (pentaerythritol), MW 20000	\$1,200

6ARM-PEG DERIVATIVES (DIPENTAERYTRITOL CORE)

6ARM PEG AMINE (DIPENTAERYTHRITOL)

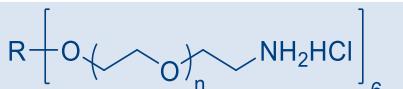

R = Dipentaerythritol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A9005-1	6ARM(DP)-NH2HCl-15K,1g	6arm PEG Amine (dipentaerythritol), HCl salt, MW 15000	\$80
A9005-10	6ARM(DP)-NH2HCl-15K,10g	6arm PEG Amine (dipentaerythritol), HCl salt, MW 15000	\$640

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A9026-1	6ARM(DP)-NH2HCl-30K,1g	6arm PEG Amine (dipentaerythritol), HCl salt, MW 30000	\$80
*A9026-10	6ARM(DP)-NH2HCl-30K,10g	6arm PEG Amine (dipentaerythritol), HCl salt, MW 30000	\$640

Please inquire if the PEG derivative that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.


8ARM-PEG DERIVATIVES (HEXAGLYCEROL CORE)

8ARM PEG AMINE (HEXAGLYCEROL)



ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A8008-1	8ARM-NH2HCl-10K,1g	8arm PEG Amine (hexaglycerol), HCl Salt, MW 10000	\$50
A8008-10	8ARM-NH2HCl-10K,10g	8arm PEG Amine (hexaglycerol), HCl Salt, MW 10000	\$400
A8009-1	8ARM-NH2HCl-20K,1g	8arm PEG Amine (hexaglycerol), HCl Salt, MW 20000	\$50
A8009-10	8ARM-NH2HCl-20K,10g	8arm PEG Amine (hexaglycerol), HCl Salt, MW 20000	\$400
A8012-1	8ARM-NH2HCl-40K,1g	8arm PEG Amine (hexaglycerol), HCl Salt, MW 40000	\$50
A8012-10	8ARM-NH2HCl-40K,10g	8arm PEG Amine (hexaglycerol), HCl Salt, MW 40000	\$400

8ARM PEG ACRYLATE (HEXAGLYCEROL)

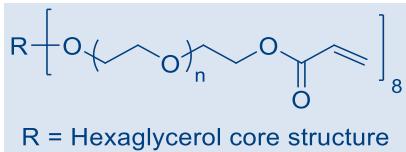

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A8002-1	8ARM-ACLT-10K,1g	8arm PEG Acrylate (hexaglycerol), MW 10000	\$80
A8002-10	8ARM-ACLT-10K,10g	8arm PEG Acrylate (hexaglycerol), MW 10000	\$640
A8005-1	8ARM-ACLT-20K,1g	8arm PEG Acrylate (hexaglycerol), MW 20000	\$80
A8005-10	8ARM-ACLT-20K,10g	8arm PEG Acrylate (hexaglycerol), MW 20000	\$640


8ARM PEG ACETIC ACID (HEXAGLYCEROL)



PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A8016-1	8ARM-CM-10K,1g	8arm PEG Acetic Acid (hexaglycerol), MW 10000	\$100
*A8016-10	8ARM-CM-10K,10g	8arm PEG Acetic Acid (hexaglycerol), MW 10000	\$800
*A8017-1	8ARM-CM-20K,1g	8arm PEG Acetic Acid (hexaglycerol), MW 20000	\$100
*A8017-10	8ARM-CM-20K,10g	8arm PEG Acetic Acid (hexaglycerol), MW 20000	\$800
*A8018-1	8ARM-CM-40K,1g	8arm PEG Acetic Acid (hexaglycerol), MW 40000	\$110
*A8018-10	8ARM-CM-40K,10g	8arm PEG Acetic Acid (hexaglycerol), MW 40000	\$880

8ARM PEG MALEIMIDE (HEXAGLYCEROL)

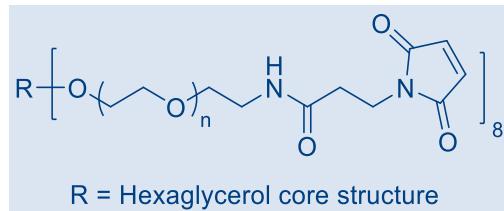


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A8027-1	8ARM-MAL-10K,1g	8arm PEG Maleimide (hexaglycerol), MW 10000	\$130
A8027-10	8ARM-MAL-10K,10g	8arm PEG Maleimide (hexaglycerol), MW 10000	\$1040
A8028-1	8ARM-MAL-20K,1g	8arm PEG Maleimide (hexaglycerol), MW 20000	\$130
A8028-10	8ARM-MAL-20K,10g	8arm PEG Maleimide (hexaglycerol), MW 20000	\$1040
A8029-1	8ARM-MAL-40K,1g	8arm PEG Maleimide (hexaglycerol), MW 40000	\$140
A8029-10	8ARM-MAL-40K,10g	8arm PEG Maleimide (hexaglycerol), MW 40000	\$1,120

8ARM PEG SUCCINIMIDYL SUCCINATE (HEXAGLYCEROL)

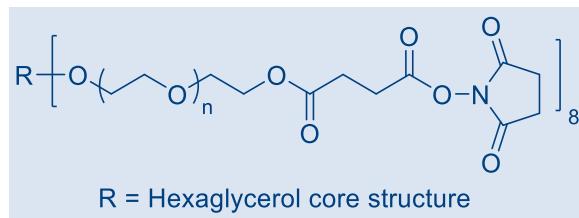


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A8019-1	8ARM-SS-10K,1g	8arm PEG Succinimidyl Succinate (hexaglycerol), MW 10000	\$80
A8019-10	8ARM-SS-10K,10g	8arm PEG Succinimidyl Succinate (hexaglycerol), MW 10000	\$640
A8010-1	8ARM-SS-15K,1g	8arm PEG Succinimidyl Succinate (hexaglycerol), MW 15000	\$80
A8010-10	8ARM-SS-15K,10g	8arm PEG Succinimidyl Succinate (hexaglycerol), MW 15000	\$640
A8020-1	8ARM-SS-20K,1g	8arm PEG Succinimidyl Succinate (hexaglycerol), MW 20000	\$80
A8020-10	8ARM-SS-20K,10g	8arm PEG Succinimidyl Succinate (hexaglycerol), MW 20000	\$640
A8021-1	8ARM-SS-40K,1g	8arm PEG Succinimidyl Succinate (hexaglycerol), MW 40000	\$90
A8021-10	8ARM-SS-40K,10g	8arm PEG Succinimidyl Succinate (hexaglycerol), MW 40000	\$720

8ARM PEG SUCCINIMIDYL GLUTARATE (HEXAGLYCEROL)

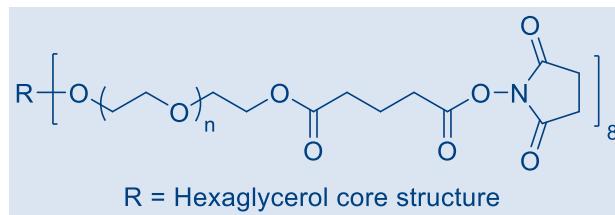


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A8013-1	8ARM-SG-10K,1g	8arm PEG Succinimidyl Glutarate (hexaglycerol), MW 10000	\$80
A8013-10	8ARM-SG-10K,10g	8arm PEG Succinimidyl Glutarate (hexaglycerol), MW 10000	\$640
A8001-1	8ARM-SG-15K,1g	8arm PEG Succinimidyl Glutarate (hexaglycerol), MW 15000	\$80

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A8001-10	8ARM-SG-15K,10g	8arm PEG Succinimidyl Glutarate (hexaglycerol), MW 15000	\$640
A8022-1	8ARM-SG-20K,1g	8arm PEG Succinimidyl Glutarate (hexaglycerol), MW 20000	\$80
A8022-10	8ARM-SG-20K,10g	8arm PEG Succinimidyl Glutarate (hexaglycerol), MW 20000	\$640
A8023-1	8ARM-SG-40K,1g	8arm PEG Succinimidyl Glutarate (hexaglycerol), MW 40000	\$100
A8023-10	8ARM-SG-40K,10g	8arm PEG Succinimidyl Glutarate (hexaglycerol), MW 40000	\$800

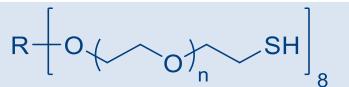8ARM PEG THIOL (HEXAGLYCEROL)


R = Hexaglycerol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A8004-1	8ARM-SH-10K,1g	8arm PEG Thiol (hexaglycerol), MW 10000	\$120
A8004-10	8ARM-SH-10K,10g	8arm PEG Thiol (hexaglycerol), MW 10000	\$960
A8007-1	8ARM-SH-20K,1g	8arm PEG Thiol (hexaglycerol), MW 20000	\$80
A8007-10	8ARM-SH-20K,10g	8arm PEG Thiol (hexaglycerol), MW 20000	\$640

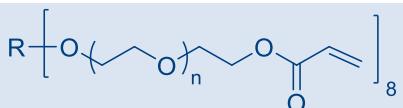8ARM-PEG DERIVATIVES (TRIPENTAERYTHRITOL CORE)**8ARM PEG ACRYLATE (TRIPENTAERYTHRITOL)**

R = Tripentaerythritol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A10021-1	8ARM(TP)-ACLT-10K,1g	8arm PEG Acrylate (tripentaerythritol), MW 10000	\$80
A10021-10	8ARM(TP)-ACLT-10K,10g	8arm PEG Acrylate (tripentaerythritol), MW 10000	\$640
A10016-1	8ARM(TP)-ACLT-20K,1g	8arm PEG Acrylate (tripentaerythritol), MW 20000	\$80
A10016-10	8ARM(TP)-ACLT-20K,10g	8arm PEG Acrylate (tripentaerythritol), MW 20000	\$640

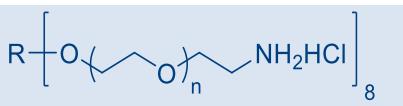
8ARM PEG AMINE (TRIPENTAERYTHRITOL)

R = Tripentaerythritol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A10001-1	8ARM(TP)-NH2HCl-10K,1g	8arm PEG Amine (tripentaerythritol), HCl Salt, MW 10000	\$50
A10001-10	8ARM(TP)-NH2HCl-10K,10g	8arm PEG Amine (tripentaerythritol), HCl Salt, MW 10000	\$400
A10002-1	8ARM(TP)-NH2HCl-20K,1g	8arm PEG Amine (tripentaerythritol), HCl Salt, MW 20000	\$50
A10002-10	8ARM(TP)-NH2HCl-20K,10g	8arm PEG Amine (tripentaerythritol), HCl Salt, MW 20000	\$400
A10003-1	8ARM(TP)-NH2HCl-40K,1g	8arm PEG Amine (tripentaerythritol), HCl Salt, MW 40000	\$50

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A10003-10	8ARM(TP)-NH2HCl-40K,10g	8arm PEG Amine (tripentaerythritol), HCl Salt, MW 40000	\$400

8ARM PEG ACETIC ACID (TRIPENTAERYTHRITOL)

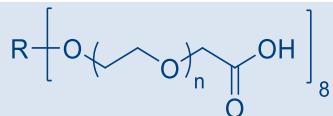


R = Tripentaerythritol core structure

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A10004-1	8ARM(TP)-CM-10K,1g	8arm PEG Acetic Acid (tripentaerythritol), MW 10000	\$100
*A10004-10	8ARM(TP)-CM-10K,10g	8arm PEG Acetic Acid (tripentaerythritol), MW 10000	\$800
*A10005-1	8ARM(TP)-CM-20K,1g	8arm PEG Acetic Acid (tripentaerythritol), MW 20000	\$100
*A10005-10	8ARM(TP)-CM-20K,10g	8arm PEG Acetic Acid (tripentaerythritol), MW 20000	\$800
*A10006-1	8ARM(TP)-CM-40K,1g	8arm PEG Acetic Acid (tripentaerythritol), MW 40000	\$110
*A10006-10	8ARM(TP)-CM-40K,10g	8arm PEG Acetic Acid (tripentaerythritol), MW 40000	\$880


8ARM PEG MALEIMIDE (TRIPENTAERYTHRITOL)

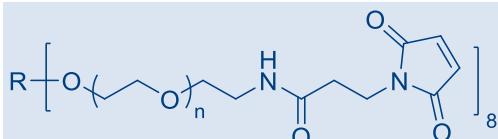

R = Tripentaerythritol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A10018-1	8ARM(TP)-MAL-10K,1g	8arm PEG Maleimide (tripentaerythritol), MW 10000	\$130
A10018-10	8ARM(TP)-MAL-10K,10g	8arm PEG Maleimide (tripentaerythritol), MW 10000	\$1,040
A10019-1	8ARM(TP)-MAL-20K,1g	8arm PEG Maleimide (tripentaerythritol), MW 20000	\$130
A10019-10	8ARM(TP)-MAL-20K,10g	8arm PEG Maleimide (tripentaerythritol), MW 20000	\$1,040
A10020-1	8ARM(TP)-MAL-40K,1g	8arm PEG Maleimide (tripentaerythritol), MW 40000	\$140
A10020-10	8ARM(TP)-MAL-40K,10g	8arm PEG Maleimide (tripentaerythritol), MW 40000	\$1,120

8ARM PEG NORBORNENE (TRIPENTAERYTHRITOL)

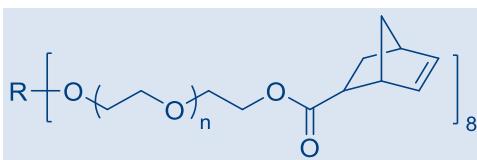


R = Tripentaerythritol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A10037-1	8ARM(TP)-NB-20K,1g	8arm PEG Norbornene (tripentaerythritol), MW 20000	\$150

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A10037-10	8ARM(TP)-NB-20K,10g	8arm PEG Norbornene (tripentaerythritol), MW 20000	\$1,200

8ARM PEG SUCCINIMIDYL SUCCINATE (TRIPENTAERYTHRITOL)

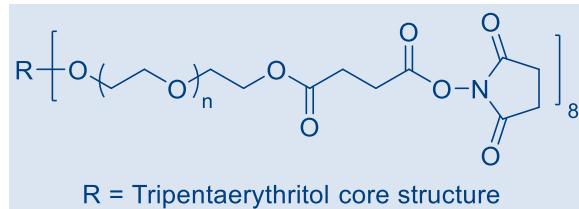


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A10007-1	8ARM(TP)-SS-10K,1g	8arm PEG Succinimidyl Succinate (tripentaerythritol), MW 10000	\$80
A10007-10	8ARM(TP)-SS-10K,10g	8arm PEG Succinimidyl Succinate (tripentaerythritol), MW 10000	\$640

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A10028-1	8ARM(TP)-SS-15K,1g	8arm PEG Succinimidyl Succinate (tripentaerythritol), MW 15000	\$80
*A10028-10	8ARM(TP)-SS-15K,10g	8arm PEG Succinimidyl Succinate (tripentaerythritol), MW 15000	\$640
*A10008-1	8ARM(TP)-SS-20K,1g	8arm PEG Succinimidyl Succinate (tripentaerythritol), MW 20000	\$80
*A10008-10	8ARM(TP)-SS-20K,10g	8arm PEG Succinimidyl Succinate (tripentaerythritol), MW 20000	\$640
*A10009-1	8ARM(TP)-SS-40K,1g	8arm PEG Succinimidyl Succinate (tripentaerythritol), MW 40000	\$90
*A10009-10	8ARM(TP)-SS-40K,10g	8arm PEG Succinimidyl Succinate (tripentaerythritol), MW 40000	\$720

8ARM PEG SUCCINIMIDYL GLUTARATE (TRIPENTAERYTHRITOL)

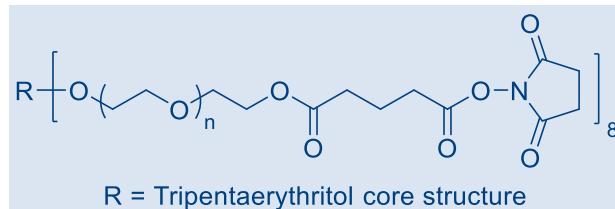


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A10010-1	8ARM(TP)-SG-10K,1g	8arm PEG Succinimidyl Glutarate (tripentaerythritol), MW 10000	\$80
A10010-10	8ARM(TP)-SG-10K,10g	8arm PEG Succinimidyl Glutarate (tripentaerythritol), MW 10000	\$640
A10012-1	8ARM(TP)-SG-20K,1g	8arm PEG Succinimidyl Glutarate (tripentaerythritol), MW 20000	\$80
A10012-10	8ARM(TP)-SG-20K,10g	8arm PEG Succinimidyl Glutarate (tripentaerythritol), MW 20000	\$640

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A10011-1	8ARM(TP)-SG-15K,1g	8arm PEG Succinimidyl Glutarate (tripentaerythritol), MW 15000	\$80
*A10011-10	8ARM(TP)-SG-15K,10g	8arm PEG Succinimidyl Glutarate (tripentaerythritol), MW 15000	\$640
*A10013-1	8ARM(TP)-SG-40K,1g	8arm PEG Succinimidyl Glutarate (tripentaerythritol), MW 40000	\$90

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A10013-10	8ARM(TP)-SG-40K,10g	8arm PEG Succinimidyl Glutarate (tripentaerythritol), MW 40000	\$720


8ARM PEG THIOL (TRIPENTAERYTHRITOL)

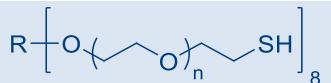

R = Tripentaerythritol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A10022-1	8ARM(TP)-SH-10K,1g	8arm PEG Thiol (tripentaerythritol), MW 10000	\$120
A10022-10	8ARM(TP)-SH-10K,10g	8arm PEG Thiol (tripentaerythritol), MW 10000	\$960
A10023-1	8ARM(TP)-SH-20K,1g	8arm PEG Thiol (tripentaerythritol), MW 20000	\$80
A10023-10	8ARM(TP)-SH-20K,10g	8arm PEG Thiol (tripentaerythritol), MW 20000	\$640

8ARM PEG VINYLSULFONE (TRIPENTAERYTHRITOL)

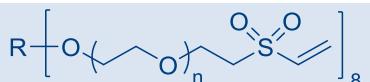

R = Tripentaerythritol core structure

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A10033-1	8ARM(TP)-VS-10K,1g	8arm PEG Vinylsulfone (tripentaerythritol), MW 10000	\$150
A10033-10	8ARM(TP)-VS-10K,10g	8arm PEG Vinylsulfone (tripentaerythritol), MW 10000	\$1,200
A10034-1	8ARM(TP)-VS-20K,1g	8arm PEG Vinylsulfone (tripentaerythritol), MW 20000	\$150
A10034-10	8ARM(TP)-VS-20K,10g	8arm PEG Vinylsulfone (tripentaerythritol), MW 20000	\$1,200
A10047-1	8ARM(TP)-VS-40K,1g	8arm PEG Vinylsulfone (tripentaerythritol), MW 40000	\$150
A10047-10	8ARM(TP)-VS-40K,10g	8arm PEG Vinylsulfone (tripentaerythritol), MW 40000	\$1,200

REFERENCES:


1. Hutanu, D., et al., Recent Applications of Polyethylene Glycols (PEGs) and PEG Derivatives. *Mod Chem Appl*, 2014.
2. Baker, A., et al., Stable oxime-crosslinked hyaluronan-based hydrogel as a biomimetic vitreous substitute, *Biomaterials*, 2021, V. 271.
3. Lee, Y.Y., et al, Long-acting nanoparticulate DNase-1 for effective suppression of SARS-CoV-2-mediated neutrophil activities and cytokine storm, *Biomaterials*, 2021, 267, 120389.
4. Ding, Y. et al, Tethering transforming growth factor β 1 to soft hydrogels guides vascular smooth muscle commitment from human mesenchymal stem cells, *Acta Biomaterialia*, 2020, V. 105, P. 68-77.
5. Newland, B., et al., Static and dynamic 3D culture of neural precursor cells on macroporous cryogel microcarriers, *MethodsX*, 2020, V.7.
6. Bachmann, D., et al., Retargeting of UniCAR T cells with an in vivo synthesized target module directed against CD19 positive tumor cells, *Oncotarget*, 2018, 9(7), p.7487.
7. Mou, C., et al., Electrochemical-mediated gelation of catechol-bearing hydrogels based on multimodal crosslinking, *Journal of Materials Chemistry B*, 2019.
8. Ding, Y., et al., Biomimetic soft fibrous hydrogels for contractile and pharmacologically responsive smooth muscle, *Acta biomaterialia*, 2018.

9. Uppal, G., et al., Tissue Failure Propagation as Mediated by Circulatory Flow, *Biophysical Journal*, 2020, 119(12), P. 2573-2583.
10. Wang, J., et al., An injectable PEG hydrogel controlling neurotrophin-3 release by affinity peptides, *Journal of Controlled Release*, 2021, 330, P. 575-586.
11. Kumar, M., et al., A fully defined matrix to support a pluripotent stem cell derived multi-cell-liver steatohepatitis and fibrosis model, *Biomaterials*, 2021, 121006.
12. Na, K.-S., et al., Effect of mesenchymal stromal cells encapsulated within polyethylene glycol-collagen hydrogels formed in situ on alkali-burned corneas in an ex vivo organ culture model, *Cytotherapy*, 2021, V. 23 (6).
13. DeForest, C.A., E.A. Sims, and K.S. Anseth, Peptide-Functionalized Click Hydrogels with Independently Tunable Mechanics and Chemical Functionality for 3D Cell Culture. *Chemistry of Materials*, 2010, 22(16): p. 4783-4790.
14. Chong, Y., et al., The in vitro and in vivo toxicity of graphene quantum dots. *Biomaterials*, 2014, 35(19).
15. Wang, M., et al., A surface convertible nanoplatform with enhanced mitochondrial targeting for tumor photothermal therapy, *Colloids and Surfaces B: Biointerfaces*, 2020, V.189.
16. Chauhan, N., et al., Dexamethasone-loaded, injectable pullulan-poly(ethylene glycol) hydrogels for bone tissue regeneration in chronic inflammatory conditions, *Materials Science and Engineering: C*, 2021, V. 130.
17. Schroeder, M. E., et al., Collagen networks within 3D PEG hydrogels support valvular interstitial cell matrix mineralization, *Acta Biomaterialia*, 2021, V. 119, P. 197-210.
18. Stock, A.A., et al., Conformal Coating of Stem Cell-Derived Islets for β Cell Replacement in Type 1 Diabetes, *Stem Cell Reports*, 2020, 14(1), P. 91-104
19. Carleton, M. M., et al., Methacrylic acid-based hydrogels enhance skeletal muscle regeneration after volumetric muscle loss in mice, *Biomaterials*, 2021, V. 275.
20. Swaminathan, G., et al., Effect of substrate stiffness on human intestinal enteroids' infectivity by enteroaggregative *Escherichia coli*, *Acta Biomaterialia*, 2021, V. 132.
21. Kaphle, P., et al., The mechanical and pharmacological regulation of glioblastoma cell migration in 3D matrices, *Journal of cellular physiology*, 2019, 234(4):3948-60.
22. Buwalda, S.J., et al., Ultrafast in situ forming poly (ethylene glycol)-poly (amido amine) hydrogels with tunable drug release properties via controllable degradation rates, *European Journal of Pharmaceutics and Biopharmaceutics*, 2019, 139:232-9.
23. Schirmer, L., et al., Glycosaminoglycan-based hydrogels with programmable host reactions, *Biomaterials*, 2020.
24. Scott, R. A., et al., Substrate stiffness directs the phenotype and polarization state of cord blood derived macrophages, *Acta Biomaterialia*, 2021, V. 122, P. 220-235.
25. Tunn, I., et al., Bioinspired Histidine-Zn²⁺ Coordination for Tuning the Mechanical Properties of Self-Healing Coiled Coil Cross-Linked Hydrogels, *Biomimetics*, 2019, 4(1):25.
26. Li, H., et al., Synthesis of thiol-terminated PEG-functionalized POSS cross-linkers and fabrication of high-strength and hydrolytic degradable hybrid hydrogels in aqueous phase, *European Polymer Journal*, 2019, 116:74-83.
27. Manzoli, V., et al., Immunoisolation of murine islet allografts in vascularized sites through conformal coating with polyethylene glycol, *American Journal of Transplantation*, 2018, 18(3):590-603.
28. Rao, V.V., et al., Rescuing mesenchymal stem cell regenerative properties on hydrogel substrates post serial expansion, *Bioengineering & translational medicine*, 2019.
29. Nguyen, M.K., et al., RNA interfering molecule delivery from in situ forming biodegradable hydrogels for enhancement of bone formation in rat calvarial bone defects, *Acta Biomaterialia*, V. 75, 2018, P. 105-114.
30. Tomaszewski, C. E., et al., Sequestered cell-secreted extracellular matrix proteins improve murine folliculogenesis and oocyte maturation for fertility preservation, *Acta Biomaterialia*, 2021.
31. Lau, C. M. L., et al., Local environment-dependent kinetics of ester hydrolysis revealed by direct ¹H NMR measurement of degrading hydrogels, *Acta Biomaterialia*, 2020, V. 101, P. 219-226.
32. Yu, W., et al., Light-Addressable Nanocomposite Hydrogels Allow Plasmonic Actuation and In Situ Temperature Monitoring in 3D Cell Matrices. *Advanced Functional Materials*. 2021, 2108234.
33. Na, K.-S., et al., Effect of mesenchymal stromal cells encapsulated within polyethylene glycol-collagen hydrogels formed in situ on alkali-burned corneas in an ex vivo organ culture model, *Cytotherapy*, 2021, V. 23 (6).
34. Sofman, M., et al., A modular polymer microbead angiogenesis scaffold to characterize the effects of adhesion ligand density on angiogenic sprouting, *Biomaterials*, 2021, 264, 120231.

35. Luo, Y., et al., Light-induced dynamic RGD pattern for sequential modulation of macrophage phenotypes, *Bioactive Materials*, 2021, V. 6(11), P. 4065-4072.
36. Khang, A., et al., Quantifying heart valve interstitial cell contractile state using highly tunable poly(ethylene glycol) hydrogels, *Acta Biomaterialia*, 2019, 96, p. 354-367.
37. Rutz, A. L., et al., Employing PEG crosslinkers to optimize cell viability in gel phase bioinks and tailor post printing mechanical properties, *Acta Biomaterialia*, 2019.
38. Dumont, C.M., et al., Aligned hydrogel tubes guide regeneration following spinal cord injury, *Acta Biomaterialia*, 2019.
39. Sun, X., et al., Three-dimensional bioprinting of multicell-laden scaffolds containing bone morphogenic protein-4 for promoting M2 macrophage polarization and accelerating bone defect repair in diabetes mellitus, *Bioactive Materials*, 6(3), 2021, P. 757-769.

MULTI-ARM HETEROBIFUNCTIONAL PEG DERIVATIVES


(HO)3-4ARMPEG-CM

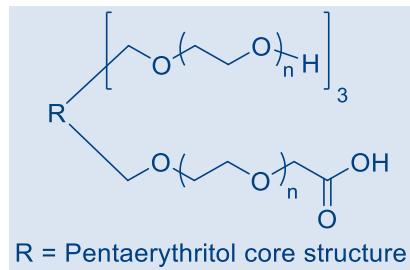

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A7147-1	(HO)3-4ARMPEG3000-CM,1g	4arm PEG, 3arm-Hydroxyl, 1arm-Acetic Acid, MW 3000	\$450
*A7147-5	(HO)3-4ARMPEG3000-CM,5g	4arm PEG, 3arm-Hydroxyl, 1arm-Acetic Acid, MW 3000	\$1,800
*A7142-1	(HO)3-4ARMPEG5000-CM,1g	4arm PEG, 3arm-Hydroxyl, 1arm-Acetic Acid, MW 5000	\$300
*A7142-5	(HO)3-4ARMPEG5000-CM,5g	4arm PEG, 3arm-Hydroxyl, 1arm-Acetic Acid, MW 5000	\$1,200
*A7047-1	(HO)3-4ARMPEG10K-CM,1g	4arm PEG, 3arm-Hydroxyl, 1arm-Acetic Acid, MW 10000	\$300
*A7047-5	(HO)3-4ARMPEG10K-CM,5g	4arm PEG, 3arm-Hydroxyl, 1arm-Acetic Acid, MW 10000	\$1,200

(HO)3-4ARMPEG-NH2HCL

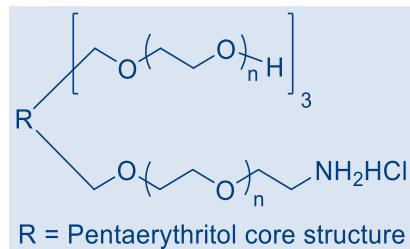


PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A7081-1	(HO)3-4ARMPEG10K-NH2HCl,1g	4arm PEG, 3arm-Hydroxyl, 1arm-Amine, HCl salt, MW 10000	\$500

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A7081-5	(HO)3-4ARMPEG10K-NH2HCl,5g	4arm PEG, 3arm-Hydroxyl, 1arm-Amine, HCl salt, MW 10000	\$2,000
*A7148-1	(HO)3-4ARMPEG20K-NH2HCl,1g	4arm PEG, 3arm-Hydroxyl, 1arm-Amine, HCl salt, MW 20000	\$550
*A7148-5	(HO)3-4ARMPEG20K-NH2HCl,5g	4arm PEG, 3arm-Hydroxyl, 1arm-Amine, HCl salt, MW 20000	\$2,200


(HO)2-4ARMPEG-(CM)2

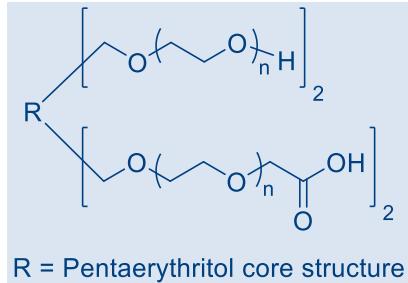

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A7048-1	(HO)2-4ARMPEG10K-(CM)2,1g	4arm PEG, 2arm-Hydroxyl, 2arm-Acetic Acid, MW 10000	\$350
*A7048-5	(HO)2-4ARMPEG10K-(CM)2,5g	4arm PEG, 2arm-Hydroxyl, 2arm-Acetic Acid, MW 10000	\$1,400


(HO)2-4ARMPEG-(NH2HCl)2

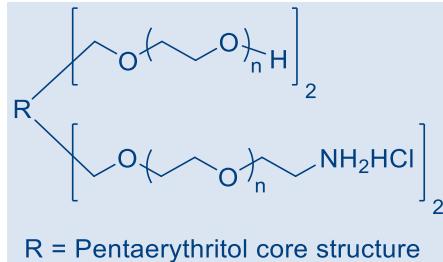

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A7077-1	(HO)2-4ARMPEG10K-(NH2HCl)2,1g	4arm-PEG, 2arm-Hydroxyl, 2arm-Amine, HCl Salt, MW 10000	\$700
*A7077-5	(HO)2-4ARMPEG10K-(NH2HCl)2,5g	4arm-PEG, 2arm-Hydroxyl, 2arm-Amine, HCl Salt, MW 10000	\$2,800
*A7149-1	(HO)2-4ARMPEG20K-(NH2HCl)2,1g	4arm PEG, 2arm-Hydroxyl, 2arm-Amine, HCl salt, MW 20000	\$750
*A7149-5	(HO)2-4ARMPEG20K-(NH2HCl)2,5g	4arm PEG, 2arm-Hydroxyl, 2arm-Amine, HCl salt, MW 20000	\$3,000


HO-4ARMPEG-(CM)3

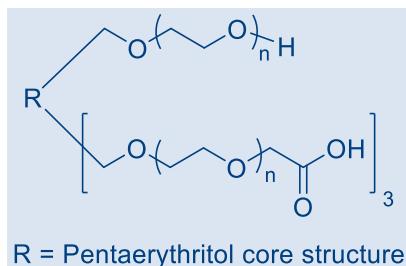

R = Pentaerythritol core structure

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A7108-1	HO-4ARMPEG5000-(CM)3,1g	4arm PEG, 1arm-Hydroxyl, 3arm-Acetic Acid, MW 5000	\$400
*A7108-5	HO-4ARMPEG5000-(CM)3,5g	4arm PEG, 1arm-Hydroxyl, 3arm-Acetic Acid, MW 5000	\$1,600

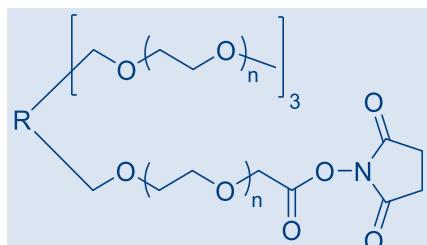
(CH₃O)3-4ARMPEG-SCM

R = Pentaerythritol core structure

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A7174-1	(CH ₃ O)3-4ARMPEG5000-SCM,1g	4arm PEG, 3arm Methoxy, 1arm Succinimidyl Carboxymethyl Ester, MW 5000	\$650
*A7174-5	(CH ₃ O)3-4ARMPEG5000-SCM,5g	4arm PEG, 3arm Methoxy, 1arm Succinimidyl Carboxymethyl Ester, MW 5000	\$2,600
*A7173-1	(CH ₃ O)3-4ARMPEG20K-SCM,1g	4arm PEG, 3arm Methoxy, 1arm Succinimidyl Carboxymethyl Ester, MW 20000	\$600
*A7173-5	(CH ₃ O)3-4ARMPEG20K-SCM,5g	4arm PEG, 3arm Methoxy, 1arm Succinimidyl Carboxymethyl Ester, MW 20000	\$2,400

(HO)7-8ARMPEG-CM

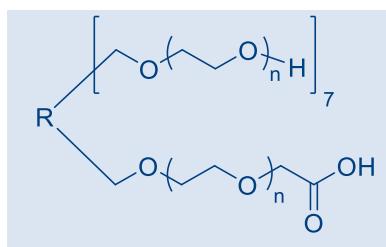


R = Hexaglycerol core structure

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A8042-1	(HO)7-8ARMPEG10K-CM,1g	8arm PEG (hexaglycerol), 7arm-Hydroxyl, 1arm-Acetic Acid, MW 10000	\$400
*A8042-5	(HO)7-8ARMPEG10K-CM,5g	8arm PEG (hexaglycerol), 7arm-Hydroxyl, 1arm-Acetic Acid, MW 10000	\$1,600
*A8043-1	(HO)7-8ARMPEG20K-CM,1g	8arm PEG (hexaglycerol), 7arm-Hydroxyl, 1arm-Acetic Acid, MW 20000	\$450
*A8043-5	(HO)7-8ARMPEG20K-CM,5g	8arm PEG (hexaglycerol), 7arm-Hydroxyl, 1arm-Acetic Acid, MW 20000	\$1,800

(HO)7-8ARMPEG(TP)-CM

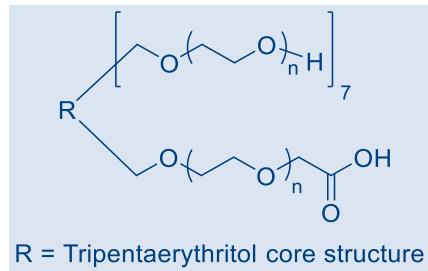

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*A10041-1	(HO)7-8ARM(TP)PEG10K-CM,1g	8arm PEG (tripentaerythritol), 7arm-Hydroxyl, 1arm-Acetic Acid, MW 10000	\$400
*A10041-5	(HO)7-8ARM(TP)PEG10K-CM,5g	8arm PEG (tripentaerythritol), 7arm-Hydroxyl, 1arm-Acetic Acid, MW 10000	\$1,600
*A10039-1	(HO)7-8ARM(TP)PEG20K-CM,1g	8arm PEG (tripentaerythritol), 7arm-Hydroxyl, 1arm-Acetic Acid, MW 20000	\$450
*A10039-5	(HO)7-8ARM(TP)PEG20K-CM,5g	8arm PEG (tripentaerythritol), 7arm-Hydroxyl, 1arm-Acetic Acid, MW 20000	\$1,800

REFERENCES:

1. Coppock, M.B., et al., A Universal Method for the Functionalization of Dyed Magnetic Microspheres with Peptides, *Methods*, 2019.
2. Coppock, M.B., et al., A novel discovery, maturation, and assay integration approach for the development of ruggedized multi-valent capture receptors exemplified against the chikungunya virus E2 protein, *Sensing and Bio-Sensing Research*, 2019, 22.
3. Liu, K.F., et al., A novel self-assembled pH-sensitive targeted nanoparticle platform based on antibody-4arm-polyethylene glycol-pterostilbene conjugates for co-delivery of anticancer drugs, *Journal of Materials Chemistry B*, 2018.
4. Dai, L., et al., Novel Multiarm Polyethylene glycol-Dihydroartemisinin Conjugates Enhancing Therapeutic Efficacy in Non-Small-Cell Lung Cancer, *Scientific Reports*, 2014, 4, 5871.

Please inquire if the PEG derivative that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.

BIODEGRADABLE PEG CO-POLYMERS

JenKem Technology offers high-quality biodegradable co-polymers. The interest in biodegradable polymers for the medical and pharmaceutical field stems from the need for extended release drug delivery systems, including for ADCs; and for degradable medical devices that do not require human intervention for removal from the body. JenKem® PEGs are ideal candidates as a polymeric backbone, as they have already been successfully applied for wound sealing and healing, controlled drug delivery, and tissue engineering.


Synthetic JenKem PEG co-polymers and PEG derivatives offer several advantages over natural biodegradable polymers:

- Highly hydrophilic
- Facile introduction of functional groups for drug attachment
- Production of JenKem PEGs back-integrated to ethylene oxide raw material
- Predictable lot-to-lot consistency and raw material sourcing
- No immunogenicity normally associated with natural polymers

Please see below a limited selection of JenKem® biodegradable PEG copolymers. Different combinations of MW and functionalities for the listed PEG co-polymers; and co-polymer derivatives with PGA (PolyGlycolide), PLL (Poly(L-lysine), or PLGA (Poly(lactic-co-glycolic acid)); are available by custom synthesis.


PEG-POLYLACTIDE BLOCK COPOLYMERS

4ARM PEG POLYLACTIDE BLOCK Co-POLYMERS

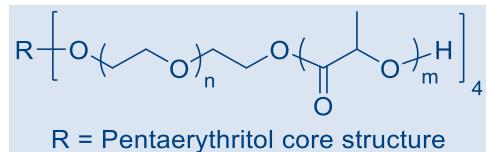


ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE(USD)
A9038-1	4ARM-PEG2500-b-PLA2500,1g	4arm PEG Polylactide Block Copolymer (pentaerythritol), PEG MW 2500, PLA MW 2500	\$500
A9038-10	4ARM-PEG2500-b-PLA2500,10g	4arm PEG Polylactide Block Copolymer (pentaerythritol), PEG MW 2500, PLA MW 2500	\$4,000

METHOXY PEG POLYLACTIDE BLOCK Co-POLYMERS

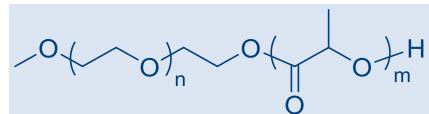

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A9019-1	MPEG5000-b-PLA2500,1g	Methoxy Polyethylene Glycol Polylactide Block Copolymer, PEG MW 5000, PLA MW 2500	\$100
A9019-10	MPEG5000-b-PLA2500,10g	Methoxy Polyethylene Glycol Polylactide Block Copolymer, PEG MW 5000, PLA MW 2500	\$800
A9037-1	MPEG5000-b-PLA5000,1g	Methoxy Polyethylene Glycol Polylactide Block Copolymer, PEG MW 5000, PLA MW 5000	\$100
A9037-10	MPEG5000-b-PLA5000,10g	Methoxy Polyethylene Glycol Polylactide Block Copolymer, PEG MW 5000, PLA MW 5000	\$800


MALEIMIDE PEG POLYLACTIDE BLOCK CO-POLYMERS

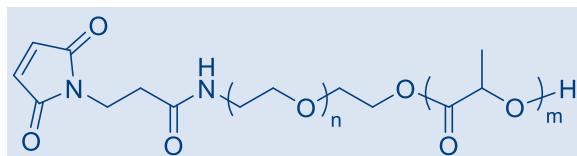

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
A9040-1	MAL-PEG2000-b-PLA2000,1g	Maleimide Polyethylene Glycol Polylactide block copolymer, PEG MW 2000, PLA MW 2000	\$450
A9040-10	MAL-PEG2000-b-PLA2000,10g	Maleimide Polyethylene Glycol Polylactide block copolymer, PEG MW 2000, PLA MW 2000	\$3,600
A9039-1	MAL-PEG5000-b-PLA20K,1g	Maleimide Polyethylene Glycol Polylactide block copolymer, PEG MW 5000, PLA MW 20000	\$500
A9039-10	MAL-PEG5000-b-PLA20K,10g	Maleimide Polyethylene Glycol Polylactide block copolymer, PEG MW 5000, PLA MW 20000	\$4,000


PEG-POLYCARROLACTONE BLOCK COPOLYMERS

4ARM PEG POLYCAPROLACTONE BLOCK CO-POLYMERS

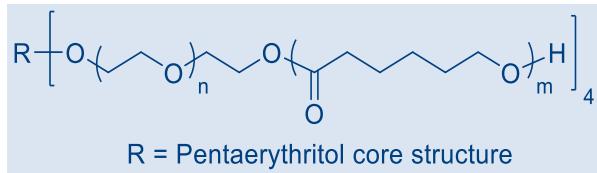

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE(USD)
A9012-1	4ARM-PEG2500-b-PCL2500,1g	4arm PEG Polycaprolactone Block Copolymer (pentaerythritol), PEG MW 2500, PCL MW 2500	\$550
A9012-10	4ARM-PEG2500-b-PCL2500,10g	4arm PEG Polycaprolactone Block Copolymer (pentaerythritol), PEG MW 2500, PCL MW 2500	\$4,400

METHOXY PEG POLYCAPROLACTONE BLOCK CO-POLYMERS

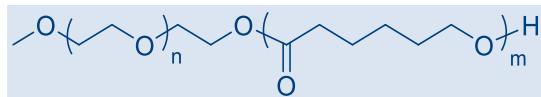

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE(USD)
A9033-1	MPEG2000-b-PCL2000,1g	Methoxy PEG Polycaprolactone Block Copolymer, PEG MW 2000, PCL MW 2000	\$200
A9033-10	MPEG2000-b-PCL2000,10g	Methoxy PEG Polycaprolactone Block Copolymer, PEG MW 2000, PCL MW 2000	\$1,600

REFERENCES:

- Wang, Y., et al., Dual micelles-loaded gelatin nanofibers and their application in lipopolysaccharide-induced periodontal disease, International journal of nanomedicine, 2019, 14:963.
- Richard, M., et al., Active cargo positioning in antiparallel transport networks, bioRxiv, 2019.

3. Xiong, Q., et al., Facile fabrication of reduction-responsive supramolecular nanoassemblies for co-delivery of doxorubicin and sorafenib towards hepatoma cells, *Frontiers in pharmacology*, 2018, 9, p.61.
4. Jin, Q., et al., Edaravone-Encapsulated Agonistic Micelles Rescue Ischemic Brain Tissue by Tuning Blood-Brain Barrier Permeability, *Theranostics*, 2017, 7(4):884.
5. Nguyen, F., et al., Enhanced Intratumoral Delivery of SN38 as a Tocopherol Oxyacetate Prodrug Using Nanoparticles in a Neuroblastoma Xenograft Model, *Clinical Cancer Research*, 2018, p.3811.
6. Liang, H., et al., Size-Shifting Micelle Nanoclusters Based on a Cross-Linked and pH-Sensitive Framework for Enhanced Tumor Targeting and Deep Penetration Features. *ACS applied materials & interfaces*, 2016, 8(16):10136-46.
7. Li, Q., et al., Development of reactive oxygen species (ROS)-responsive nanoplatform for targeted oral cancer therapy. *Journal of Materials Chemistry B*, 2016.
8. Schaedel, L., et al., Microtubules self-repair in response to mechanical stress, *Nature Mater.*, 2015, 14(11).
9. Bhajun, R., et al., A statistically inferred microRNA network identifies breast cancer target miR-940 as an actin cytoskeleton regulator, *Scientific Reports*, 2015, 5, 8336.
10. Portran, D., et al., Micropatterning Microtubules, *Methods in Cell Biology*, Matthieu P, Manuel T (Editors) Academic Press, 2014, 39-51.
11. Boujemaa-Paterski, R., et al., Directed actin assembly and motility, *Methods Enzymol*, 2014, 540: 283-300.
12. Vignaud, T., et al., Polyacrylamide hydrogel micropatterning, *Methods Cell Biol*, 2014, 120: 93-116.
13. Reymann, A.-C., et al, Nucleation geometry governs ordered actin networks structures, *Nature Materials*, 2010.
14. Kiss, A., et al., Nuclear Motility in Glioma Cells Reveals a Cell-Line Dependent Role of Various Cytoskeletal Components, *PLoS ONE*, 2014, 9(4): e93431.
15. Portran, D., et al., Quantification of MAP and molecular motor activities on geometrically controlled microtubule networks, *Cytoskeleton*, 2013, 70: 12-23.
16. Schiller, H.B., et al., β 1-and av-class integrins cooperate to regulate myosin II during rigidity sensing of fibronectin-based microenvironments, *Nature cell biology*, 2013, 15.6: 625-636.
17. Galland, R., et al., Fabrication of three-dimensional electrical connections by means of directed actin self-organization, *Nature materials*, 2013, 12.5 : 416-421.
18. Pitaval, A., et al., Cell shape and contractility regulate ciliogenesis in cell cycle-arrested cells, *JCB*, 2010, 191.

PEG RAW MATERIALS

JenKem Technology manufactures high quality polyethylene glycol raw materials with high purity and low polydispersity. JenKem Technology's in-house PEG production is back-integrated to the ethylene oxide starting material. JenKem Technology's PEG raw materials (PEG alcohols) are suitable for further derivatization into PEG derivatives.

JenKem Technology provides high quality linear methoxy and benzyl polyethylene glycol raw materials with high purity, low polydispersity, and low to no diol content. These PEG raw materials are end-capped with a methoxy group or benzyl group. Linear PEG raw materials are provided as bulk orders of 100g or more, with molecular weights of 5kDa, 10kDa, 20kDa, 30kDa, and 40kDa.

JenKem Technology now provides novel 8ARM(TP)-PEG raw materials with tripentaerythritol core. 8ARM(TP)-PEGs with tripentaerythritol core have a higher purity as evidenced by MALDI compared to the generic 8ARM-PEGs with a hexaglycerin core.

JenKem Technology reports the total molecular weight of multi-arm PEG raw materials as the sum of the molecular weights of all the arms. 4ARM-PEG-OH (4arm PEG hydroxyl) is prepared by ethoxylation of pentaerythritol. 6ARM-PEG-OH (6arm PEG hydroxyl) is prepared by ethoxylation of dipentaerythritol. 8ARM-PEG-OH (8arm PEG hydroxyl) is prepared either by ethoxylation of tripentaerythritol (TP), or hexaglycerol.

PEG raw materials with molecular weights and branching not listed in our catalog may be available by custom synthesis. Please inquire at tech@jenkemusa.com about pricing and availability.


LINEAR METHOXY PEG RAW MATERIALS



PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*P1005R-100	M-PEG-5000R,100g	Methoxy Polyethylene Glycol, MW 5000	\$775
*P1006R-100	M-PEG-10KR,100g	Methoxy Polyethylene Glycol, MW 10000	\$775
*P1009R-100	M-PEG-20KR,100g	Methoxy Polyethylene Glycol, MW 20000	\$925
*P1010R-100	M-PEG-30KR,100g	Methoxy Polyethylene Glycol, MW 30000	\$1075
*P1012R-100	M-PEG-40KR,100g	Methoxy Polyethylene Glycol, MW 40000	\$1,150

LINEAR BENZYL PEG RAW MATERIALS

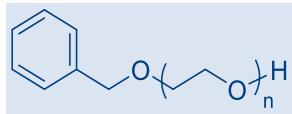


PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*P10001R-100	BENZYL-PEG-2000R,100g	Benzyl Polyethylene Glycol, MW 2000	\$775
*P10002R-100	BENZYL-PEG-5000R,100g	Benzyl Polyethylene Glycol, MW 5000	\$775

MULTI-ARM PEG RAW MATERIALS


4ARM-PEG RAW MATERIALS (PENTAERYTHRITOL CORE)

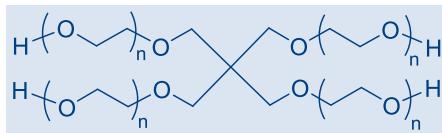

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
P2001-100	4ARM-PEG-2000,100g	4arm Polyethylene Glycol, MW 2000	\$450
*P2002R-100	4ARM-PEG-5000R,100g	4arm Polyethylene Glycol, MW 5000	\$625
*P2003R-100	4ARM-PEG-10KR,100g	4arm Polyethylene Glycol, MW 10000	\$625
*P2004R-100	4ARM-PEG-15KR,100g	4arm Polyethylene Glycol, MW 15000	\$625
*P2005R-100	4ARM-PEG-20KR,100g	4arm Polyethylene Glycol, MW 20000	\$625


6ARM-PEG RAW MATERIALS (DIPENTAERYTHRITOL CORE)

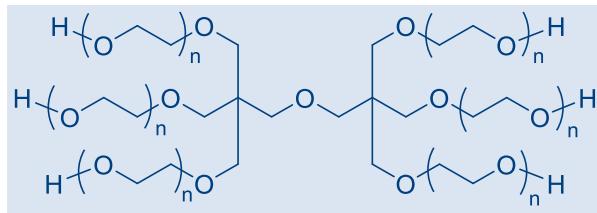

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*P6001R-100	6ARM(DP)-PEG-15KR,100g	6arm Poly(ethylene glycol) (dipentaerythritol), MW 15000	\$675
*P6003R-100	6ARM(DP)-PEG-30KR,100g	6arm Poly(ethylene glycol) (dipentaerythritol), MW 30000	\$675


8ARM-PEG RAW MATERIALS (HEXAGLYCEROL CORE)

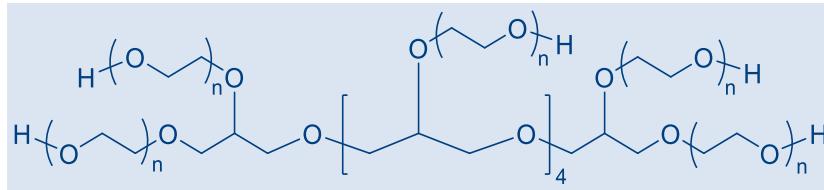

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*P3001R-100	8ARM-PEG-10KR,100g	8arm Polyethylene Glycol (hexaglycerol), MW 10000	\$775
*P3003R-100	8ARM-PEG-15KR,100g	8arm Polyethylene Glycol (hexaglycerol), MW 15000	\$775
*P3004R-100	8ARM-PEG-20KR,100g	8arm Polyethylene Glycol (hexaglycerol), MW 20000	\$775
*P3006R-100	8ARM-PEG-40KR,100g	8arm Polyethylene Glycol (hexaglycerol), MW 40000	\$775

8ARM-PEG RAW MATERIALS (TRIPENTAERYTHRITOL CORE)

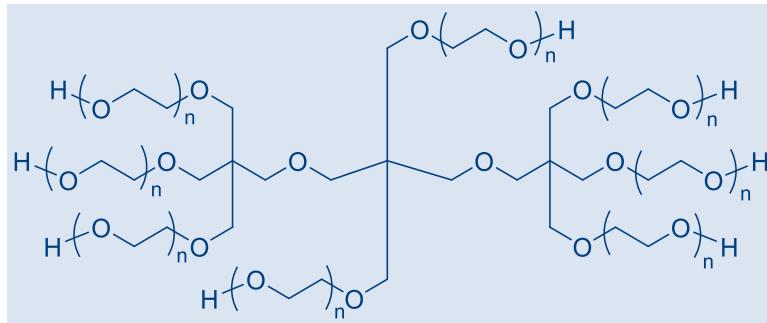

PRODUCTS MARKED WITH * ARE MADE TO ORDER AND MAY NOT BE IN STOCK. THE LEAD TIME FOR SYNTHESIS OF THESE PRODUCTS MAY RANGE FROM 4-10 WEEKS FROM THE DATE WHEN THE ORDER IS RECEIVED.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
*P4001R-100	8ARM(TP)-PEG-10KR,100g	8arm Polyethylene Glycol (tripentaerythritol), MW 10000	\$475
*P4003R-100	8ARM(TP)-PEG-15KR,100g	8arm Polyethylene Glycol (tripentaerythritol), MW 15000	\$475
*P4004R-100	8ARM(TP)-PEG-20KR,100g	8arm Polyethylene Glycol (tripentaerythritol), MW 20000	\$475
*P4006R-100	8ARM(TP)-PEG-40KR,100g	8arm Polyethylene Glycol (tripentaerythritol), MW 40000	\$625
*P4006R-100	8ARM(TP)-PEG-40KR,100g	8arm Polyethylene Glycol (tripentaerythritol), MW 40000	\$625

REFERENCES:

1. Yang, Z., et al., Defeating relapsed and refractory malignancies through a nano-enabled mitochondria-mediated respiratory inhibition and damage pathway, *Biomaterials*, 2020, V. 229.
2. Lee, S., et al., Inflammatory responses of macrophage-like RAW264.7 cells in a 3D hydrogel matrix to ultrasonicated schizophyllan, *Carbohydrate Polymers*, 2020, V. 229.
3. Kwak, H., et al., Colorimetric assay of tyrosinase inhibition using melanocyte laden hydrogel fabricated by digital light processing printing, *Journal of Industrial and Engineering Chemistry*, 2020, V. 84, P. 252-259.
4. Ma, H., et al., Bioorthogonal click chemistries enable simultaneous spatial patterning of multiple proteins to probe synergistic protein effects on fibroblast function, *Biomaterials*, 2020, 255:120205.
5. Eles, J.R., et al., Meningeal inflammatory response and fibrous tissue remodeling around intracortical implants: An in vivo two-photon imaging study, *Biomaterials*, 2019, 195, P. 111-12.
6. Xin, S., et al., Clickable PEG Hydrogel Microspheres as Building Blocks for 3D Bioprinting, *Biomaterials science*, 2019.
7. Tae, H., et al., β -glucan hybridized poly (ethylene glycol) microgels for macrophage-targeted protein delivery, *Journal of Industrial and Engineering Chemistry*, 2019.
8. Jones, C.E., et al., Stromal PTEN Regulates Extracellular Matrix Organization in the Mammary Gland, *Neoplasia*, 2019.
9. Shagan, A., et al., Near-Infrared Light Induced Phase Transition of Biodegradable Composites for On-Demand Healing and Drug Release, *ACS applied materials & interfaces*, 2018, 10(4), p.4131-4139.
10. Carthew, J., et al., Polyethylene glycol-gelatin hydrogels with tuneable stiffness prepared by horseradish peroxidase-activated tetrazine-norbornene ligation, *Journal of Materials Chemistry B*, 2018, 6:9.
11. Croitoru-Sadger, T., et al., Two-component cross-linkable gels for fabrication of solid oral dosage forms, *Journal of Controlled Release*, 2019.
12. Khodamoradi, M., et al., An electro-conductive hybrid scaffold as an artificial Bruch's membrane, *Materials Science and Engineering: C*, 2021, V. 126.
13. Machado Cruz, R., et al., Impact of polyethylene glycol polymers on the physicochemical properties and mucoadhesivity of itraconazole nanoparticles, *European Journal of Pharmaceutics and Biopharmaceutics*, 2019, 144.
14. Kwak, H., et al., Formation of a keratin layer with silk fibroin-polyethylene glycol composite hydrogel fabricated by digital light processing 3D printing, *Journal of Industrial and Engineering Chemistry*, 2018.
15. Saghiri, M.A., et al., Hydrogel Arrays and Choroidal Neovascularization Models for Evaluation of Angiogenic Activity of Vital Pulp Therapy Biomaterials, *Journal of endodontics*, 2018, 44(5), pp.773-779.
16. Pushparajan, C., et al., A mechanically strengthened polyacrylamide gel matrix fully compatible with electrophoresis of proteins and nucleic acids, *Electrophoresis*, 2018, 39(5-6), pp.824-832.
17. Gangolphe, L., et al., Degradable multi (aryl azide) star copolymer as universal photo-crosslinker for elastomeric scaffolds, *Materials Today Chemistry*, 2019, 12:209-21.
18. Peuler, K. et al., Clickable modular polysaccharide nanoparticles for selective cell-targeting. *Carbohydrate Polymers*. 2020; 234:115901.
19. Graf, M., et al., Hydrogel microspheres evading alveolar macrophages for sustained pulmonary protein delivery, *International Journal of Pharmaceutics*, 2019, 566, p.652-661.
20. Chu, S., et al., Cell encapsulation spatially alters crosslink density of poly(ethylene glycol) hydrogels formed from free-radical polymerizations, *Acta Biomaterialia*, 2020, 109, p. 37-50.

21. Wang, C., et al., A comparative study of brain tumor cells from different age and anatomical locations using 3D biomimetic hydrogels. *Acta Biomaterialia*, 2020.
22. Li, Y., et al., Matrix metalloproteinase (MMP)-degradable tissue engineered periosteum coordinates allograft healing via early stage recruitment and support of host neurovasculature, *Biomaterials*, 2021, V. 268.
23. Ziegler, C. E., et al., A novel anhydrous preparation of PEG hydrogels enables high drug loading with biologics for controlled release applications, *European Polymer Journal*, 2021, V. 147.
24. Crocini, C., et al., Three-dimensional encapsulation of adult mouse cardiomyocytes in hydrogels with tunable stiffness, *Progress in Biophysics and Molecular Biology*, 2019.
25. Peuler, K., et al., Clickable modular polysaccharide nanoparticles for selective cell-targeting, *Carbohydrate Polymers*, 2020, V. 234.
26. Wang, C., et al., Mimicking brain tumor-vasculature microanatomical architecture via co-culture of brain tumor and endothelial cells in 3D hydrogels, *Biomaterials*, 2019.

PEG GPC STANDARDS

JenKem Technology provides a Methoxy PEG GPC Calibration Standards Kit (STDMPG) containing linear methoxy PEGs with MW of 2kDa, 5kDa, 10kDa, 20kDa and 40kDa. JenKem Technology also provides a 4ARM PEG GPC Calibration Standards Kit (STD4ARMPEG) that includes 200 mg each of 4ARM-PEG-2000, 4ARM-PEG-5000, 4ARM-PEG-10K, 4ARM-PEG-20K, and 4ARM-PEG-40K.

The quality of the PEG standards is very important for GPC characterization of PEGs and PEGylated products. PEG GPC Calibration Standard Kits from JenKem Technology contain inactive hydroxyl-ended methoxy PEG and 4arm PEG raw material standards, with narrower polydispersity and more accurate molecular weight compared to other PEG standards on the market. The molecular weight of the JenKem® PEG standards is determined by MALDI-TOF. The polydispersity of the JenKem® PEG standards is determined via GPC.

ITEM CODE	ITEM NUMBER	ITEM DESCRIPTION	PRICE (USD)
STD001	STDMPG	Methoxy PEG GPC Standards Kit (200mg each MW 2000, 5000, 10000, 20000, 40000)	\$500
STD002	STD4ARMPEG	4ARM PEG GPC Standards Kit (200mg each MW 2000, 5000, 10000, 20000, 40000)	\$500

REFERENCES:

1. Croitoru-Sadger, T., et al., Two-component cross-linkable gels for fabrication of solid oral dosage forms, *Journal of Controlled Release*, 2019.
2. Abolmaali, S.S., et al., Nano-hydrogels of methoxy polyethylene glycol-grafted branched polyethyleneimine via biodegradable cross-linking of Zn²⁺-ionomer micelle template, *Journal of nanoparticle research*, 2013.
3. Prednisolone- α -cyclodextrin-star poly(ethylene glycol) polypseudorotaxane with delayed pH-sensitivity as a targeted drug delivery system, *International Journal of Pharmaceutics*, 2011, 414:1-2, 29, P. 42-47.
4. Prednisolone- α -cyclodextrin-star PEG polypseudorotaxanes with controlled drug delivery properties, *Org. Biomol. Chem.*, 2010, 8, 5423–5430.
5. Targeted antifungal delivery system: β -Glucosidase sensitive nystatin-star poly(ethylene glycol) conjugate, *International Journal of Pharmaceutics*, 2010, 386(1-2), p: 1–5.

ORDERING INFORMATION

To Place an Order

We encourage you to visit our website to place an online order for our catalog PEGs:
<http://www.jenkemusa.com/shop>

You can also place an order directly with us for catalog and custom PEGs, and PEGylation services, by phone, fax, email or mail. Please provide the following information in your purchase order: your name, phone number, and email address, complete shipping and billing addresses with contact phone name, phone number, purchase order number, payment method, product name, grade, and quantity.

Phone numbers:

866.953.6536 (866.9JENKEM)
1.972.673.0603

Fax:

469.327.2615

Email:

sales@jenkemusa.com

Address:

JenKem Technology USA, Inc., 4105 W. Spring Creek Pkwy Suite 606B, Plano, TX 75024, U. S. A.

Payment Information

Payments via postal service please remit to:

JenKem Technology USA, 4105 W. Spring Creek Pkwy Suite 606B, Plano, TX 75024, U. S. A.

Payment by bank draft or check must be made in U.S. dollars and must mention your purchase order number and invoice number. Please contact us for wire transfer instructions.

Visa, Mastercard, and American Express credit cards are accepted.

Custom Synthesis and PEGylation Services

We have the capability to provide many PEG compounds that are not listed in this catalog and invite your inquiries for custom synthesis. PEGylation services are provided for conjugation of our PEGs to proteins, peptides, oligonucleotides, and small molecules. Please call us or email us at tech@jenkemusa.com to discuss or request a quote for any custom synthesis and PEGylation services. Strict confidentiality can be arranged.

Bulk Orders and Commercial GMP Manufacture

Larger than catalog, bulk orders up through commercial scale manufacturing can be supplied at a lower price per gram than our catalog prices. Both non-GMP and GMP grade can be provided. Please contact us at tech@jenkemusa.com for quotations on bulk orders.

Technical and Customer Support

Our Technical Support Team is ready to answer your PEG-related questions. Please do not hesitate to contact us at tech@jenkemusa.com.

STORAGE AND HANDLING OF PEG DERIVATIVES

All PEG derivatives should be stored under an inert atmosphere (e.g., dry Nitrogen or dry Argon), in the dark, at or below -20°C.

PEG NHS Esters and Maleimides react with moisture in the atmosphere. PEG Maleimide, PEG Thiol, and PEG Acrylates are very sensitive to light, so it is important to keep these products in the dark in dry conditions all the times to avoid loss of activity. PEG aldehydes and PEG thiol can be deactivated by exposure to atmospheric oxygen. In addition, all PEG and PEG derivatives react with atmospheric oxygen to form peroxides along the backbone, a very slow process facilitated by light and heavy metals. The peroxides formed in this process can lead to chain cleavage, increasing the polydispersity of PEG and PEG derivatives.

Our PEG derivatives are thoroughly dried and packaged under an inert atmosphere. Generally, shipment of PEG derivatives at ambient temperature will not affect their activity. For long term storage, PEG derivatives should be stored at or below -20°C. Our standard shipping is at ambient temperature. Low temperature shipping (blue ice and dry ice) is also available upon request at an additional charge.

When using a PEG derivative after storage at low temperature, the material should be warmed up slowly to ambient temperature. After part of the material is taken out, the container with the remaining PEG derivative should be blanketed with an inert gas (e.g., dry Nitrogen or dry Argon) and stored again at or below -20°C.

For multiple uses, it is suggested that PEG containing compounds be repackaged into several containers to reduce the risk of loss of activity due to repetitive freeze-thaw cycles. We offer repacking services for our PEG products for an additional charge.

Please inquire if the PEG product that you need is not listed in this catalog. We may still have your compound, as not all accessible functional groups and molecular weight combinations are listed in the catalog. In addition, we have chemists dedicated to developing new products and processes on a custom synthesis basis, under confidentiality if needed, and we would be glad to provide a quote for you for this service. Please call or email us at tech@jenkemusa.com.